
Data Acquisition Toolbox™

Reference





How to Contact MathWorks

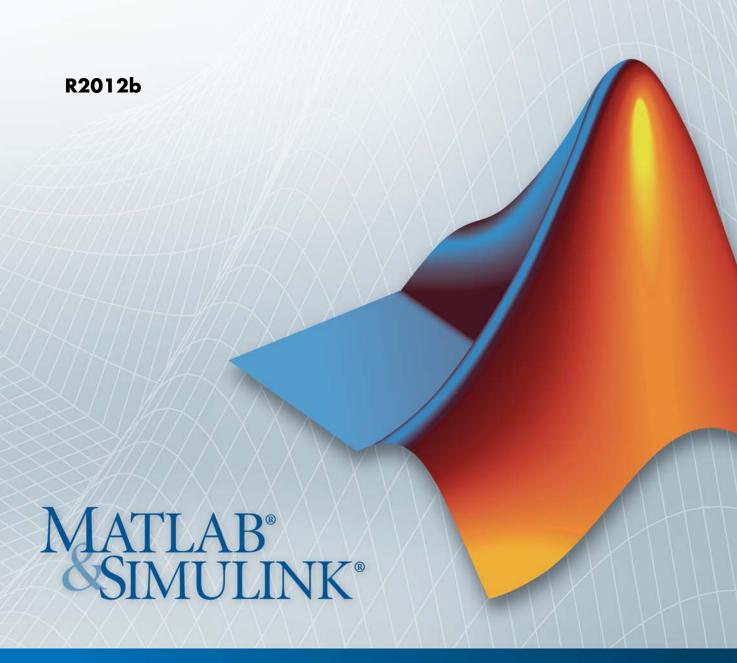

www.mathworks.com

Newsgroup comp.soft-sys.matlab

www.mathworks.com/contact TS.html Technical Support

suggest@mathworks.com bugs@mathworks.com

doc@mathworks.com

service@mathworks.com info@mathworks.com

Order status, license renewals, passcodes Sales, pricing, and general information

508-647-7000 (Phone)



508-647-7001 (Fax)



The MathWorks, Inc. 3 Apple Hill Drive Natick. MA 01760-2098

For contact information about worldwide offices, see the MathWorks Web site.

Data Acquisition ToolboxTM Reference

© COPYRIGHT 2005–2010 by The MathWorks, Inc.

The software described in this document is furnished under a license agreement. The software may be used or copied only under the terms of the license agreement. No part of this manual may be photocopied or reproduced in any form without prior written consent from The MathWorks, Inc.

Web

Bug reports

Product enhancement suggestions

Documentation error reports

FEDERAL ACQUISITION: This provision applies to all acquisitions of the Program and Documentation by, for, or through the federal government of the United States. By accepting delivery of the Program or Documentation, the government hereby agrees that this software or documentation qualifies as commercial computer software or commercial computer software documentation as such terms are used or defined in FAR 12.212, DFARS Part 227.72, and DFARS 252.227-7014. Accordingly, the terms and conditions of this Agreement and only those rights specified in this Agreement, shall pertain to and govern the use, modification, reproduction, release, performance, display, and disclosure of the Program and Documentation by the federal government (or other entity acquiring for or through the federal government) and shall supersede any conflicting contractual terms or conditions. If this License fails to meet the government's needs or is inconsistent in any respect with federal procurement law, the government agrees to return the Program and Documentation, unused, to The MathWorks, Inc.

Trademarks

MATLAB and Simulink are registered trademarks of The MathWorks, Inc. See www.mathworks.com/trademarks for a list of additional trademarks. Other product or brand names may be trademarks or registered trademarks of their respective holders.

Patents

MathWorks products are protected by one or more U.S. patents. Please see www.mathworks.com/patents for more information.

Revision History

September 2010	Online only	Revised for Version 2.17 (Release 2010b)
April 2011	Online only	Revised for Version 2.18 (Release 2011a)
September 2011	Online only	Revised for Version 3.0 (Release 2011b)
March 2012	Online only	Revised for Version 3.1(Release 2012a)
September 2012	Online only	Revised for Version 3.0 (Release 2012b)

Base Property Reference

1		
	Analog Input and Output and Digital I/O Properties	1-2
	Analog Input Properties	1-3
	Common Properties	1-3 1-6
	Analog Output Properties	1-8
	Common Properties	1-8 1-10
	Digital I/O Properties	1-12
	Common Properties	1-12
	Line Properties	1-12 1-13
	Session-Based Interface	1-14
_ [Base Properties — Alphabetical	List
2		
) [Device-Specific Proper	rties
	Groups by Vendor	3-2
	Advantagh	2_2

	Measurement Computing	3-4
	National Instruments	3-5
	Parallel Port	3-7
	Sound Card	3-8 3-8
a 1	Device-Specific Properties — Alphabetical 1	List
7		
	Block Refere	nce
5		
- 1	Function Refere	nce
5		
	Device Object Creation	6-2
	Channels and Lines	6-3
	Property Getting and Setting	6-4
	Object Execution	6-5
	Data Manipulation	6-6
	Information and Help	6-7

Utilities 6-	6-8
Session-Based Interface 6-1	10
Class Reference	ce
7	_
Functions — Alphahetical Lis	st
Functions — Alphabetical Lis	<u>st</u>
Functions — Alphabetical Lis	

Base Property Reference

Analog Input and Output and Digital I/O Properties

Base properties apply to all supported hardware subsystems of a given type (analog input, analog output, or digital I/O). For example, the SampleRate property is supported for all analog input subsystems regardless of the vendor.

The properties are categorized according to these subsystems:

Depending on the hardware device you are using, additional property names or property values might be present. For example, only analog input and analog output objects associated with a sound card have a BitsPerSample property. The additional property values are also device-specific but are included in this chapter. For example, all supported devices have an InputType property, but the value AC-Coupled is unique to analog input objects associated with a sound card.

Analog Input Properties

Common Properties

The analog input common properties are grouped into the following categories based on usage.

Analog Input Basic Setup Properties

SampleRate Specify per-channel rate at which

analog data is converted to digital

data, or vice versa

SamplesPerTrigger Specify number of samples to acquire

for each channel group member for

each trigger that occurs

TriggerType Specify type of trigger to execute

Analog Input Logging Properties

LogFileName Specify name of disk file information

is logged to

Logging Indicate whether data is being

logged to memory or disk file

LoggingMode Specify destination for acquired data

LogToDiskMode Specify whether data, events, and

hardware information are saved to

one or more disk files

Analog Input Trigger Properties

InitialTriggerTime Absolute time of first trigger

ManualTriggerHw0n Specify hardware device starts at

manual trigger

TriggerChannel Specify channel serving as trigger

source

TriggerCondition Specify condition that must be

satisfied before trigger executes

TriggerConditionValue Specify voltage value(s) that must be

satisfied before trigger executes

TriggerDelay Specify delay value for data logging
TriggerDelayUnits Specify units in which trigger delay

data is measured

TriggerFcn Specify callback function to execute

when trigger occurs

TriggerRepeat Specify number of additional times

trigger executes

TriggersExecuted Indicate number of triggers that

execute

TriggerType Specify type of trigger to execute

Analog Input Status Properties

Logging Indicate whether data is being

logged to memory or disk file

Running Indicate whether device object is

running

SamplesAcquired Indicate number of samples acquired

per channel

Samples Available Indicate number of samples available

per channel in engine

Analog Input Hardware Configuration Properties

ChannelSkew Specify time between consecutive

scanned hardware channels

ChannelSkewMode Specify how channel skew is

determined

ClockSource Specify clock that governs hardware

conversion rate

InputType Specify analog input hardware

channel configuration

SampleRate Specify per-channel rate at which

analog data is converted to digital

data, or vice versa

Analog Input Callback Properties

DataMissedFcn Specify callback function to execute

when data is missed

InputOverRangeFcn Specify callback function to execute

when acquired data exceeds valid

hardware range

RuntimeErrorFcn Specify callback function to execute

when run-time error occurs

SamplesAcquired Indicate number of samples acquired

per channel

SamplesAcquiredFcn Specify callback function to execute

when predefined number of samples is acquired for each channel group

member

SamplesAcquiredFcnCount Specify number of samples to acquire

for each channel group member before samples acquired event is

generated

StartFcn Specify callback function to execute

before device object runs

StopFcn Specify callback function to execute

after device object runs

TimerFcn Specify callback function to execute

when predefined time period passes

TimerPeriod Specify time period between timer

events

TriggerFcn Specify callback function to execute

when trigger occurs

Analog Input General Purpose Properties

BufferingConfig Specify per-channel allocated

memory

BufferingMode Specify how memory is allocated

Channel Contain hardware channels added to

device object

EventLog Store information for specific events

Name Specify descriptive name for the

channel

Tag Specify device object label

Timeout Specify additional waiting time to

extract or queue data

Type Indicate device object type, channel,

or line

UserData Store data to associate with device

object

Channel Properties

The analog input channel properties are given below.

Analog Input Channel Properties

Channel Name Specify descriptive channel name

HwChannel Specify hardware channel ID
HwLine Specify hardware line ID

Index MATLAB® index of hardware

channel or line

InputRange Specify range of analog input

subsystem

NativeOffset Indicate offset to use when

converting between native data

format and doubles

NativeScaling Indicate scaling to use when

converting between native data

format and doubles

Parent Indicate parent (device object) of

channel or line

SensorRange Specify range of data expected from

sensor

Type Indicate device object type, channel,

or line

Units Specify engineering units label

UnitsRange Specify range of data as engineering

units

Analog Output Properties

Common Properties

The analog output common properties are grouped into the following categories based on usage.

Analog Output Basic Setup Properties

SampleRate Specify per-channel rate at which

analog data is converted to digital

data, or vice versa

TriggerType Specify type of trigger to execute

Analog Output Trigger Properties

InitialTriggerTime Absolute time of first trigger

TriggerFcn Specify callback function to execute

when trigger occurs

TriggersExecuted Indicate number of triggers that

execute

TriggerType Specify type of trigger to execute

Analog Output Status Properties

Running Indicate whether device object is

running

Samples Available Indicate number of samples available

per channel in engine

SamplesOutput Indicate number of samples output

per channel from engine

Sending Indicate whether data is being sent

to hardware device

Analog Output Hardware Configuration Properties

ClockSource Specify clock that governs hardware

conversion rate

SampleRate Specify per-channel rate at which

analog data is converted to digital

data, or vice versa

Analog Output Data Management Properties

MaxSamplesQueued Indicate maximum number of

samples that can be gueued in

engine

RepeatOutput Specify number of additional times

queued data is output

Timeout Specify additional waiting time to

extract or queue data

Analog Output Callback Properties

RuntimeErrorFcn Specify callback function to execute

when run-time error occurs

SamplesOutputFcn Specify callback function to execute

when predefined number of samples is output for each channel group

member

SamplesOutputFcnCount Specify number of samples to output

for each channel group member before samples output event is

generated

StartFcn Specify callback function to execute

before device object runs

StopFcn Specify callback function to execute

after device object runs

TimerFcn Specify callback function to execute

when predefined time period passes

TimerPeriod Specify time period between timer

events

TriggerFon Specify callback function to execute

when trigger occurs

Analog Output General Purpose Properties

BufferingConfig Specify per-channel allocated

memory

BufferingMode Specify how memory is allocated

Channel Contain hardware channels added to

device object

EventLog Store information for specific events

Name Specify descriptive name for the

channel

OutOfDataMode Specify how value held by analog

output subsystem is determined

Tag Specify device object label

Type Indicate device object type, channel,

or line

UserData Store data to associate with device

object

Channel Properties

The analog output channel properties are given below.

Analog Output Channel Properties

ChannelName Specify descriptive channel name

DefaultChannelValue Specify value held by analog output

subsystem

HwChannel Specify hardware channel ID

Index MATLAB index of hardware channel

or line

NativeOffset Indicate offset to use when

converting between native data

format and doubles

NativeScaling Indicate scaling to use when

converting between native data

format and doubles

OutputRange Specify range of analog output

hardware subsystem

Parent Indicate parent (device object) of

channel or line

Type Indicate device object type, channel,

or line

Units Specify engineering units label

UnitsRange Specify range of data as engineering

units

Digital I/O Properties

Common Properties

The digital I/O common properties are given below.

Digital I/O Common Properties

Line Contain hardware lines added to

device object

Name Specify descriptive name for the

channel

Running Indicate whether device object is

running

Tag Specify device object label

TimerFcn Specify callback function to execute

when predefined time period passes

TimerPeriod Specify time period between timer

events

Type Indicate device object type, channel,

or line

UserData Store data to associate with device

object

Line Properties

The digital I/O line properties are given below.

Digital I/O Line Properties

Direction Specify whether line is for input or

output

HwLine Specify hardware line ID

Index MATLAB index of hardware channel

or line

LineName Specify descriptive line name

Parent Indicate parent (device object) of

channel or line

Port Specify port ID

Type Indicate device object type, channel,

or line

Getting Command-Line Property Help

To get command-line property help, use the daqhelp function. For example, to get help for the SampleRate property, type:

daqhelp SampleRate

Note You can use daghelp without creating a device object.

You can also get property characteristics, such as the default property value, using the propinfo function. For example, suppose you create the analog input object ai for a sound card and you want to find the default value for the SampleRate property.

```
ai = analoginput('winsound');
out = propinfo(ai, 'SampleRate');
out.DefaultValue
ans =
    8000
```

Session-Based Interface

Properties methods used with the session-based interface are listed below.

ActiveEdge Rising or falling edges of EdgeCount

signals

Active Pulse Massurement of

 $Pulse Width\ counter\ channel$

ADCTimingMode Set channel timing mode

AutoSyncDSA Automatically Synchronize DSA

devices

BridgeMode Specify analog input device bridge

mode

Channels Array of channel objects associated

with session object

Connections Array of connections in a session

CountDirection Specify direction of counter channel
Coupling Specify input coupling mode

Destination Indicates trigger destination

terminal

Device annel device

Direction Specify digital channel direction
DurationInSeconds Specify duration of acquisition

DutyCycle Duty cycle of counter output channel
EncoderType Encoding type of counter channel

ExcitationCurrent Voltage of external source of

excitation

ExcitationSource External source of excitation

ExcitationVoltage Voltage of excitation source

ExternalTriggerTimeout Indicate if external trigger timed out

Frequency of generated pulses on

counter output channel

ID ID

IdleState Default state of counter output

channel

InitialCount Specify initial count point

InitialDelay Delay until output channel generates

pulses

IsContinuous Specify if operation continues until

manually stopped

IsDone Indicate if operation is complete

IsLogging Indicate if hardware is acquiring or

generating data

IsNotifyWhenDataAvailableExceedsAutontrol if is set automatically

IsNotifyWhenScansQueuedBelowAuto Control if is set automatically

IsRunning Indicate if operation is still in

progress

IsSimulated Indicate if device is simulated

IsWaitingForExternalTrigger Indicates if synchronization is

waiting for an external trigger

MaxSoundPressureLevel Sound pressure level for microphone

channels

MeasurementType Type counter channel measurement

Name Specify descriptive name for the

channel

Name Specify descriptive name for the

channel

NominalBridgeResistance Resistance of sensor

NotifyWhenDataAvailableExceeds Control firing of DataAvailable

event

NotifyWhenScansQueuedBelow Control firing of DataRequired event

Number Of Scans for operation when

starting

R0 Specify resistance value

Range Specify channel measurement range
Rate Rate of operation in scans per second
RateLimit Limit of rate of operation based on

hardware configuration

RTDConfiguration Specify wiring configuration of RTD

device

RTDType Specify sensor sensitivity

ScansAcquired Number of scans acquired during

operation

ScansOutputByHardware Indicate number of scans output by

hardware

ScansQueued Indicate number of scans queued for

output

Sensitivity Sensitivity of an channel

ShuntLocation Indicate location of channel's shunt

resistor

ShuntResistance Resistance value of channel's shunt

resistor

Source Indicates trigger source terminal

Terminal PFI terminal

TerminalConfig Specify terminal configuration

Terminals Terminals available on device or

CompactDAQ chassis

ThermocoupleType Select thermocouple type

TriggerCondition Specify condition that must be

satisfied before trigger executes

TriggersPerRun Indicate the number of times the

trigger executes in an operation

TriggersRemaining Indicates the number of trigger to

execute in an operation

TriggerType Type of trigger executed

Type Display synchronization trigger type
Units Specify unit of RTD measurement
Vendor Vendor information associated with

session object

ZResetCondition Reset condition for Z-indexing
ZResetEnable Enable reset for Z-indexing
ZResetValue Reset value for Z-indexing

Base Properties — Alphabetical List

ActiveEdge

Purpose Rising or falling edges of EdgeCount signals **Description** When working with the session-based interface, use the ActiveEdge property to represent rising or falling edges of a EdgeCount signal. **Values** You can set the Active edge of a counter input channel to Rising or Falling. **Examples** s = daq.createSession('ni'); ch = s.addCounterInputChannel ('cDAQ1Mod5', 0, 'EdgeCount') ch = Data acquisition counter input edge count channel 'ctr0' on device 'Dev2' ActiveEdge: Rising CountDirection: Increment InitialCount: 0 Terminal: 'PFI8' Name: empty ID: 'ctr0' Device: [1x1 daq.ni.DeviceInfo] MeasurementType: 'EdgeCount' Change the Active Edge property to 'Falling': ch.ActiveEdge = 'Falling' ch = Data acquisition counter input edge count channel 'ctr0' on device 'Dev2': ActiveEdge: Falling CountDirection: Increment InitialCount: 0 Terminal: 'PFI8' Name: empty

ID: 'ctr0'

Device: [1x1 daq.ni.DeviceInfo]

MeasurementType: 'EdgeCount'

See Also Methods

daq.Session.addCounterInputChannel,
daq.Session.addCounterOutputChannel

Class

daq.Session

Purpose

Active pulse measurement of PulseWidth counter channel

Description

When working with the session-based interface, the ActivePulse property displays the pulse width measurement in seconds of your counter channel, with PulseWidth measurement type.

Values

Active pulse measurement values include:

- 'High'
- 'Low'

Examples

Create a session object, add a counter input channel, with the 'EdgeCount' MeasurementType.

```
s = daq.createSession('ni');
ch = s.addCounterInputChannel ('cDAQ1Mod5', 0, 'PulseWidth')

ch =

Data acquisition counter input pulse width channel 'ctr0' on device 'cDAQ1Mod5':

        ActivePulse: High
        Terminal: 'PFI4'
        Name: empty
        ID: 'ctr1'
        Device: [1x1 daq.ni.DeviceInfo]

MeasurementType: 'PulseWidth

Change the ActiveEdge property to Low.

ch.ActivePulse = 'Low'

ch =

Data acquisition counter input pulse width channel 'ctr0' on device 'cDAQ1Mod5':
```

ActivePulse: Low

ActivePulse

Terminal: 'PFI4'
Name: empty
ID: 'ctr1'

Device: [1x1 daq.ni.DeviceInfo]

MeasurementType: 'PulseWidth'

See Also Class

 ${\tt daq.Session.addCounterInputChannel}$

ADCTimingMode

Purpose

Set channel timing mode

Description

When working with the session-based interface, use the ADCTimingMode property to specify if the timing mode in of all channels in the device is high resolution or high speed.

Note The ADCTimingMode must be the same for all channels on the device.

Values

You can set the ADCTimingMode to:

s = daq.createSession('ni');

- 'HighResolution'
- 'HighSpeed'
- 'Best50HzRejection'
- 'Best60HzRejection'

Examples

Create a session and add an analog input channel:

```
s.addAnalogInputChannel('cDAQ1Mod1','ai1','Voltage');
s.Channels
ans =
```

Data acquisition analog input voltage channel 'ai1' on device 'cDAQ1Mod1'

```
Coupling: DC
TerminalConfig: SingleEnded
Range: -10 to +10 Volts
Name: ''
ID: 'ai1'
Device: [1x1 daq.ni.CompactDAQModule]
MeasurementType: 'Voltage'
```

ADCTimingMode

ADCTimingMode: ''

Set the ADCTimingMode property to 'HighResolution':

s.Channels.ADCTimingMode = 'HighResolution';

See Also Methods

daq.Session.addAnalogInputChannel

Class

daq.Session

AutoSyncDSA

Purpose Automatically Synchronize DSA devices

Description Use this property to enable or disable automatic synchronization

between DSA (PXI or PCI) devices in the same session. By default the

sessions automatic synchronization capability is disabled.

Examples To enable automatic synchronization, create a session and add channels

from a DSA device:

```
s=daq.createSession('ni')
s.addAnalogInputChannel('PXI1Slot2',0,'Voltage');
s.addAnalogInputChannel('PXI1Slot3',1,'Voltage');
```

Enable automatic synchronization and acquire data"

```
s.AutoSyncDSA=true;
s.startForeground;
```

See Also daq.Session.addAnalogInputChannel

Purpose

Specify analog input device bridge mode

Description

Use this property in the session-based interface to specify the bridge mode, which represents the active gauge of the analog input channel.

The bridge mode is 'Unknown' when you add a bridge channel to the session. Change this value to a valid mode to use the channel. Valid bridge modes are:

- 'Full' All four gauges are active.
- 'Half'— Only two bridges are active.
- 'Quarter'— Only one bridge is active.

See Also Class

daq.Session

BufferingConfig

Purpose

Specify per-channel allocated memory

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

BufferingConfig is a two-element vector that specifies the per-channel allocated memory. The first element of the vector specifies the block size, while the second element of the vector specifies the number of blocks. The total allocated memory (in bytes) is given by

(block size).(number of blocks).(number of channels).(native data type)

You can determine the native data type with daghwinfo.

You can allocate memory automatically or manually. If BufferingMode is Auto, the BufferingConfig values are automatically set by the engine. If BufferingMode is Manual, then you must manually set the BufferingConfig values. If you change the BufferingConfig values, BufferingMode is automatically set to Manual.

When memory is automatically allocated by the engine, the block-size value depends on the sampling rate and is typically a binary number. The number of blocks is initially set to a value of 30 but can dynamically increase to accommodate the memory requirements. In most cases, the number of blocks used results in a per-channel memory that is somewhat greater than the SamplesPerTrigger value. When you manually allocate memory, the number of blocks is not dynamic and care must be taken to ensure there is sufficient memory to store the acquired data. If the number of samples acquired or queued exceeds the allocated memory, then an error is returned.

You can easily determine the memory allocated and available memory for each device object with the dagmem function.

Characteristics

Usage AI, AO, common to all channels

Access Read/write

Data type Two-element vector of doubles

Yes

Read-only when

running

Values

The default value is determined by the engine, and is based on the number of channels contained by the device object and the sampling rate. The BufferingMode value determines whether the values are automatically updated as data is acquired. For analog output objects, the default number of blocks is two.

Note If you change the BufferingConfig property for an analog output object, all previously queued output data will get discarded.

Examples

Create the analog input object ai for a sound card and add two channels to it.

```
ai = analoginput('winsound');
addchannel(ai,1:2);
```

The block size and number of blocks are given by BufferingConfig, while the native data type for the sound card is given by daghwinfo.

```
ai.BufferingConfig
ans =
    512    30
out = daqhwinfo(ai);
out.NativeDataType
ans =
int16
```

BufferingConfig

With this information, the total allocated memory is calculated to be 61,440 bytes. This number is stored by dagmem.

```
out = daqmem(ai);
out.UsedBytes
ans =
61440
```

The allocated memory is more than sufficient to store 8000 two-byte samples for two channels. If more memory was required, then the number of blocks would dynamically grow because BufferingMode is set to Auto.

See Also Functions

daghwinfo, dagmem

Properties

BufferingMode, SampleRate, SamplesPerTrigger

Specify how memory is allocated

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

BufferingMode can be set to Auto or Manual. If BufferingMode is set to Auto, the data acquisition engine automatically allocates the required memory. If BufferingMode is set to Manual, you must manually allocate memory with the BufferingConfig property.

If BufferingMode is set to Auto and the SampleRate value is changed, then the BufferingConfig values might be recalculated by the engine. Specifically, you can increase (decrease) the block size if SampleRate is increased (decreased). If BufferingMode is set to Auto and you change the BufferingConfig values, then BufferingMode is automatically set to Manual. If BufferingMode is set to Manual, then you cannot set the number of blocks to a value less than three.

For most data acquisition applications, you should set BufferingMode to Auto and have memory allocated by the engine because this minimizes the chance of an out-of-memory condition.

Characteristics

Usage AI, AO, common to all channels

Access Read/write

Data type String
Read-only when Yes

running

Values

{Auto} Memory is allocated by the data acquisition engine.

Manual Memory is allocated manually.

BufferingMode

See Also Functions

daqmem

Properties

BufferingConfig

Contain hardware channels added to device object

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Channel is a vector of all the hardware channels contained by an analog input (AI) or analog output (AO) object. Because a newly created AI or AO object does not contain hardware channels, Channel is initially an empty vector. The size of Channel increases as channels are added with the addchannel function, and decreases as channels are removed using the delete function.

Channel is used to reference one or more individual channels. To reference a channel, you must know its MATLAB index, which is given by the Index property. For example, you must use Channel with the appropriate indices when configuring channel property values.

For scanning hardware, the scan order follows the MATLAB index. Therefore, the hardware channel associated with index 1 is sampled first, the hardware channel associated with index 2 is sampled second, and so on. To change the scan order, you can specify a permutation of the indices with Channel.

Characteristics

Usage AI, AO

Access Read/write

Data type Vector of channels

Read-only when Yes

running

Values

Values are automatically defined when channels are added to the device object with the addchannel function. The default value is an empty column vector.

Examples

Create the analog input object ai for a National Instruments® card and add three hardware channels to it.

```
ai = analoginput('nidaq','Dev1');
addchannel(ai,0:2);
```

To set a property value for the first channel added (ID = 0), you must reference the channel by its index using the Channel property.

```
chans = ai.Channel(1);
set(chans, 'InputRange', [-10 10])
```

Based on the current configuration, the hardware channels are scanned in order from 0 to 2. To swap the scan order of channels 0 and 1, you can specify the appropriate permutation of the MATLAB indices with Channel.

```
ai.Channel([1 2 3]) = ai.Channel([2 1 3]);
```

See Also

Functions

addchannel, delete

Properties

HwChannel, Index

Specify descriptive channel name

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

ChannelName specifies a descriptive name for a hardware channel. If a channel name is defined, then you can reference that channel by its name. If a channel name is not defined, then the channel must be referenced by its index. Channel names are not required to be unique.

You can also define descriptive channel names when channels are added to a device object with the addchannel function.

Characteristics

Usage AI, AO, per channel

Access Read/write
Data type String
Read-only when Yes

running

Values

The default value is an empty string. To reference a channel by name, it must contain only letters, numbers, and underscores and must begin with a letter.

Examples

Create the analog input object ai for a sound card and add two channels to it.

```
ai = analoginput('winsound');
addchannel(ai,1:2);
```

To assign a descriptive name to the first channel contained by ai:

```
Chan1 = ai.Channel(1)
set(Chan1, 'ChannelName', 'Joe')
```

ChannelName

You can now reference this channel by name instead of by index.

set(ai.Joe, 'Units', 'Decibels')

See Also Functions

addchannel

Array of channel objects associated with session object

Description

This session object property contains and displays an array of channels added to the session. For more information on the session-based interface, see "Session-Based Interface".

Tip You cannot directly add or remove channels using the Channels object properties. Use daq.Session.addAnalogInputChannel and daq.Session.addAnalogOutputChannel to add channels. Use daq.Session.removeChannel to remove channels.

Values

The value is determined by the channels you add to the session object.

Examples

Create a session object, add an analog input channel, and display the session Channels property.

```
s = daq.createSession('ni');
s.addAnalogInputChannel('cDAQ1Mod1','ai1','Voltage');
s.Channels
ans =
Data acquisition analog input channel 'ai1' on device 'cDAQ1Mod1':
```

```
Coupling: DC
InputType: Differential
Range: -10 to +10 Volts
Name: empty
ID: 'ai1'
Device: [1x1 daq.ni.CompactDAQModule]
ADCTimingMode: empty
```

Add an analog output channel and view the Channels property:

```
s.addAnalogOutputChannel('cDAQ1Mod2', 'ao1', 'Voltage');
```

Channels

Change the InputType property of the input channel to SingleEnded:

s.Channels(1).InputType = 'SingleEnded';

See Also Methods

daq.Session.addAnalogInputChannel,
daq.Session.addAnalogOutputChannel

Class

daq.Session

Specify time between consecutive scanned hardware channels

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

ChannelSkew applies only to scanning hardware and not to simultaneous sample and hold (SS/H) hardware.

If ChannelSkewMode is set to Minimum or Equisample, then ChannelSkew is automatically set to the appropriate device-specific read-only value. For SS/H hardware, the only valid ChannelSkew value is zero. For some vendors, ChannelSkewMode is automatically set to Manual if you first set ChannelSkew to a valid value.

Characteristics

Usage AI, common to all channels

Access Read/write (depends on ChannelSkewMode value)

Data type Double Read-only when Yes

running

Values

For SS/H hardware, the only valid value is zero. For scanning hardware, the value depends on ChannelSkewMode. ChannelSkew is specified in seconds.

See Also Properties

ChannelSkewMode

ChannelSkewMode

Purpose

Specify how channel skew is determined

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

For simultaneous sample and hold (SS/H) hardware, ChannelSkewMode is None. For scanning hardware, ChannelSkewMode can be Minimum, Equisample, or Manual (National Instruments only). SS/H hardware includes sound cards, while scanning hardware includes most Measurement ComputingTM and NI boards. Note that some supported boards from these vendors are SS/H, such as Measurement Computing's PCI-DAS4020/12.

If ChannelSkewMode is Minimum, then the minimum channel skew supported by the hardware is used. Some vendors refer to this as burst mode. If ChannelSkewMode is Equisample, the channel skew is given by [(sampling rate)(number of channels)]-1. If ChannelSkewMode is Manual, then you must specify the channel skew with the ChannelSkew property. For some vendors, ChannelSkewMode is automatically set to Manual if you first set ChannelSkew to a valid value.

Notes If you want to use the maximum sampling rate of your hardware, you should set ChannelSkewMode to Equisample.

Large loads on the input device, especially if you are using multiple channels with scanning hardware, can increase the settling time. To improve the settling time, set ChannelSkewMode to Equisample and lower your sample rate.

Characteristics

Usage AI, common to all channels

Access Read/write

ChannelSkewMode

Data type String
Read-only when Yes

running

Values Advantech®

{Equisample} The channel skew is given by [(sampling

rate)(number of channels)]-1.

Measurement Computing

{Minimum} The channel skew is set to the minimum supported

value.

Equisample The channel skew is given by [(sampling

rate)(number of channels)]⁻¹.

National Instruments

{Minimum} The channel skew is set to the minimum supported

value.

Equisample The channel skew is given by [(sampling

rate)(number of channels)]-1.

Manual The channel skew is given by Channel Skew.

Sound Cards

{None} This is the only supported value for SS/H

hardware.

Examples Create an analog input object for an MCC device and add eight channels.

ai = analoginput('mcc',1);
addchannel(ai,0:7);

ChannelSkewMode

Using the default ChannelSkewMode value of Min and the default SampleRate value of 1000, the corresponding ChannelSkew value is

```
ai.ChannelSkew
ans =
  1.0000e-005
```

To use the maximum sampling rate, set ChannelSkewMode to Equisample.

```
ai.ChannelskewMode = 'Equisample';
ai.Samplerate = 100000/8;
```

See Also Properties

ChannelSkew, SampleRate

Specify clock that governs hardware conversion rate

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

For all supported hardware except Measurement Computing analog output subsystems, ClockSource can be set to Internal, which specifies that the acquisition rate is governed by the internal hardware clock.

Use this table to map to the National Instruments terminology.

Data Acquisition Toolbox™	NI_DAQmx
Scan Clock	Sample Clock
Sample Clock	Convert Clock

For subsystems without a hardware clock, you must use software clocking to govern the sampling rate. Software clocking allows a maximum sampling rate of 500 Hz and a minimum sampling rate of 0.0002 Hz. An error is returned if more than 1 sample of jitter is detected. Note that you might not be able to attain rates over 100 Hz on all systems.

Characteristics

Usage AI, AO, common to all channels

Access Read/write

Data type String
Read-only when Yes

running

ClockSource

Values Advantech

{Internal} The internal hardware clock is used

(AI only).

External Externally control the channel clock

(AI only).

Software The computer clock is used.

Measurement Computing

{Internal}ExternalExternalExternally control the channel clock.

Software The computer clock is used.

National Instruments

{Internal} The internal hardware clock is used.

External Externally control the channel clock (AO only).

ExternalSampleCtrlExternally control the channel clock. This value overrides the ChannelSkew property value (AI only). This value does not apply to cards with

simultaneous sample and hold.

Note If you set ClockSource to ExternalSampleCtrl then the value of ExternalSampleClockSource specifies the pin whose signal is used as the channel clock for conversions on each channel.

ExternalScanCtrl Externally control the scan clock. This value overrides the SampleRate property value (AI only).

> Note If you set ClockSource to ExternalScanCtrl then the value of ExternalScanClockSource specifies the pin whose signal is used as the scan clock to initiate conversions across a group of channels.

ExternalSampleAndSExtectrally control the channel and scan clocks.

This value overrides the ChannelSkew and SampleRate property values (AI only). This value does not apply to cards with simultaneous sample and hold.

Note If you set ClockSource to ExternalSampleAndScanCtrl then the value of ExternalSampleClockSource specifies the pin whose signal is used as the channel clock for conversions on each channel, and the value of ExternalScanClockSource specifies the pin whose signal is used as the scan clock to initiate conversions across a group of channels.

Note If you set the ClockSource property to one of the External options, you must also set the SampleRate property to a value close to the external clock rate. SampleRate does not directly affect the external device, and the device will not use SampleRate if you have set an external clock rate, but this ensures that the toolbox configures itself correctly for expected data rates.

ClockSource

Sound Cards

{Internal}

The internal hardware clock is used.

See Also Properties

ChannelSkew, SampleRate

Array of connections in a session

Description

This session property contains and displays all connections added to the session.

Tip You cannot directly add or remove connections using the Connections object properties. Use daq.Session.addTriggerConnection and daq.Session.addClockConnection to add connections. Use daq.Session.removeConnection to remove connections.

Values

The value is determined by the connections you add to the session.

Examples

This example shows you how to remove a synchronization connection.

Create a session and add analog input channels and trigger and clock connections.

```
s = daq.createSession('ni')
s.addAnalogInputChannel('Dev1', 0, 'voltage');
s.addAnalogInputChannel('Dev2', 0, 'voltage');
s.addAnalogInputChannel('Dev3', 0, 'voltage');
s.addTriggerConnection('Dev1/PFI4','Dev2/PFI0','StartTrigger');
s.addTriggerConnection('Dev1/PFI4','Dev3/PFI0','StartTrigger');
s.addClockConnection('Dev1/PFI5','Dev2/PFI1','ScanClock');
```

Examine the session Connections property.

s.Connections

Remove the last clock connection at index 4 and display the session connections.

- s.removeConnection(4)
- s.Connections

```
ans =
```

Scan Clock is provided by 'Dev1' at 'PFI5' and will be received by 'Dev2'

```
index    Type    Source    Destination

1    StartTrigger    Dev1/PFI4    Dev2/PFI0
2    StartTrigger    Dev1/PFI4    Dev3/PFI0
3    ScanClock    Dev1/PFI5    Dev2/PFI1
```

See Also Methods

daq.Session.addTriggerConnection
daq.Session.addClockConnection

Class

daq.Session

Specify direction of counter channel

Description

When working with the session-based interface, use the CountDirection property to set the direction of the counter. Count direction can be 'Increment', in which case the counter operates in incremental order, or 'Decrement', in which the counter operates in decremental order.

Examples

Create a session object, add a counter input channel, and change the CountDirection.

```
s = daq.createSession('ni');
ch = s.addCounterInputChannel ('cDAQ1Mod5', 0, 'EdgeCount')
ch =
Data acquisition counter input edge count channel 'ctr0' on device 'cDAQ1Mod5':
     ActiveEdge: Rising
 CountDirection: Increment
   InitialCount: 0
       Terminal: 'PFI8'
          Name: empty
            ID: 'ctr0'
        Device: [1x1 daq.ni.DeviceInfo]
MeasurementType: 'EdgeCount'
Change CountDirection to 'Decrement':
ch.CountDirection = 'Decrement'
ch =
Data acquisition counter input edge count channel 'ctr0' on device 'cl
       ActiveEdge: Rising
  CountDirection: Decrement
     InitialCount: 0
```

CountDirection

Terminal: 'PFI8'
Name: empty
ID: 'ctr0'

Device: [1x1 daq.ni.DeviceInfo]

MeasurementType: 'EdgeCount'

See Also Class

daq.Session, daq.Session.addCounterInputChannel

Specify callback function to execute when data is missed

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

A data missed event is generated immediately after acquired data is missed. This event executes the callback function specified for DataMissedFcn. The default value for DataMissedFcn is daqcallback, which displays the event type and the device object name.

In most cases, data is missed because:

- The engine cannot keep up with the rate of acquisition.
- The driver wrote new data into the hardware's FIFO buffer before the previously acquired data was read. You can usually avoid this problem by increasing the size of the memory block with the BufferingConfig property.

Data missed event information is stored in the Type and Data fields of the EventLog property. The Type field value is DataMissed. The Data field values are given below.

Data Field Value	Description
AbsTime	The absolute time (as a clock vector) the event occurred.
RelSample	The acquired sample number when the event occurred.

When a data missed event occurs, the analog input object is automatically stopped.

DataMissedFcn

 $\begin{array}{ccc} \textbf{Characteristics} & U_{sage} & & AI, \ common \ to \ all \ channels \end{array}$

Access Read/write

Data type String
Read-only when No

running

Values The default value is daqcallback.

See Also Functions

daqcallback

Properties

EventLog

Specify value held by analog output subsystem

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

DefaultChannelValue specifies the value to write to the analog output (AO) subsystem when data is finished being output from the engine.

DefaultChannelValue is used only when OutOfDataMode is set to DefaultValue. This property guarantees that a known value is held by the AO subsystem if a run-time error occurs. Note that sound cards do not have an OutOfDataMode property.

Characteristics

Usage AO, per channel

Access Read/write
Data type Double
Read-only when Yes

running

Values

The default value is zero.

Examples

Create the analog output object ao and add two channels to it.

```
ao = analogoutput('nidaq','Dev1');
addchannel(ao,0:1);
```

You can configure **ao** so that when it stops outputting data, a value of 1 volt is held for both channels.

```
ao.OutOfDataMode = 'DefaultValue';
ao.Channel.DefaultChannelValue = 1.0;
```

DefaultChannelValue

See Also Properties

OutOfDataMode

Destination

Purpose Indicates trigger destination terminal

Description When working with the session-based interface, the Source property

indicates the device and terminal to which you connect a trigger.

See Also Sourcedaq.Session.addTriggerConnection

Channel device information

Description

When working with the session-based interface, the read-only Device property displays device information for the channel.

Examples

Create a session object, add a counter input channel, and view the Device property.

```
s = daq.createSession('ni');
ch = s.addCounterInputChannel('cDAQ1Mod5', 0, 'EdgeCount');
ch.Device

ans =

ni cDAQ1Mod5: National Instruments NI 9402
   Counter input subsystem supports:
    Rates from 0.1 to 80000000.0 scans/sec
   2 channels
   'EdgeCount','PulseWidth','Frequency','Position' measurement types

Counter output subsystem supports:
   Rates from 0.1 to 80000000.0 scans/sec
   3 channels
   'PulseGeneration' measurement type
This module is in chassis 'cDAQ1', slot 5
```

See Also

Class

daq.Session, daq.Session.addCounterInputChannel,
daq.Session.addCounterOutputChannel

Specify digital channel direction

Description

When you add a digital channel or a group to a session, you can specify the measurement type to be:

- Input
- Output
- Unknown

When you specify the MeasurementType as Bidirectional, you can use the channel to input and output messges. By default the channel is set to Unknown. Change the direction to output singnal on the channel.

Example

To change the direction of a bidirectional signal on a digital channel in the session **s**, type:

s.Channels(1).Direction='Output';

Direction

Purpose

Specify whether line is for input or output

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

When adding hardware lines to a digital I/O object with addline, you must configure the line direction. The line direction can be In or Out, and is automatically stored in Direction. If a line direction is In, you can only read a value from that line. If a line direction is Out, you can write or read a line value.

For line-configurable devices, you can change individual line directions using Direction. For port-configurable devices, you cannot change individual line directions.

Characteristics

Usage DIO, per line
Access Read/write
Data type String

Read-only when

running

Yes

Values

{In} The line can be read from.

Out The line can be read from or written to.

Examples

Create the digital I/O object dio and add two input lines and two output lines to it.

```
dio = digitalio('nidaq','Dev1');
addline(dio,0:3,{'In','In','Out','Out'});
```

Direction

```
To configure all lines for output:
```

dio.Line(1:2).Direction = 'Out';

See Also Functions

addline

DurationInSeconds

Purpose

Specify duration of acquisition

Description

When working with the session-based interface, use the DurationInSeconds property to change the duration of an acquisition. When the session contains output channels, DurationInSeconds becomes a read only property and its value is determined by

```
\frac{s.ScansQueued}{s.Rate}.
```

Values

In a session with only input channels, you can enter a value in seconds for the length of the acquisition. Changing the duration changes the number of scans accordingly. By default, the DurationInSeconds is set to 1 second.

Examples

Create a session object, add an analog input channel, and change the duration:

See Also Properties

NumberOfScans, Rate

Class

daq.Session, daq.Session.addCounterInputChannel

Duty cycle of counter output channel

Description

When working with the session-based interface, use the DutyCycle property to specify the fraction of time that the generated pulse is in active state.

Duty cycle is the ratio between the duration of the pulse and the pulse period. For example, if a pulse duration is 1 microsecond and the pulse period is 4 microseconds, the duty cycle is 0.25. In a square wave, you will see that the time the signal is high is equal to the time the signal is low.

Examples

Create a session object and add a 'PulseGeneration' counter output channel:

```
s = daq.createSession('ni');
s.addCounterOutputChannel('cDAQ1Mod5', 'ctr0', 'PulseGeneration')
```

Change the DutyCycle to 0.25 and display the channel:

```
s.Channels.Frequency = 200;
```

s.Channels

ans =

Data acquisition counter output pulse generation channel 'ctr0' on de

```
IdleState: Low
InitialDelay: 2.5e-008
Frequency: 100
DutyCycle: 0.25
Terminal: 'PFI12'
Name: empty
ID: 'ctr0'
Device: [1x1 daq.ni.DeviceInfo]
MeasurementType: 'PulseGeneration'
```

DutyCycle

See Also Class

daq.Session, daq.Session.addCounterOutputChannel

Encoding type of counter channel

Description

When working with the session-based interface, use the EncoderType property to specify the encoding type of the counter input 'Position' channel.

Encoder types include:

- 'X1'
- 'X2'
- 'X4'
- 'TwoPulse'

See Also Class

daq.Session, daq.Session.addCounterInputChannel

EnhancedAliasRejectionEnable

Purpose Set enhanced alias rejection mode

Description Enable or disable the enhanced alias rejection on your DSA device's

analog channel. See "Synchronize DSA Devices" for more information. Enhanced alias reject is disabled by default. This property only takes

logical values. To enable it, type:

s.Channels(1).EnhancedAliasRejectionEnable = 1

You cannot modify enhanced rejection mode if you are synchronizing

your DSA device using ${\tt AutoSyncDSA}.$

See Also AutoSyncDSA

Store information for specific events

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Eventlog is a structure array that stores information related to specific analog input (AI) or analog output (AO) events. Event information is stored in the Type and Data fields of Eventlog. Type stores the event type. The logged event types are shown below.

Event Type	Description	Al	AO
Data missed	Data is missed by the engine.	✓	
Input overrange	A signal exceeds the hardware input range.	✓	
Run-time error	A run-time error is encountered. Run-time errors include timeouts and hardware errors.	✓	√
Start	The start function is issued.	✓	✓
Stop	The device object stops executing.	√	✓
Trigger	A trigger executes.	✓	✓

Timer events, samples available events (AI), and samples output events (AO) are not logged.

Data stores event-specific information associated with the event type in several fields. For all stored events, Data contains the RelSample field, which returns the input or output sample number at the time the event occurred. For the start, stop, run-time error, and trigger events, Data contains the AbsTime field, which returns the absolute time (as a clock vector) the event occurred. Other event-specific fields are included in Data. For a description of these fields, refer to "Events and Callbacks"

for analog input objects and analog output objects, or the appropriate reference pages in this chapter.

EventLog can store a maximum of 1000 events. If this value is exceeded, then the most recent 1000 events are stored. You can use the showdagevents function to easily display stored event information.

Characteristics

Usage AI, AO, common to all channels

Access Read-only

Data type Structure array

Read-only when N/A

running

Values

Values are automatically added as events occur. The default value is an empty structure array.

Examples

Create the analog input object ai and add four channels to it.

```
ai = analoginput('nidaq','Dev1');
chans = addchannel(ai,0:3);
```

Acquire 1 second of data and display the logged event types.

```
start(ai)
events = ai.EventLog;
{events.Type}
ans =
    'Start' 'Trigger' 'Stop
```

To examine the data associated with the trigger event:

```
events(2).Data
ans =
        AbsTime: [1999 2 12 14 54 52.5456]
    RelSample: 0
```

EventLog

Channel: []
Trigger: 1

See Also Functions

showdaqevents

ExcitationCurrent

Purpose Voltage of external source of excitation

Description When working with the session-based interface, the ExcitationCurrent

property indicates the voltage of the external source of excitation of an

accelerometer, IEPE, or microphone analog input channel.

See Also Properties

ExcitationSource

Class

daq.Session.addAnalogInputChannel

External source of excitation

Description

When working with the session-based interface, the ExcitationSource property indicates the source of excitation of an analog input channel, with accelerometer, IEPE, microphone, or bridge measurement type. Excitation source can be:

- Internal
- External
- None
- Unknown

By default, ExcitationSource is set to Unknown.

See Also

Properties

ExcitationCurrent

ExcitationVoltage

Class

daq.Session.addAnalogInputChannel

ExcitationVoltage

Purpose Voltage of excitation source

Description When working with RTD measurements in the session-based interface,

the ${\tt ExcitationVoltage}$ property indicates the excitation voltage value

to apply to the bridge measurement.

You can specify any positive voltage value. The default value is 0 volts.

See Also Properties

ExcitationSource

Class

ExternalTriggerTimeout

Purpose Indicate if external trigger timed out

Description When working with the session-based interface, the

ExternalTriggerTimeout property indicates if the if an

external trigger timed out.

See Also daq.Session.addTriggerConnection

Frequency

Purpose Frequency of generated pulses on counter output channel **Description** When working with the session-based interface, use the Frequency property to set the pulse repetition rate of a counter input channel. **Values** Specify the frequency in hertz. **Examples** Create a session object and add a 'PulseGeneration' counter output channel: s = daq.createSession('ni'); s.addCounterOutputChannel('cDAQ1Mod5', 'ctr0', 'PulseGeneration') Change the Frequency to 200 and display the channel: s.Channels.Frequency = 200; s.Channels ans = Data acquisition counter output pulse generation channel 'ctr0' on device IdleState: Low InitialDelay: 2.5e-008 Frequency: 200 DutyCycle: 0.5 Terminal: 'PFI12' Name: empty ID: 'ctr0' Device: [1x1 daq.ni.DeviceInfo] MeasurementType: 'PulseGeneration' See Also Class

daq.Session, daq.Session.addCounterInputChannel

Specify hardware channel ID

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

All channels contained by a device object have a hardware channel ID and an associated MATLAB index. The channel ID is given by HwChannel and the MATLAB index is given by the Index property. The HwChannel value is defined when hardware channels are added to a device object with the addchannel function.

The beginning channel ID value depends on the hardware device. For National Instruments hardware, channel IDs are zero-based (begin at zero). For sound cards, channel IDs are one-based (begin at one).

For scanning hardware, the scan order follows the MATLAB index. Therefore, the hardware channel associated with index 1 is sampled first, the hardware channel associated with index 2 is sampled second, and so on. To change the scan order, you can assign the channel IDs to different indices using HwChannel.

Characteristics

Usage AI, AO, per channel

Access Read/write
Data type Double
Read-only when Yes

running

Values

Values are automatically defined when channels are added to the device object with the addchannel function. The default value is one.

Examples

Create the analog input object ai for a National Instruments board and add the first three hardware channels to it.

HwChannel

```
ai = analoginput('nidaq','Dev1');
addchannel(ai,0:2);
```

Based on the current configuration, the hardware channels are scanned in order from 0 to 2. To swap the scan order of channels 0 and 1, you can assign these channels to the appropriate indices using HwChannel.

```
ai.Channel(1).HwChannel = 1;
ai.Channel(2).HwChannel = 0;
```

See Also Functions

addchannel

Properties

Channel, Index

Specify hardware line ID

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

All lines contained by a digital I/O object have a hardware ID and an associated MATLAB index. The hardware ID is given by HwLine and the MATLAB index is given by the Index property. The HwLine value is defined when hardware lines are added to a digital I/O object with the addline function.

The beginning line ID value depends on the hardware device. For National Instruments hardware, line IDs are zero-based (begin at zero).

Characteristics

Usage DIO, per line

Access Read/write
Data type Double
Read-only when running Yes

Values

Values are automatically defined when lines are added to the digital I/O object with the addline function. The default value is one.

Examples

Suppose you create the digital I/O object dio and add four hardware lines to it.

```
dio = digitalio('nidaq','Dev1');
addline(dio,0:3,'out');
```

addline automatically assigns the indices 1-4 to these hardware lines. You can swap the hardware lines associated with index 1 and index 2 with HwLine.

```
dio.Line(1).HwLine = 1;
```

HwLine

dio.Line(2).HwLine = 0;

See Also Functions

addline

Properties

Line, Index

ID of channel in session

Description

When working with the session-based interface, the ID property displays the ID of the channel. You set the channel ID when you add the channel to a session object.

Values

Examples

Create a session object, add a counter input channel, with the ID 'ctr0'.

```
s = daq.createSession('ni');
ch = s.addCounterInputChannel ('cDAQ1Mod5', 'ctr0', 'EdgeCount')
ch =
Data acquisition counter input edge count channel 'ctr0' on device 'cDAQ1Mod5':
     ActiveEdge: Rising
 CountDirection: Increment
   InitialCount: 0
       Terminal: 'PFI8'
          Name: empty
            ID: 'ctr0'
        Device: [1x1 daq.ni.DeviceInfo]
MeasurementType: 'EdgeCount'
Change CountDirection to 'Decrement':
ch.CountDirection = 'Decrement'
ch =
Data acquisition counter input edge count channel 'ctr0' on device 'cl
       ActiveEdge: Rising
```

CountDirection: Decrement

InitialCount: 0

Terminal: 'PFI8'
Name: empty
ID: 'ctr0'

Device: [1x1 daq.ni.DeviceInfo]

MeasurementType: 'EdgeCount'

See Also Class

Purpose Default state of counter output channel

Description When working with the session-based interface, the IdleState

property indicates the default state of the counter output channel with a 'PulseGeneration' measurement type when the counter is

not running.

Values IdleState is either 'High' or 'Low'.

Examples Create a session object and add a 'PulseGeneration' counter output channel:

```
s = daq.createSession('ni');
s.addCounterOutputChannel('cDAQ1Mod5', 'ctr0', 'PulseGeneration')
```

Change the IdleState property to 'High' and display the channel:

```
s.Channels.IdleState = 'High';
```

s.Channels

ans =

Data acquisition counter output pulse generation channel 'ctr0' on de

```
IdleState: High
InitialDelay: 2.5e-008
Frequency: 100
DutyCycle: 0.5
Terminal: 'PFI12'
Name: empty
ID: 'ctr0'
Device: [1x1 daq.ni.DeviceInfo]
MeasurementType: 'PulseGeneration'
```

IdleState

See Also Class

daq.Session, daq.Session.addCounterOutputChannel

MATLAB index of hardware channel or line

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Every hardware channel (line) contained by a device object has an associated MATLAB index that is used to reference that channel (line). For example, to configure property values for an individual channel, you must reference the channel through the Channel property using the appropriate Index value. Likewise, to configure property values for an individual line, you must reference the line through the Line property using the appropriate Index value.

For channels (lines), you can assign indices automatically with the addchannel (addline) function. Channel (line) indices always begin at 1 and increase monotonically up to the number of channels (lines) contained by the device object. For channels, index assignments can also be made manually with the addchannel function.

For scanning hardware, the scan order follows the MATLAB index. Therefore, the hardware channel associated with index 1 is sampled first, the hardware channel associated with index 2 is sampled second, and so on. To change the scan order, you can assign the channel IDs to different indices using the HwChannel or Channel property.

Index provides a convenient way to access channels and lines programmatically.

Characteristics

Usage AI, AO, per channel; DIO, per line

Access Read-only
Data type Double
Read-only when running N/A

Index

Values

Values are automatically defined when channels (lines) are added to the device object with the addchannel (addline) function. The default value is one.

Examples

Create the analog input object ai for a sound card and add two hardware channels to it.

```
ai = analoginput('winsound');
chans = addchannel(ai,1:2);
```

You can access the MATLAB indices for these channels with Index.

```
Index1 = chans(1).Index;
Index2 = chans(2).Index;
```

See Also

Functions

addchannel, addline

Properties

Channel, HwChannel, HwLine, Line

InitialDelay

Purpose Delay until output channel generates pulses

Description When working with the session-based interface, use the InitialDelay

property to set an initial delay on the counter output channel in which

the counter is running but does not generate any pulse.

See Also Class

Specify initial count point

Description

When working with the session-based interface, use the InitialCount property to set the point from which the device starts the counter.

Values

Examples

Create a session object, add counter input channel, and change the InitialCount.

```
s = daq.createSession('ni');
ch = s.addCounterInputChannel ('cDAQ1Mod5', 0, 'EdgeCount')
ch =
Data acquisition counter input edge count channel 'ctr0' on device 'cDAQ1Mod5':
     ActiveEdge: Rising
 CountDirection: Increment
   InitialCount: 0
       Terminal: 'PFI8'
          Name: empty
            ID: 'ctr0'
         Device: [1x1 dag.ni.DeviceInfo]
MeasurementType: 'EdgeCount'
Change InitalCount to 15:
 ch.InitialCount = 15
ch =
Data acquisition counter input edge count channel 'ctr0' on device 'cDAQ1
```

ActiveEdge: Rising CountDirection: Increment

InitialCount: 15

InitialCount

Terminal: 'PFI8'
Name: empty
ID: 'ctr0'

Device: [1x1 daq.ni.DeviceInfo]

MeasurementType: 'EdgeCount'

See Also Properties

Class

InitialTriggerTime

Purpose

Absolute time of first trigger

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

For all trigger types, InitialTriggerTime records the time when Logging or Sending is set to On. The absolute time is recorded as a clock vector.

You can return the InitialTriggerTime value with the getdata function, or with the Data.AbsTime field of the EventLog property.

If you synchronize multiple analoginput and analogoutput objects by setting TriggerType to HwDigitalTrigger and use the same digital trigger signal for all of the subsystems, the InitialTriggerTime property will not show the exact identical time for all subsystems.

Although the actual trigger events occurred simultaneously across all subsystems, the InitialTriggerTime events are recorded serially on a single thread. This causes the discrepancy of a few milliseconds. The time difference between InitialTriggerTime for multiple Data Acquisition Toolbox objects will not be consistent due to operating system process scheduling algorithms.

Characteristics

Usage AI, AO, common to all channels

Access Read-only

Data type Six-element vector of doubles

Read-only when N/A

running

Values

The value is automatically updated when the trigger executes. The default value is a vector of zeros.

Examples

Create the analog input object ai for a sound card and add two hardware channels to it.

```
ai = analoginput('winsound');
chans = addchannel(ai,1:2);
```

After starting ai, the trigger immediately executes and the trigger time is recorded.

```
start(ai)
abstime = ai.InitialTriggerTime
abstime =
1.0e+003 *
    1.9990    0.0020    0.0190    0.0130    0.0260    0.0208
```

To convert the clock vector to a more convenient form:

```
t = fix(abstime);
sprintf('%d:%d:%d', t(4),t(5),t(6))
ans =
13:26:20
```

See Also

Functions

getdata

Properties

EventLog, Logging, Sending

InputOverRangeFcn

Purpose

Specify callback function to execute when acquired data exceeds valid hardware range

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

An input overrange event is generated immediately after an overrange condition is detected for any channel group member. This event executes the callback function specified for InputOverRangeFcn.

An overrange condition occurs when an input signal exceeds the range specified by the SensorRange property. Overrange detection is enabled only if the analog input object is running and a callback function is specified for InputOverRangeFon.

Input overrange event information is stored in the Type and Data fields of the EventLog property. The Type field value is OverRange. The Data field values are given below.

Note The input overrange event is not generated if a signal begins outside the range and then goes into the range.

Data Field Value	Description
AbsTime	The absolute time (as a clock vector) the event occurred.
Channel	The index of the channel that experienced an overrange signal.

Data Field Value	Description
OverRange	The OverRange value, Off indicates that the channel went from overrange to in range, and On indicates that it went from in range to overrange.
RelSample	The acquired sample immediately before the moment when the overrange transition occurs.

Note The input signal values will not exceed the values set by the InputRange property. If you set InputRange and SensorRange to the same value, the OverRange event is never received. To receive OverRange events set the value of SensorRange within, and not equal to, the InputRange value.

Characteristics

Usage AI, common to all channels

Access Read/write

Data type String

Read-only when No
running

Values

The default value is an empty string.

See Also Properties

EventLog, SensorRange

InputRange

Purpose

Specify range of analog input subsystem

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

InputRange is a two-element vector that specifies the range of voltages that can be accepted by the analog input (AI) subsystem. You should configure InputRange so that the maximum dynamic range of your hardware is utilized.

If an input signal exceeds the InputRange value, then an overrange condition occurs. Overrange detection is enabled only if the analog input object is running and a value is specified for the InputOverRangeFcn property. For many devices, the input range is expressed in terms of the gain and polarity.

AI subsystems have a finite number of InputRange values that you can set. If an input range is specified but does not match a valid range, then the next highest supported range is automatically selected by the engine. If InputRange exceeds the range of valid values, then an error is returned. Use the daqhwinfo function to return the input ranges supported by your board.

Because the engine can set the input range to a value that differs from the value you specify, you should return the actual input range for each channel using the get function or the device object display summary. Alternatively, you can use the setverify function, which sets the InputRange value and then returns the actual value that is set.

Note If your hardware supports a channel gain list, then you can configure InputRange for individual channels. Otherwise, InputRange must have the same value for all channels contained by the analog input object.

You should use InputRange in conjunction with the SensorRange property. These two properties should be configured such that the maximum precision is obtained and the full dynamic range of the sensor signal is covered.

Characteristics

Usage AI, per channel

Access Read/write

Data type Two-element vector of doubles

Read-only when Yes

running

Values

The default value is supplied by the hardware driver.

Examples

Create the analog input object ai for a National Instruments board, and add two hardware channels to it.

```
ai = analoginput('nidaq','Dev1');
addchannel(ai,0:1);
```

You can return the input ranges supported by the board with the InputRanges field of the daqhwinfo function.

```
out = daqhwinfo(ai);
out.InputRanges
ans =
    -0.0500     0.0500
    -0.5000     0.5000
    -5.0000     5.0000
    -10.0000     10.0000
```

To configure both channels contained by ai to accept input signals between -10 volts and 10 volts:

```
ai.Channel.InputRange = [-10 10];
```

InputRange

Some devices allow you to set each channel's InputRange property independently:

```
ai.Channel(1).InputRange = [-0.05     0.05];
ai.Channel(2).InputRange = [-10 10];
Alternatively, you can use the setverify function.
```

ActualRange = setverify(ai.Channel, 'InputRange', [-10 10]);

See Also Functions

daqhwinfo, setverify

Properties

InputOverRangeFcn, SensorRange, Units, UnitsRange

Specify analog input hardware channel configuration

Description

For National Instruments devices, InputType can be SingleEnded, Differential, NonReferencedSingleEnded, or PseudoDifferential. For Measurement Computing devices, InputType can be SingleEnded, or Differential. For sound cards, InputType can only be AC-Coupled.

If channels have been added to a National Instruments or Measurement Computing analog input object and you change the InputType value, then the channels are automatically deleted if the hardware reduces the number of available channels.

Characteristics

Usage AI, common to all channels

Access Read/write

Data type String
Read-only when Yes

running

Values

Advantech and Measurement Computing

Differential Channels are configured for differential input.

SingleEnded Channels are configured for single-ended

input.

The value for InputType on Advantech and MCC boards is always read-only in MATLAB. For Advantech boards, the setting is made in the Advantech Device Manager. For Measurement Computing boards, the setting is made in InstaCal.

InputType

National Instruments

{Differential} Channels are configured for

differential input.

SingleEnded Channels are configured for

single-ended input.

NonReferencedSingleEnded This channel configuration is

used when the input signal has its own ground reference, which is tied to the negative input of the instrumentation amplifier.

PseudoDifferential Channels are configured for

pseudodifferential input, which are all referred to a common ground but this ground is not connected to the computer

ground.

Sound Cards

{AC-Coupled} The input is coupled so that constant (DC)

signal levels are suppressed.

Specify if operation continues until manually stopped

Description

When working with the session-based interface, use IsContinuous to specify that the session operation runs until you execute daq.Session.stop. When set to true, the session will run continuously, acquiring or generating data until stopped.

Values

{false}

Set the IsContinuous property to false to make the session operation stop automatically. This property is set to false by default.

true

Set the IsContinuous property to true to make the session operation run until you execute daq. Session.stop.

Examples

Create a session object, add an analog input channel, and set the session to run until manually stopped:

See Also Properties

IsDone

IsContinuous

Methods

 ${\tt daq.Session.stop}, {\tt daq.Session.startBackground}$

Class

Indicate if operation is complete

Description

When working with the session-based interface, the read-only IsDone property indicates if the session operation is complete.

Tip IsDone indicates if the session object has completed acquiring or generating data. IsRunning indicates if the operation is in progress, but the hardware may not be acquiring or generating data. IsLogging indicates that the hardware is acquiring or generating data.

Values

true

Value is true if the operation is complete.

false

Value is false if the operation is not complete.

Examples

Create an acquisition session and see if the operation is complete:

```
s = daq.createSession('ni');
s.addAnalogOutputChannel('cDAQ1Mod2', 'ao1', 'vVoltage');
s.queueOutputData (linspace(-1, 1, 1000)');
s.startBackground();
s.IsDone
ans =
   0
```

Issue a wait and see if the operation is complete:

```
s.wait()
s.IsDone
ans =
```

1

See Also Methods

daq.Session.startBackground

Class

Indicate if hardware is acquiring or generating data

Description

When working with the session-based interface, the status of the read-only IsLogging property indicates if the hardware is acquiring or generating data.

Tip IsLogging indicates that the hardware is acquiring or generating data. IsRunning indicates if the operation is in progress, but the hardware may not be acquiring or generating data. IsDone indicates if the session object has completed acquiring or generating data.

Values

true

Value is true if the device is acquiring or generating data.

false

Value is false if the device is not acquiring or generating data.

Examples

Create a session and see if the operation is logging:

```
s = daq.createSession('ni');
s.addAnalogOutputChannel('cDAQ1Mod2', 'ao1', 'Voltage');
s.queueOutputData (linspace(-1, 1, 1000)');
s.startBackground();
s.IsLogging
ans =
1
```

Wait until the operation is complete:

```
s.wait()
s.IsLogging
ans =
```

IsLogging

0

See Also Properties

IsRunning, IsDone

Methods

daq.Session.startBackground

Class

IsNotifyWhenDataAvailableExceedsAuto

Purpose

Control if NotifyWhenDataAvailableExceeds is set automatically

Description

When working with the session-based interface, the IsNotifyWhenDataAvailableExceedsAuto property indicates if the NotifyWhenDataAvailableExceeds property is set automatically, or you have set a specific value.

Tip This property is typically used to set NotifyWhenDataAvailableExceeds back to its default behavior.

Values

{true}

When the value is true, then the NotifyWhenDataAvailableExceeds property is set automatically.

false

When the value is false, when you have set the NotifyWhenDataAvailableExceeds property to a specific value.

See Also

Properties

NotifyWhenDataAvailableExceeds

Events

DataAvailable

Class

IsNotifyWhenScansQueuedBelowAuto

Purpose Control if NotifyWhenScansQueuedBelow is set automatically

Description When working with the session-based interface, the

IsNotifyWhenScansQueuedBelowAuto property indicates if

the NotifyWhenScansQueuedBelow property is set automatically, or you

have set a specific value.

Values {true}

When the value is true, then NotifyWhenScansQueuedBelow is

set automatically.

false

When the value is false, you have set

NotifyWhenScansQueuedBelow property to a specific value.

See Also Properties

NotifyWhenScansQueuedBelow, ScansQueued

Events

DataRequired

Class

Indicate if operation is still in progress

Description

When working with the session-based interface, the IsRunning status indicates if the operation is still in progress.

Tip IsRunning indicates if the operation is in progress, but the hardware may not be acquiring or generating data. IsLogging indicates that the hardware is acquiring or generating data. IsDone indicate is if the session object has completed acquiring or generating.

Values

true

When the value is true if the operation is in progress.

false

When the value is false if the operation is not in progress.

Examples

Create an acquisition session, add a DataAvailable event listener and start the acquisition.

IsRunning

Wait until operation completes and see if session is in progress:

```
s.wait()
s.IsRunning
ans =
0
```

See Also

Properties

IsLogging, IsDone

Methods

daq.Session.startBackground

Class

daq.Session

Indicate if device is simulated

Description

When working with the session-based interface, the IsSimulated property indicates if the session is using a simulated device.

Values

true

When the value is true if the operation is in progress.

false

When the value is false if the operation is not in progress.

Examples

Discover available devices.

>> d = daq.getDevices

d =

Data acquisition devices:

index	Vendor	Device ID		Description	
1	ni	cDAQ1Mod1	National	Instruments	NI 9201
2	ni	cDAQ2Mod1	National	Instruments	NI 9201
3	ni	Dev1	National	Instruments	USB-6211
4	ni	Dev2	National	Instruments	USB-6218
5	ni	Dev3	National	Instruments	USB-6255
6	ni	Dev4	National	Instruments	USB-6363
7	ni	PXI1Slot2	National	Instruments	PXI-4461
8	ni	PXI1Slot3	National	Instruments	PXI-4461

Examine properties of NI 9201, with the device id cDAQ1Mod1 with the index 1.

>> d(1)

ans =

```
ni: National Instruments NI 9201 (Device ID: 'cDAQ1Mod1')
   Analog input subsystem supports:
      -10 to +10 Volts range
      Rates from 0.1 to 800000.0 scans/sec
      8 channels ('ai0', 'ai1', 'ai2', 'ai3', 'ai4', 'ai5', 'ai6', 'ai7')
      'Voltage' measurement type
This module is in slot 4 of the 'cDAQ-9178' chassis with the name 'cDAQ1'
Properties, Methods, Events
Click the Properties link to see the properties of the device.
  ChassisName: 'cDAQ1'
    ChassisModel: 'cDAQ-9178'
      SlotNumber: 4
     IsSimulated: true
       Terminals: [48x1 cell]
          Vendor: National Instruments
               ID: 'cDAQ1Mod1'
           Model: 'NI 9201'
      Subsystems: [1x1 daq.ni.AnalogInputInfo]
     Description: 'National Instruments NI 9201'
RecognizedDevice: true
Note that the IsSimulated value is true, indicating that this device is
simulated.
Properties
IsLogging, IsDone
Methods
dag.Session.startBackground
Class
daq.Session
```

See Also

IsWaitingForExternalTrigger

Purpose Indicates if synchronization is waiting for an external trigger

Description When working with the session-based interface, the

read-onlyIsWaitingForExternalTrigger property indicates

if the acquisition or generation session is waiting for a trigger from an external device. If you have added an external trigger, this property

displays true, if not, it displays false.

See Also daq.Session.addTriggerConnection

Contain hardware lines added to device object

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Line is a vector of all the hardware lines contained by a digital I/O (DIO) object. Because a newly created DIO object does not contain hardware lines, Line is initially an empty vector. The size of Line increases as lines are added with the addline function, and decreases as lines are removed with the delete function.

You can use Line to reference one or more individual lines. To reference a line, you must know its MATLAB index and hardware ID. The MATLAB index is given by the Index property, while the hardware ID is given by the HwLine property.

Characteristics

Usage DIO

Access Read/write

Data type Vector of lines

Read-only when

running

Values

Values are automatically defined when lines are added to the DIO object with the addline function. The default value is an empty column vector.

Examples

Create the digital I/O object dio and add four input lines to it.

Yes

```
dio = digitalio('nidaq','Dev1');
addline(dio,0:3,'In');
```

To set a property value for the first line added (ID = 0), you can reference the line by its index using the Line property.

```
line1 = dio.Line(1);
set(line1, 'Direction', 'Out')
```

See Also Functions

addline, delete

Properties

HwLine, Index

Specify descriptive line name

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

LineName specifies a descriptive name for a hardware line. If a line name is defined, then you can reference that line by its name. If a line name is not defined, then the line must be referenced by its index. Line names are not required to be unique.

You can also define descriptive line names when lines are added to a digital I/O object with the addline function.

Characteristics

Usage DIO, per line

Access Read/write

Data type String

Read-only when Yes

running

Values

The default value is an empty string. To reference a line by name, it must contain only letters, numbers, and underscores and must begin with a letter.

Examples

Create the digital I/O object dio and add four hardware lines to it.

```
dio = digitalio('nidaq','Dev1');
addline(dio,0:3,'out');
```

To assign a descriptive name to the first line contained by dio:

```
line1 = dio.Line(1);
set(line1, 'LineName', 'Joe')
```

You can now reference this line by name instead of index.

set(dio.Joe, 'Direction', 'In')

See Also Functions

addline

LogFileName

Purpose

Specify name of disk file information is logged to

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

You can log acquired data, device object property values and event information, and hardware information to a disk file by setting the LoggingMode property to Disk or Disk&Memory.

You can specify any value for LogFileName as long as it conforms to the MATLAB software naming conventions: the name cannot start with a number and cannot contain spaces. If no extension is specified as part of LogFileName, then daq is used. The default value for LogFileName is logfile.daq.

You can choose whether an output file is overwritten or if multiple log files are created with the LogToDiskMode property. Setting LogToDiskMode to Overwrite causes the output file to be overwritten. Setting LogToDiskMode to Index causes new data files to be created, each with an indexed name based on the value of LogFileName.

Characteristics

Usage AI, common to all channels

Yes

Access Read/write

Data type String

Read-only when

running

Values

The default value is logfile.daq.

See Also Properties

Logging, LoggingMode, LogToDiskMode

Indicate whether data is being logged to memory or disk file

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Along with the Running property, Logging reflects the state of an analog input object. Logging can be On or Off.

Logging is automatically set to On when a trigger occurs. When Logging is On, acquired data is being stored in memory or to a disk file.

Logging is automatically set to Off when the requested samples are acquired, an error occurs, or a stop function is issued. When Logging is Off, you can still preview data with the peekdata function provided Running is On. However, peekdata does not guarantee that all the requested data is returned.

To guarantee that acquired data contains no gaps, is must be logged to memory or to a disk file. Data stored in memory is extracted with the getdata function, while data stored to disk is returned with the daqread function. The destination for logged data is controlled with the LoggingMode property.

Characteristics

Usage AI, common to all channels

Access Read-only
Data type String
Read-only when N/A

running

Values

{Off} Data is not logged to memory or a disk file.

On Data is logged to memory or a disk file.

Logging

See Also Functions

daqread, getdata, peekdata, stop

Properties

LoggingMode, Running

Specify destination for acquired data

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

You can set LoggingMode to Disk, Memory, or Disk&Memory. When you set LoggingMode to Disk, then acquired data (as well as device object and hardware information) is streamed to a disk file. If LoggingMode is set to Memory, then acquired data is stored in the data acquisition engine. If LoggingMode is set to Disk&Memory, then acquired data is stored in the data acquisition engine and is streamed to a disk file.

When logging to the engine, you must extract the data with the getdata function. If you do not extract this data, and the amount of data stored in memory reaches the limit for the data acquisition object (see dagmem(obj)), a **DataMissed** event occurs. At this point, the acquisition stops.

When logging to disk, you can specify the log filename with the LogFileName property, and you can control the number of log files created with the LogToDiskMode property. You can return data stored in a disk file to the MATLAB workspace with the dagread function.

Characteristics

Usage AI, common to all channels

Read/write Access

Data type String Yes

Read-only when

running

LoggingMode

Values Disk Acquired data is logged to a disk file.

{Memory} Acquired data is logged to memory.

Disk&Memory Acquired data is logged to a disk file and to

memory.

See Also Functions

daqread, getdata

Properties

LogFileName, LogToDiskMode

Specify whether data, events, and hardware information are saved to one or more disk files

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

LogToDiskMode can be set to Overwrite or Index. If LogToDiskMode is set to Overwrite, then the log file is overwritten each time start is issued. If LogToDiskMode is set to Index, a different disk file is created each time start is issued and these rules are followed:

- The first log filename is specified by the initial value of LogFileName.
- If the specified file already exists, it is overwritten and no warning is issued.
- LogFileName is automatically updated with a numeric identifier after each file is written. For example, if LogFileName is initially specified as data.daq, then data.daq is the first filename, dataO1.daq is the second filename, and so on.

Separate analog input objects are logged to separate files. You can return data stored in a disk file to the MATLAB workspace with the daqread function. If an error occurs during data logging, an error message is returned and data logging is stopped.

Characteristics

Usage AI, common to all channels

Access Read/write

Data type String
Read-only when Yes

running

${\bf LogToDiskMode}$

Values Index Multiple log files are written, each with an

indexed filename based on the LogFileName

property.

{Overwrite} The log file is overwritten.

See Also Functions

dagread

Properties

LogFileName, LoggingMode

ManualTriggerHwOn

Purpose

Specify hardware device starts at manual trigger

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

You can set ManualTriggerHwOn to Start or Trigger, and it has an effect only when the TriggerType property value is Manual. If ManualTriggerHwOn is Start, then the hardware device associated with your device object starts running after you issue the start function. If ManualTriggerHwOn is Trigger, then the hardware device associated with your device object starts acquiring after you issue both the start function and you execute a manual trigger with the trigger function. You can use trigger only when you configure the TriggerType property to Manual.

You should configure ManualTriggerHwOn to Trigger when you want to synchronize the input and output of data, or you require more control over when your hardware starts. Note that you cannot use peekdata or acquire pretrigger data when you use this value. Additionally, you should not use this value with repeated triggers because the subsequent behavior is undefined.

Characteristics

Usage AI, common to all channels

Access Read/write

Data type String

Read-only when Yes

running

ManualTriggerHwOn

Values

{Start} Start the hardware after the start function is issued.

Trigger Start the hardware after the trigger function is

issued.

Examples

Create the analog input object ai and the analog output object ao for a sound card and add two channels to each device object.

```
ai = analoginput('winsound');
addchannel(ai,1:2);
ao = analogoutput('winsound');
addchannel(ao,1:2);
```

To operate the sound card in full duplex mode, and to minimize the time between when ai starts and ao starts, you configure ManualTriggerHwOn to Trigger for ai and TriggerType to Manual for both ai and ao.

```
set([ai ao],'TriggerType','Manual')
ai.ManualTriggerHwOn = 'Trigger';
```

The analog input and analog output hardware devices will both start after you issue the trigger function. For a detailed example that uses ManualTriggerHwOn, refer to "Start Multiple Device Objects".

See Also

Functions

peekdata, start, trigger

Properties

TriggerType

Indicate maximum number of samples that can be queued in engine

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

MaxSamplesQueued indicates the maximum number of samples allowed in the analog output queue.

If the BufferingMode is set to Auto, the default value is calculated by the engine, and is based on the memory resources of your system. You can override the default value of MaxSamplesQueued with the dagmem function.

If the BufferingMode is set to Manual, MaxSamplesQueued is updated to indicate the maximum number of samples allowed in the analog output queue based on the number of buffers selected in BufferingConfig.

The value of MaxSamplesQueued can affect the behavior of putdata. For example, if the queued data exceeds the value of MaxSamplesQueued, then putdata becomes a blocking function until there is enough space in the queue to add the additional data.

Characteristics

Usage AO, common to all channels

Access Read-only
Data type Double
Read-only when N/A

running

Values

The value is calculated by the data acquisition engine.

See Also Functions

daqmem, putdata

MaxSoundPressureLevel

Purpose Sound pressure level for microphone channels

Description When working with the session-based interface, use the

MaxSoundPressureLevel set the maximum sound pressure of the

microphone channel in decibels.

Values The maximum sound pressure level is based on the sensitivity and the

voltage range of your device. When you sent your device Sensitivity, the MaxSoundPressure value is automatically corrected to match the specified sensitivity value and the device voltage range. You can also specify any acceptable pressure level in decibels. Refer to your

microphone specifications for more information.

Type counter channel measurement

Description

When working with the session-based interface, the MeasurementType property displays the selected measurement type for your channel.

Values

Counter measurement types include:

- 'EdgeCount' (input)
- 'PulseWidth' (input)
- 'Frequency'(input)
- 'Position'(input)
- 'PulseGeneration' (output)

Analog measurement types include:

- 'Voltage' (input and output)
- 'Thermocouple' (input)
- 'Current' (input and output)
- 'Accelerometer' (input)
- 'RTD' (input)
- 'Bridge' (input)
- 'Microphone' (input)
- 'IEPE' (input)

Examples

Create a session object, add a counter input channel, with the 'EdgeCount' MeasurementType.

```
s = daq.createSession('ni');
ch = s.addCounterInputChannel ('cDAQ1Mod5', 0, 'EdgeCount')
ch =
```

MeasurementType

```
Data acquisition counter input edge count channel 'ctr0' on device 'cDAQ1Mod5':

ActiveEdge: Rising
CountDirection: Increment
InitialCount: 0
Terminal: 'PF18'
Name: empty
ID: 'ctr0'
Device: [1x1 daq.ni.DeviceInfo]
MeasurementType: 'EdgeCount'
```

See Also

daq.Session.addAnalogInputChannel, daq.Session.addAnalogOutputChannel, daq.Session.addCounterInputChannel, daq.Session.addCounterOutputChannel

Class

daq.Session

Specify descriptive name for the channel

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

When you add a channel, a descriptive name is automatically generated and stored in Name. The name is a concatenation the name of the adaptor, the device ID, and the device object type. You can change the value of Name at any time.

Values

The value is defined when you add the channel.

Examples

Create the analog input object ai for a sound card.

```
ai = analoginput('winsound');
```

The descriptive name for ai is given by

ai.Name

```
ans =
winsound0-AI
```

Change the name to WindowsSoundChannel and access the name

```
ai.Name='WindowsSoundChannel'
```

Specify descriptive name for the channel

Description

When you add a channel, a descriptive name is stored in Name. By default there is no name assigned to the channel. You can change the value of Name at any time.

Values

You can specify a string value for the name.

Examples

Change the name of an analog input channel

Create a session and add an analog input channel.

Create a session object, add a counter input channel, and change the Name.

```
s = daq.createSession('ni');
ch = s.addCounterInputChannel ('cDAQ1Mod5', 0, 'EdgeCount')

ch =

Data acquisition counter input edge count channel 'ctr0' on device 'cDAQ1Mod5':

    ActiveEdge: Rising
    CountDirection: Increment
    InitialCount: 0
```

```
Name: empty
ID: 'ctr0'
Device: [1x1 daq.ni.DeviceInfo]

MeasurementType: 'EdgeCount'

Change Name to 'AI-Voltage':
s.Channels(1).Name='AI-Voltage'

s =

Data acquisition session using National Instruments hardware:
Will run for 1 second (1000 scans) at 1000 scans/second.
Number of channels: 1
index Type Device Channel MeasurementType Range National Instruments hardware:

1 ai Dev1 ai0 Voltage (Diff) -10 to +10 Volts AI-Voltage
```

See Also Class

daq.Session

Terminal: 'PFI8'

NativeOffset

Purpose

Indicate offset to use when converting between native data format and doubles

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

NativeOffset, along with NativeScaling, is used to convert data between the native hardware format and doubles.

For analog input objects, you return native data from the engine with the getdata function. Additionally, if you log native data to a .daq file, then you can read back that data using the daqread function. The formula for converting from native data to doubles is

doubles data = (native data)(native scaling) + native offset

For analog output objects, you queue native data in the engine with the putdata function. The formula for converting from doubles to native data is

native data = (doubles data)(native scaling) + native offset

You return the native data type of your hardware device with the daqhwinfo function. Note that the NativeScaling value for a given channel might change if you change its InputRange (AI) or OutputRange (AO) property value.

You might want to return or queue data in native format to conserve memory and to increase data acquisition or data output speed.

Characteristics

Usage AI, AO, per channel

Access Read-only

Data type Double
Read-only when N/A
running

Values

The default value is device-specific.

Examples

Create the analog input object ai for a National Instruments board, and add eight channels to it.

```
ai = analoginput('nidaq','Dev1');
addchannel(ai,0:7);
```

Start ai, collect one second of data for each channel, and extract the data from the engine using the native format of the device.

```
start(ai)
nativedata = getdata(ai,1000, 'native');
```

You can return the native data type of the board with the daqhwinfo function.

```
out = daqhwinfo(ai);
out.NativeDataType
ans =
double
```

Convert the data to doubles using the NativeScaling and NativeOffset properties.

```
scaling = get(ai.Channel(1),'NativeScaling');
offset = get(ai.Channel(1),'NativeOffset');
data = double(nativedata)*scaling + offset;
```

See Also Functions

daqhwinfo, daqread, getdata, putdata

NativeOffset

Properties

InputRange, NativeScaling, OutputRange

Indicate scaling to use when converting between native data format and doubles

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

NativeScaling, along with NativeOffset, is used to convert data between the native hardware format and doubles.

For analog input objects, you return native data from the engine with the getdata function. Additionally, if you log native data to a .daq file, then you can read back that data using the daqread function. The formula for converting from native data to doubles is

doubles data = (native data)(native scaling) + native offset

For analog output objects, you queue native data in the engine with the putdata function. The formula for converting from doubles to native data is

native data = (doubles data)(native scaling) + native offset

You return the native data type of your hardware device with the daqhwinfo function. Note that the NativeScaling value for a given channel might change if you change its InputRange (AI) or OutputRange (AO) property value.

You might want to return or queue data in native format to conserve memory and to increase data acquisition or data output speed.

Characteristics

Usage AI, AO, per channel

Access Read-only

NativeScaling

Data type Double

Read-only when N/A

running

Values The default value is device-specific.

See Also Functions

daqhwinfo, daqread, getdata, putdata

Properties

InputRange, NativeOffset, OutputRange

NominalBridgeResistance

Purpose Resistance of sensor

Description When working with the session-based interface, the

NominalBridgeResistance property displays the resistance of a bridge—based sensor in ohms. This value is used to calculate

voltage.

You can specify any accepted positive value in ohms. The default value is 0 until you change it. You must set the resistance to use the channel.

See Also Class

daq.Session

NotifyWhenDataAvailableExceeds

Purpose

Control firing of DataAvailable event

Description

When working with the session-based interface the DataAvailable event is fired when the scans available to the session object exceeds the value specified in the NotifyWhenDataAvailableExceeds property.

Values

By default the DataAvailable event fires when 1/10 second worth of data is available for analysis. To specify a different threshold change this property to control when DataAvailable fires.

Examples

Create the session and add an analog input voltage channel. Add an event listener to display the total number of scans acquired for this operation:

```
s = daq.createSession('ni');
s.addAnalogInputChannel('cDAQ1Mod4', 1, 'Voltage');
lh = s.addlistener('DataAvailable', ...
@(src, event) disp(s.ScansAcquired));
```

The default the Rate is 1000 scans per second. The session is automatically configured to fire the DataAvailable notification 10 times per second.

Increase the Rate to 800,000 scans per second and the DataAvailable notification automatically fires 10 times per second:

```
s.Rate = 800000;
```

NotifyWhenDataAvailableExceeds

Running the acquisition causes the number of scans acquired to be displayed by the callback 10 times:

```
data = s.startForeground;
80000
160000
240000
320000
400000
480000
560000
640000
720000
800000
```

Increase NotifyWhenDataAvailableExceeds to 160,000. NotifyWhenDataAvailableExceeds no longer configured automatically when the Rate changes.

The DataAvailable event is fired only five times per second.

NotifyWhenDataAvailableExceeds

Set IsNotifyWhenDataAvailableExceedsAuto back to true. This causes NotifyWhenDataAvailableExceeds to set automatically when Rate changes.

See Also Properties

Is Notify When Data Available Exceeds Auto

Events

DataAvailable

Class

daq.Session

NotifyWhenScansQueuedBelow

Purpose Control firing of DataRequired event

Description When working with the session-based interface to generate output

signals continuously, the DataRequired event is fired when you need to queue more data. This occurs when the ScansQueued property drops below the value specified in the NotifyWhenScansQueuedBelow

property.

Values By default the DataRequired event fires when 1/2 second worth of data

remains in the queue. To specify a different threshold, change the this

property to control when DataRequired is fired.

See Also Properties

ScansQueued, IsNotifyWhenScansQueuedBelowAuto

Events

DataRequired

Class

daq.Session

Number of scans for operation when starting

Description

When working with the session-based interface, use the NumberOfScans property to specify the number of scans the session will acquire during the operation. Changing the number of scans changes the duration of an acquisition. When the session contains output channels, NumberOfScans becomes a read only property and the number of scans in a session is determined by the amount of data queued.

Tips

- To specify length of the acquisition, use DurationInSeconds.
- To control length of the output operation, use daq.Session.queueOutputData.

Values

You can change the value only when you use input channels.

Examples

Create an acquisition session, add an analog input channel, and display the NumberOfScans.

Change the NumberOfScans property.

```
s.NumberOfScans = 2000
s =
```

NumberOfScans

```
Data acquisition session using National Instruments hardware:

Will run for 2000 scans (2 seconds) at 1000 scans/second.

Operation starts immediately.

Number of channels: 1

index Type Device Channel InputType Range Name

1 ai cDAQ1Mod1 ai0 Diff -10 to +10 Volts
```

Changing NumberOfScans changes the duration.

See Also Properties

ScansQueued, DurationInSeconds

Methods

daq.Session.startForeground, daq.Session.startBackground,
daq.Session.queueOutputData

Class

daq.Session

OutputRange

Purpose

Specify range of analog output hardware subsystem

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

OutputRange is a two-element vector that specifies the range of voltages that can be output by the analog output (AO) subsystem. You should configure OutputRange so that the maximum dynamic range of your hardware is utilized. For many devices, the output range is expressed in terms of the gain and polarity.

AO subsystems have a finite number of OutputRange values that you can set. If an output range is specified but does not match a valid range, then the next highest supported range is automatically selected by the engine. If OutputRange exceeds the range of valid values, then an error is returned. Use the daqhwinfo function to return the output ranges supported by your board.

Because the engine can set the output range to a value that differs from the value you specify, you should return the actual output range for each channel using the get function or the device object display summary. Alternatively, you can use the setverify function, which sets the OutputRange value and then returns the actual value that is set.

Characteristics

Usage AO, per channel

Access Read/write

Data type Two-element vector of doubles

Read-only when Yes

running

Values

The default value is determined by the hardware driver.

Examples

Create the analog output object ao for a National Instruments board and add two hardware channels to it.

```
ao = analogoutput('nidaq','Dev1');
addchannel(ao,0:1);
```

You can return the output ranges supported by the board with the OutputRanges field of the daghwinfo function.

```
out = daqhwinfo(ao);
out.OutputRanges
ans =
    0.0000    10.0000
    -10.0000    10.0000
```

To configure both channels contained by **ao** to output signals between -10 volts and 10 volts:

```
ao.Channel.OutputRange = [-10 10];
```

Alternatively, you can use the setverify function to configure and return the OutputRange value.

ActualRange = setverify(ao.Channel, 'OutputRange',[-10 10]);

See Also Functions

daqhwinfo, setverify

Properties

Units, UnitsRange

Indicate parent (device object) of channel or line

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

The parent of a channel (line) is defined as the device object that contains the channel (line).

You can create a copy of the device object containing a particular channel or line by returning the value of Parent. You can treat this copy like any other device object. For example, you can configure property values, add channels or lines to it, and so on.

Characteristics

Usage AI, AO, per channel; DIO, per line

Access Read-only

Data type Device object

Read-only when N/A

running

Values

The value is defined when channels or lines are added to the device object.

Examples

Create the analog input object ai for a National Instruments board and add three hardware channels to it.

```
ai = analoginput('nidaq','Dev1');
chans = addchannel(ai,0:2);
```

To return the parent for channel 2:

parent = ai.Channel(2).Parent;

parent is an exact copy of ai.

```
isequal(ai,parent)
ans =
    1
```

Specify port ID

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Hardware lines are often grouped together as a port. Digital I/O subsystems can consist of multiple ports and typically have eight lines per port. When adding hardware lines to a digital I/O object with addline, you can specify the port ID. The port ID is stored in the Port property. If the port ID is not specified, then the smallest port ID value is automatically used.

Characteristics

Usage DIO, per line

Access Read-only
Data type Double
Read-only when N/A

running

Values

The port ID is defined when line are added to the digital I/O object with addline.

Examples

Create the digital I/O object dio and add two hardware channels to it.

```
dio = digitalio('nidaq','Dev1');
addline(dio,0:1,'In');
```

You can use Port property to return the port IDs associated with the lines contained by dio.

```
dio.Line.Port
ans =
  [0]
```

[0]

See Also Functions

addline

Purpose Specify resistance value

Description Use this property to specify the resistance of the device.

You can specify any acceptable value in ohms. When you add an RTD Channel, the resistance is unknown and the RO property displays Unknown. You must change this value to set the resistance of this device.

See Also Class

daq.Session

Properties

RTDConfiguration, RTDType

Purpose Specify channel measurement range

Description When working with the session-based interface, use the Range to

indicate the measurement range of a channel.

Values Range is not applicable for counter channels. For analog channels,

value is dependent on the measurement type. This property is read-only for all measurement types except 'Voltage'. You can specify a range in

volts for analog channels.

Rate of operation in scans per second

Description

When working with the session-based interface, use the Rate property to set the number of scans per second.

Note Many hardware devices accept fractional rates.

Tip On most devices, the hardware limits the exact rates that you can set. When you set the rate, Data Acquisition Toolbox sets the rate to the next higher rate supported by the hardware. If the exact rate affects your analysis of the acquired data, obtain the actual rate after you set it, and then use that in your analysis.

Values

You can set the rate to any positive nonzero scalar value supported by the hardware in its current configuration.

Examples

Create a session object and add an analog input channel and change the session rate.

```
s = daq.createSession('ni');
s.addAnalogInputChannel('cDAQ1Mod1','ai1','Voltage');
s.Rate = 10000
s =
```

Data acquisition session using National Instruments hardware: Will run for 1 second (10000 scans) at 10000 scans/second. Operation starts immediately.

```
Number of channels: 1

index Type Device Channel InputType Range Name

1 ai cDAQ1Mod1 ai1 Diff -10 to +10 Volts
```

Properties, Methods, Events

See Also Properties

DurationInSeconds, NumberOfScans, RateLimit

Class

daq.Session

RateLimit

Purpose

Limit of rate of operation based on hardware configuration

Description

In the session-based interface, the read-only RateLimit property displays the minimum and maximum rates that the session supports, based on the device configuration for the session.

Tip RateLimit changes dynamically as the session configuration changes.

See Also

Properties

Rate

Class

daq.Session

Specify number of additional times queued data is output

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

To send data to an analog output subsystem, it must first be queued in the data acquisition engine with the putdata function. If you want to continuously output the same data, you can use multiple calls to putdata. However, because each putdata call consumes memory, a long output sequence can quickly bring your system to halt.

As an alternative to putdata, you can continuously output previously queued data using RepeatOutput. Because RepeatOutput requeues the data, additional memory resources are not consumed. While the data is being output, you cannot add additional data to the queue.

Characteristics

Usage AO, common to all channels

Access Read/write
Data type Double
Read-only when Yes
running

Values

The default value is zero.

Examples

Create the analog output object **ao** for a sound card and add one channel to it.

```
ao = analogoutput('winsound');
chans = addchannel(ao,1);

To queue one second of data:
data = sin(linspace(0,10,8000))';
```

RepeatOutput

putdata(ao,data)

To continuously output data for 10 seconds:

set(ao, 'RepeatOutput',9)

See Also Functions

putdata

RTDConfiguration

Purpose

Specify wiring configuration of RTD device

Description

Use this property to specify the wiring configuration for measuring resistance.

When you create an RTD channel, the wiring configuration is unknown and the RTDConfiguration property displays Unknown. You must change this to one of the following valid configurations:

- TwoWire
- ThreeWire
- FourWire

See Also

Class

daq.Session

Properties

RO, RTDType

RTDType

Purpose

Specify sensor sensitivity

Description

Use this property to specify the sensitivity of a standard RTD sensor in the session-based interface. A standard RTD sensor is defined as a 100–ohm platinum sensor.

When you create an RTD channel, the sensitivity is unknown and the RTDType property displays Unknown. You must change this to one of these valid values:

- Pt3750
- Pt3851
- Pt3911
- Pt3916
- Pt3920
- Pt3928

See Also

daq.Session.addAnalogInputChannel

Class

daq.Session

Properties

RTDConfiguration, RTDType

Indicate whether device object is running

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Along with the Logging or Sending property, Running reflects the state of an analog input or analog output object. Running can be On or Off.

Running is automatically set to On once the start function is issued. When Running is On, you can acquire data from an analog input device or send data to an analog output device after the trigger occurs. For digital I/O objects, Running is typically used to indicate if time-based events are being generated.

Running is automatically set to Off once the stop function is issued, the specified data is acquired or sent, or a run-time error occurs. When Running is Off, you cannot acquire or send data. However, you can acquire one sample with the getsample function, or send one sample with the putsample function.

Characteristics

Usage AI, AO, DIO, common to all channels and lines

Access Read-only

Data type String
Read-only when N/A

running

Values

{Off} The device object is not running.

On The device object is running.

Running

See Also Functions

getsample, putsample, start

Properties

Logging, Sending

Specify callback function to execute when run-time error occurs

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

A run-time error event is generated immediately after a run-time error occurs. This event executes the callback function specified for RuntimeErrorFcn. Additionally, a toolbox error message is automatically displayed to the MATLAB workspace. If an error occurs that is not explicitly handled by the toolbox, then the hardware-specific error message is displayed.

The default value for RunTimeErrorFcn is daqcallback, which displays the event type, the time the event occurred, and the device object name along with the error message.

Run-time error event information is stored in the Type and Data fields of the EventLog property. The Type field value is Error. The Data field values are given below.

Data Field Value	Description
AbsTime	The absolute time (as a clock vector) the event occurred.
RelSample	The acquired (AI) or output (AO) sample number when the event occurred.
String	The descriptive error message.

Run-time errors include hardware errors and timeouts. Run-time errors do not include configuration errors such as setting an invalid property value.

RuntimeErrorFcn

 $\begin{array}{ccc} \textbf{Characteristics} & U_{sage} & & AI,\,AO,\,common\;to\;all\;channels \end{array}$

Access Read/write

Data type String

Read-only when No

running

Values The default value is daqcallback.

See Also Functions

daqcallback

Properties

EventLog, Timeout

Specify per-channel rate at which analog data is converted to digital data, or vice versa

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

SampleRate specifies the per-channel rate (in samples/second) that an analog input (AI) or analog output (AO) subsystem converts data. AI subsystems convert analog data to digital data, while AO subsystems convert digital data to analog data.

AI and AO subsystems have a finite (though often large) number of valid sampling rates. If you specify a sampling rate that does not match one of the valid values, the data acquisition engine automatically selects the nearest available sampling rate. In most data acquisition hardware, some valid sample rates can be non integers. See The Sampling Rate for more info about valid sample rates.

Because the engine can set the sampling rate to a value that differs from the value you specify, you should return the actual sampling rate using the get function or the device object display summary. Alternatively, you can use the setverify function, which sets the SampleRate value and then returns the actual value that is set. To find out the range of sampling rates supported by your board, use the propinfo function. Additionally, because the actual sampling rate depends on the number of channels contained by the device object and the ChannelSkew property value (AI only), SampleRate should be the last property you set before starting the device object.

Characteristics

Usage AI, AO, common to all channels

Access Read/write

SampleRate

Data type Double
Read-only when Yes
running

Values

The default value is obtained from the hardware driver.

Examples

Create the analog input object ai for a sound card and add two channels to it.

```
ai = analoginput('winsound');
addchannel(ai,1:2);
```

You can find out the range of valid sampling rates with the ConstraintValue field of the propinto function.

```
rates = propinfo(ai, 'SampleRate');
rates.ConstraintValue
ans =
    8000    48000
```

To configure the per-channel sampling rate to 48 kHz:

```
set(ai, 'SampleRate', 48000)
```

Alternatively, you can use the setverify function to configure and return the SampleRate value.

```
ActualRate = setverify(ai, 'SampleRate', 48000);
```

See Also Functions

propinfo, setverify

Properties

ChannelSkew

Indicate number of samples acquired per channel

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

SamplesAcquired is continuously updated to reflect the current number of samples acquired by an analog input object. It is reset to zero after a start function is issued.

Use the SamplesAvailable property to find out how many samples are available to be extracted from the engine.

Characteristics

Usage AI, common to all channels

Access Read-only
Data type Double
Read-only when N/A

running

Values

The value is continuously updated to reflect the current number of samples acquired. The default value is zero.

See Also

Functions

start

Properties

SamplesAvailable

SamplesAcquiredFcn

Purpose

Specify callback function to execute when predefined number of samples is acquired for each channel group member

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

A samples acquired event is generated immediately after the number of samples specified by the SamplesAcquiredFcnCount property is acquired for each channel group member. This event executes the callback function specified for SamplesAcquiredFcn.

The samples acquired event is executed regardless of its waiting time in the queue.

Use SamplesAcquiredFcn to trigger an event each time a specified number of samples is acquired. To process samples at regular time intervals, use the TimerFcn property.

Samples acquired event information is not stored in the EventLog property. When the callback function is executed, the second argument is a structure containing two fields. The Type field value is set to the string 'SamplesAcquired', and the Data field values are given below.

Data Field Value	Description
AbsTime	The absolute time (as a clock vector) the event occurred.
RelSample	The acquired sample number when the event occurred.

Characteristics Usage

Usage AI, common to all channels

Access Read/write

${\bf Samples Acquired Fcn}$

Data type String

Read-only when No

running

Values The default value is an empty string.

See Also Properties

 ${\tt EventLog}, {\tt SamplesAcquiredFcnCount}, {\tt TimerFcn}$

SamplesAcquiredFcnCount

Purpose

Specify number of samples to acquire for each channel group member before samples acquired event is generated

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

A samples acquired event is generated immediately after the number of samples specified by SamplesAcquiredFcnCount is acquired for each channel group member. This event executes the callback function specified by the SamplesAcquiredFcn property.

Characteristics

Usage AI, common to all channels

Access Read/write
Data type Double
Read-only when Yes

running

Values The default value is 1024.

See Also Properties

 ${\tt SamplesAcquiredFcn}$

Indicate number of samples available per channel in engine

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

For analog input (AI) objects, SamplesAvailable indicates the number of samples that can be extracted from the engine for each channel group member with the getdata function. For analog output (AO) objects, SamplesAvailable indicates the number of samples that have been queued with the putdata function, and can be sent (output) to each channel group member.

After data has been extracted (AI) or output (AO), the SamplesAvailable value is reduced by the appropriate number of samples. For AI objects, SamplesAvailable is reset to zero after a start function is issued.

For AI objects, use the SamplesAcquired property to find out how many samples have been acquired since the start function was issued. For AO objects, use the SamplesOutput property to find out how many samples have been output since the start function was issued.

Characteristics

Usage AI, AO, common to all channels

Access Read-only
Data type Double

Read-only when N/A

running

Values

The value is automatically updated based on the number of samples acquired (analog input) or sent (analog output). The default value is zero.

Samples Available

See Also Functions

start

Properties

SamplesAcquired, SamplesOutput

Indicate number of samples output per channel from engine

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

SamplesOutput is continuously updated to reflect the current number of samples output by an analog output object. It is reset to zero after the device objects stops and data has been queued with the putdata function.

Use the SamplesAvailable property to find out how many samples are available to be output from the engine.

Characteristics

Usage AO, common to all channels

Access Read-only
Data type Double
Read-only when N/A

running

Values

The value is continuously updated to reflect the current number of samples output. The default value is zero.

See Also Functions

putdata

Properties

SamplesAvailable

SamplesOutputFcn

Purpose

Specify callback function to execute when predefined number of samples is output for each channel group member

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

A samples output event is generated immediately after the number of samples specified by the SamplesOutputFcnCount property is output for each channel group member. This event executes the callback function specified for SamplesOutputFcn.

Use SamplesOutputFcn to trigger an event each time a specified number of samples is output. To process samples at regular time intervals, use the TimerFcn property.

Samples output event information is not stored in the EventLog property. When the callback function is executed, the second argument is a structure containing two fields. The Type field value is set to the string 'SamplesOutput', and the event Data field values are given below.

Data Field Value	Description
AbsTime	The absolute time (as a clock vector) the event occurred.
RelSample	The output sample number when the event occurred.

Characteristics Usage

Usage AO, common to all channels

Access Read/write

SamplesOutputFcn

Data type String

Read-only when No

running

Values The default value is an empty string.

See Also Properties

EventLog, SamplesOutputFcnCount

SamplesOutputFcnCount

Purpose

Specify number of samples to output for each channel group member before samples output event is generated

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

A samples output event is generated immediately after the number of samples specified by SamplesOutputFcnCount is output for each channel group member. This event executes the callback function specified by the SamplesOutputFcn property.

Characteristics

Usage AO, common to all channels

Access Read/write
Data type Double
Read-only when Yes

running

Values The default value is 1024.

See Also Properties

 ${\tt SamplesOutputFcn}$

Specify number of samples to acquire for each channel group member for each trigger that occurs

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

SamplesPerTrigger specifies the number of samples to acquire for each analog input channel group member for each trigger that occurs. If SamplesPerTrigger is set to Inf, then the analog input object continually acquires data until a stop function is issued or an error occurs.

The default value of SamplesPerTrigger is calculated by the data acquisition engine such that one second of data is acquired. This calculation is based on the value of SampleRate.

Characteristics

Usage AI, common to all channels

Access Read/write
Data type Double
Read-only when running Yes

Values

The default value is set by the engine such that one second of data is acquired.

Examples

Create the analog input object ai for a sound card and add two channels to it.

```
ai = analoginput('winsound');
addchannel(ai,1:2);
```

SamplesPerTrigger

By default, a one second acquisition in which 8000 samples are acquired for each channel is defined. To define a two second acquisition at the same sampling rate:

set(ai, 'SamplesPerTrigger',16000)

See Also Functions

stop

Properties

SampleRate

ScansAcquired

Purpose Number of scans acquired during operation

Description In the session-based interface, the ScansAcquired property displays

the number of scans acquired after you start the operation using

daq.Session.startBackground.

Values The read-only value represents the number of scans acquired

by the hardware. This value is reset each time you call

daq.Session.startBackground.

See Also Properties

NumberOfScans, ScansOutputByHardware

Methods

daq.Session.startBackground

Class

daq.Session

ScansOutputByHardware

Purpose

Indicate number of scans output by hardware

Description

In the session-based interface, the ScansOutputByHardware property displays the number of scans output by the hardware after you start the operation using daq.Session.startBackground.

Tip The value depends on information from the hardware.

Values

This read-only value is based on the output of the hardware configured for your session.

See Also

Properties

ScansAcquired, ScansQueued

Methods

daq.Session.queueOutputData, daq.Session.startBackground

Class

daq.Session

Purpose Indicate number of scans queued for output

Description In the session-based interface, the ScansQueued property displays the

number of scans queued for output daq.Session.queueOutputData. The ScansQueued property increases when you successfully call daq.Session.queueOutputData. The ScansQueued property decreases

when the hardware reports that it has successfully output data.

Values This read-only value is based on the number of scans queued.

See Also Properties

ScansOutputByHardware

Methods

daq.Session.queueOutputData

Class

daq.Session

Sending

Purpose

Indicate whether data is being sent to hardware device

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Along with the Running property, Sending reflects the state of an analog output object. Sending can be On or Off.

Sending is automatically set to On when a trigger occurs. When Sending is On, queued data is being output to the analog output subsystem.

Sending is automatically set to Off when the queued data has been output, an error occurs, or a stop function is issued. When Sending is Off, data is not being output to the analog output subsystem although you can output a single sample with the putsample function.

Characteristics

Usage AO, common to all channels

Access Read-only
Data type String
Read-only when running N/A

Values

{Off} Data is not being sent to the analog output hardware.

On Data is being sent to the analog output hardware.

See Also Functions

putsample

Properties

Running

Sensitivity of an analog channel

Description

When working with the session-based interface, the Sensitivity property to set the accelerometer or microphone sensor channel.

Sensitivity in an accelerometer channel is expressed as $\frac{v}{g}$, or volts per gravity.

Sensitivity in a microphone channel is expressed as $\frac{v}{pa}$, or volts per pascal.

Examples

Create a session object, add an analog input channel, with the 'accelerometer' MeasurementType.

```
s = daq.createSession('ni');
s.addAnalogInputChannel('Dev4', 'ai0', 'accelerometer')

Data acquisition session using National Instruments hardware:

Will run for 1 second (2000 scans) at 2000 scans/second.

Number of channels: 1

index Type Device Channel MeasurementType Range Name

1 ai Dev4 ai0 Accelerometer (PseudoDiff) -5.0 to +5.0 Volts
```

Change the Sensitivity to 10.2e-3 V/G:

```
s.Channels(1).Sensitivity = 10.2e-3

s =

Data acquisition session using National Instruments hardware:
    Will run for 1 second (2000 scans) at 2000 scans/second.
    Number of channels: 1
    index Type Device Channel MeasurementType Range Name
```

Sensitivity

1 ai Dev4 ai0 Accelerometer (PseudoDiff) -490 to +490 Gravities

See Also Class

daq.Session, daq.Session.addAnalogInputChannel

Specify range of data expected from sensor

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

You use SensorRange to scale your data to reflect the range you expect from your sensor. You can find the appropriate sensor range from your sensor's specification sheet.

The data is scaled while it is extracted from the engine with the getdata function according to the formula

$$scaled\ value = \frac{(A\ /\ D\ value)(units\ range)}{(sensor\ range)}$$

The A/D value is constrained by the InputRange property, which reflects the gain and polarity of your hardware channels. The units range is given by the UnitsRange property.

Characteristics

Usage AI, per channel

Access Read/write

Data type Two-element vector of doubles

Read-only when running No

Values

The default value is determined by the default value of the InputRange property.

See Also Functions

getdata

SensorRange

Properties

InputRange, Units, UnitsRange

ShuntLocation property

Purpose

Indicate location of channel's shunt resistor

Description

When working with the session-based interface, ShuntLocation on the analog input current channel indicates if the shunt resistor is located internally on the device or externally. Values are:

- 'Internal': when the shunt resistor is located internally.
- 'External': when the shunt resistor is located externally.

If your device supports an internal shunt resistor, this property is set to Internal by default. If the shunt location is external, you must specify the shunt resistance value.

See Also ShuntResistance |

ShuntResistance

Purpose

Resistance value of channel's shunt resistor

Description

When working with the session-based interface, the analog input current channel's ShuntResistance property indicates resistance in ohms. This value is automatically set if the shunt resistor is located internally on the device and is read only.

Note Before starting an analog output channel with an external shunt resistor, specify the shunt resistance value.

See Also ShuntLocation |

Purpose Indicates trigger source terminal

Description When working with the session-based interface, the Source property

indicates the device and terminal to which you added a trigger.

See Also Destinationdaq.Session.addTriggerConnection

Specify callback function to execute before device object runs

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

A start event is generated immediately after the start function is issued. This event executes the callback function specified for StartFcn. When the callback function has finished executing, Running is automatically set to On and the device object and hardware device begin executing. Note that the device object is not started if an error occurs while executing the callback function.

Start event information is stored in the Type and Data fields of the EventLog property. The Type field value is Start. The Data field values are given below.

Data Field Value	Description
AbsTime	The absolute time (as a clock vector) the event occurred.
RelSample	The acquired (AI) or output (AO) sample number when the event occurred.

Characteristics

Usage AI, AO, common to all channels

Access Read/write
Data type String
Read-only when No

running

Values The default value is an empty string.

See Also Functions

start

Properties

EventLog, Running

Specify callback function to execute after device object runs

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

A stop event is generated immediately after the device object and hardware device stop executing. This occurs when

- A stop function is issued.
- For analog input (AI) objects, the requested number of samples to acquire was reached or data was missed. For analog output (AO) objects, the requested number of samples to output was reached.
- A run-time error occurred.

A stop event executes the callback function specified for StopFcn. Under most circumstances, the callback function is not guaranteed to complete execution until sometime after the device object and hardware device stop, and the Running property is set to Off.

Stop event information is stored in the Type and Data fields of the EventLog property. The Type field value is Stop. The Data field values are given below.

Data Field Value	Description
AbsTime	The absolute time (as a clock vector) the event occurred.
RelSample	The acquired (AI) or output (AO) sample number when the event occurred.

Characteristics Usage AI, AO, common to all channels

Access Read/write

Data type String
Read-only when No

running

Values The default value is an empty string.

See Also Functions

stop

Properties

EventLog, Running

Specify device object label

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Tag provides a means to identify device objects with a label. Using the daqfind function and the Tag value, you can identify and retrieve a device object that was cleared from the MATLAB workspace.

Characteristics

Usage AI, AO, DIO, common to all channels and lines

Access Read/write

Data type String

Read-only when No

running

Values

The default value is an empty string.

Examples

Create the analog input object ai for a sound card and add two channels to it.

```
ai = analoginput('winsound');
addchannel(ai,1:2);
```

Assign ai a label using Tag.

```
set(ai, 'Tag', 'Sound')
```

If ai is cleared from the workspace, you can use daqfind and the Tag value to identify and retrieve the device object.

```
clear ai
aicell = daqfind('Tag','Sound');
ai = aicell{1};
```

See Also Functions

daqfind

Terminal

Purpose PFI terminal of counter subsystem

Description When working with the session-based interface, the Terminal property

indicates the counter subsystem's corresponding PFI terminal.

See Also Class

daq.Session, daq.Session.AddCounterInputChannel,

 ${\tt daq.Session.addCounterOutputChannel},$

Purpose Specia

Specify terminal configuration

Description

Use the TerminalConfig to change the configuration of your analog channel. The property displays the hardware default configuration. You can change this to

- SingleEnded
- NonReferencedSingleEnded
- Differential
- PseudoDifferential

See Also

TerminalConfigsAvailable | daq.Session.addAnalogInputChannel
| daq_ref.Session.addAnalogOutputChannel

Terminals

Purpose

Terminals available on device or CompactDAQ chassis

Description

When working with the session-based interface, the Terminals on the device or the CompactDAQ chassis lists all available terminals. The list includes terminals available for trigger and clock connections. When you access the Terminals property on modules on a CompactDAQ chassis, the terminals are on the chassis, not on the module.

Examples

Discover available devices:

d=daq.getDevices

d =

Data acquisition devices:

index	Vendor	Device ID		Description		
1	ni	cDAQ1Mod1	National	Instruments	ΝI	9205
2	ni	cDAQ1Mod2	National	Instruments	ΝI	9263
3	ni	cDAQ1Mod3	National	Instruments	ΝI	9234
4	ni	cDAQ1Mod4	National	Instruments	ΝI	9201
5	ni	cDAQ1Mod5	National	Instruments	ΝI	9402
6	ni	cDAQ1Mod6	National	Instruments	ΝI	9213
7	ni	cDAQ1Mod7	National	Instruments	ΝI	9219
8	ni	cDAQ1Mod8	National	Instruments	ΝI	9265

Access the Terminals property of NI 9205 with index 1:

```
d(1).Terminals
```

```
ans =
   'cDAQ1/PFI0'
   'cDAQ1/PFI1'
   'cDAQ1/20MHzTimebase'
   'cDAQ1/80MHzTimebase'
```

Terminals

```
'cDAQ1/ChangeDetectionEvent'
'cDAQ1/AnalogComparisonEvent'
'cDAQ1/100kHzTimebase'
'cDAQ1/SyncPulse0'
'cDAQ1/SyncPulse1'
.
```

See Also Class

daq.Session

Functions

daq.getDevices,
daq.Session.addTriggerConnection,daq.Session.addClockConnection

ThermocoupleType

Purpose

Select thermocouple type

Description

When working with the session-based interface, use the ThermocoupleType property to select the type of thermocouple you will use to make your measurements. Select the type based on the temperature range and sensitivity you need.

Values

You can set the ThermocoupleType to:

- 'J'
- 'K'
- 'N'
- 'R'
- 'S'
- 'T'
- 'B'
- 'E'

By default the thermocouple type is 'Unknown'.

Examples

Create a session and add an analog input channel with 'Thermocouple' measurement type:

```
s = daq.createSession('ni');
s.addAnalogInputChannel('cDAQ1Mod6','ai1','Thermocouple');
ans =
```

Data acquisition analog input voltage channel 'ai0' on device 'cDAQ1Mod6'

Units: Celsius ThermocoupleType: Unknown Coupling: DC

TerminalConfig: Differential

ThermocoupleType

```
Range: 0 to +750 Celsius
Name: ''
ID: 'ai0'
Device: [1x1 daq.ni.CompactDAQModule]
MeasurementType: 'Voltage'

Set the ThermocoupleType property to 'J':
s.Channels.Thermocoupletype = 'J';
```

See Also Methods

daq.Session.addAnalogInputChannel

Class

daq.Session

Timeout

Purpose

Specify additional waiting time to extract or queue data

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

The Timeout value (in seconds) is added to the time required to extract data from the engine or queue data to the engine. Because data is extracted with the getdata function, and queued with the putdata function, Timeout is associated only with these two "blocking" functions.

If the requested data is not extracted or queued after waiting the required time, then a time-out condition occurs and control is immediately returned to the MATLAB workspace. A time-out is one of the conditions for stopping an acquisition. When a time-out occurs, the callback function specified by RuntimeErrorFcn is called.

Timeout is not associated with hardware time-out conditions. Possible hardware time-out conditions include

- Triggering on a voltage level and that level never occurs
- Externally clocking an acquisition and the external clock signal never occurs
- Losing the hardware connection

To check for hardware timeouts, you might need to poll the appropriate property value.

Characteristics

Usage AI, AO, common to all channels

Access Read/write

Data type Double

Read-only when Yes

running

Timeout

Values The default value is one second.

See Also Functions

getdata, putdata

Properties

RuntimeErrorFcn

TimerFcn

Purpose

Specify callback function to execute when predefined time period passes

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

A timer event is generated whenever the time specified by the TimerPeriod property passes. This event executes the callback function specified for TimerFcn. Time is measured relative to when the device object starts running.

Some timer events might not be processed if your system is significantly slow or if the TimerPeriod value is too small. The time taken to process an event depends on the sample rate, the performance of your system, and the data itself.

There can only be one timer event waiting in the queue at a given time. The callback function must process all available data to ensure that it keeps up with the inflow of data. Alternatively, you can use the SamplesAcquiredFcn (analog input) or SamplesOutputFcn (analog output) property to process the data when a specified number of samples is acquired or output.

Note For analog input objects, use the SamplesAvailable property inside a callback function to determine the number of samples available in the queue.

For digital I/O objects, timer events are typically used to update and display the state of the device object.

Timer event information is not stored in the EventLog property. When the callback function is executed, the second argument is a structure containing two fields. The Type field value is set to the string 'Timer', and the event Data field value is given below.

Data Field Value	Description
AbsTime	The absolute time (as a clock vector) the event occurred.

Characteristics Usage AI, AO, DIO, common to all channels and lines

Access Read/write
Data type String

Read-only when No running

Values The default value is an empty string.

See Also Properties

 ${\bf EventLog, Samples Acquired Fcn, Samples Output Fcn, Timer Period}$

TimerPeriod

Purpose

Specify time period between timer events

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

TimerPeriod specifies the time, in seconds, that must pass before the callback function specified for TimerFcn is called. Time is measured relative to when the hardware device starts running.

Some timer events might not be processed if your system is significantly slowed or if the TimerPeriod value is too small. For example, a common application for timer events is to display data. However, because displaying data is a CPU-intensive task, some of these events might be dropped.

Characteristics

Usage AI, AO, DIO, common to all channels and lines

Access Read/write
Data type Double

Read-only when

running

No

Values

The default value is 0.1 second.

See Also

Properties

TimerFcn

Specify channel serving as trigger source

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

TriggerChannel specifies the channel serving as the trigger source. The trigger channel must be specified before the trigger type. You might need to configure the TriggerCondition and TriggerConditionValue properties in conjunction with TriggerChannel.

For all supported vendors, if TriggerType is Software, then you must acquire data from the channel being used for the trigger source. For National Instruments hardware, if TriggerType is HwAnalogChannel, then TriggerChannel must be the first element of the channel group. The exception is if you are using simultaneous acquisition devices such as the S-series boards, with which you can specify any channel for the TriggerChannel value.

Characteristics

Usage AI, common to all channels

Access Read/write

Data type Vector or scalar

Read-only when Yes

running

Values

The data type can be either vector or scalar, representing one channel. The default value is an empty vector.

Examples

Create the analog input object ai, add two channels, and define the trigger source as channel 2.

```
ai = analoginput('winsound');
ch = addchannel(ai,1:2);
```

TriggerChannel

```
set(ai, 'TriggerChannel',ch(2))
set(ai, 'TriggerType', 'Software')
```

See Also Properties

 ${\tt TriggerCondition, TriggerConditionValue, TriggerType}$

Specify condition that must be satisfied before trigger executes

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

When working with the session-based interface, use the TriggerCondition property to specify the signal condition that executes the trigger, which synchronizes operations on devices in a session. For more information, see .

The available trigger conditions depend on the value of TriggerType. If TriggerType is Immediate or Manual, the only available TriggerCondition is None. If TriggerType is Software, then TriggerCondition can be Rising, Falling, Leaving, or Entering. These trigger conditions require one or more voltage values to be specified for the TriggerConditionValue property.

Based on the hardware you are using, additional trigger conditions might be available. Refer to the values listed below.

Values

All Supported Hardware

The following trigger condition is used when TriggerType is Immediate or Manual.

{None} No trigger condition is required.

The following trigger conditions are available when TriggerType is Software.

TriggerCondition

{Rising} The trigger occurs when the signal has a

positive slope when passing through the

specified value.

Falling The trigger occurs when the signal has a

negative slope when passing through the

specified value.

Leaving The trigger occurs when the signal leaves the

specified range of values.

Entering The trigger occurs when the signal enters the

specified range of values.

Measurement Computing

The following trigger conditions are available when TriggerType is HwDigital.

GateHigh The trigger occurs as long as the digital signal

is high.

GateLow The trigger occurs as long as the digital signal

is low.

TrigHigh The trigger occurs when the digital signal is

high.

TrigLow The trigger occurs when the digital signal is

low.

TrigPosEdge The trigger occurs when the positive (rising)

edge of the digital signal is detected.

{TrigNegEdge} The trigger occurs when the negative (falling)

edge of the digital signal is detected.

The following trigger conditions are available when TriggerType is HwAnalog.

TriggerCondition

{TrigAbove} The trigger occurs when the analog signal

makes a transition from below the specified

value to above.

TrigBelow The trigger occurs when the analog signal

makes a transition from above the specified

value to below.

GateNegHys The trigger occurs when the analog signal

is more than the specified high value. The acquisition stops if the analog signal is less

than the specified low value.

GatePosHys The trigger occurs when the analog signal is less

than the specified low value. The acquisition stops if the analog signal is more than the

specified high value.

GateAbove The trigger occurs as long as the analog signal

is more than the specified value.

GateBelow The trigger occurs as long as the analog signal

is less than the specified value.

GateInWindow The trigger occurs as long as the analog signal

is within the specified range of values.

GateOutWindow The trigger occurs as long as the analog signal

is outside the specified range of values.

National Instruments

The following trigger conditions are available for AI objects when TriggerType is HwDigital.

PositiveEdge The trigger occurs when the positive (rising)

edge of a digital signal is detected.

{NegativeEdge} The trigger occurs when the negative (falling)

edge of a digital signal is detected.

TriggerCondition

The following trigger conditions are available for AO objects on NI-DAQmx devices when TriggerType is HwDigital.

PositiveEdge The trigger occurs when the positive (rising)

edge of a digital signal is detected.

{NegativeEdge} The trigger occurs when the negative (falling)

edge of a digital signal is detected.

The following trigger conditions are available when TriggerType is HwAnalogChannel or HwAnalogPin.

{AboveHighLevel} The trigger occurs when the analog signal is

above the specified value.

BelowLowLevel The trigger occurs when the analog signal is

below the specified value.

InsideRegion The trigger occurs when the analog signal is

inside the specified region.

LowHysteresis The trigger occurs when the analog signal is

less than the specified low value with hysteresis

given by the specified high value.

HighHysteresis The trigger occurs when the analog signal

is greater than the specified high value with hysteresis given by the specified low value.

See Also Properties

TriggerChannel, TriggerConditionValue, TriggerType

Specify condition that must be satisfied before trigger executes

Description

When working with the session-based interface, use the TriggerCondition property to specify the signal condition that executes the trigger, which synchronizes operations on devices in a session. For more information, see .

Values

Set the trigger condition to RisingEdge or FallingEdge.

Examples

Create a session and add channels and trigger to the session.

```
s = daq.createSession('ni');
s.addAnalogInputChannel('Dev1', 0, 'voltage');
s.addAnalogInputChannel('Dev2', 0, 'voltage');
s.addTriggerConnection('Dev1/PFI4','Dev2/PFI0','StartTrigger');
```

Change the trigger condition to FallingEdge.

```
s.Connections(1).TriggerCondition='FallingEdge'

s =

Data acquisition session using National Instruments hardware:
  Will run for 1 second (1000 scans) at 1000 scans/second.

Trigger Connection added. (Details)

Number of channels: 2
  index Type Device Channel MeasurementType Range Name

1 ai Dev1 ai0 Voltage (Diff) -10 to +10 Volts
2 ai Dev2 ai0 Voltage (Diff) -10 to +10 Volts
```

See Also

dag.Session.addTiggerConnection

Properties

TriggerType

TriggerConditionValue

Purpose

Specify voltage value(s) that must be satisfied before trigger executes

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

For all hardware TriggerConditionValue is used when TriggerType is Software, and is ignored when TriggerCondition is None. For vendor specific triggers, refer to the TriggerCondition and the TriggerType properties.

To execute a software trigger, the values specified for TriggerCondition and TriggerConditionValue must be satisfied. When TriggerCondition is Rising or Falling, TriggerConditionValue accepts a single value. When TriggerCondition is Entering or Leaving, TriggerConditionValue accepts a two-element vector of values. For vendor specific values, refer to the TriggerCondition property.

Characteristics

Usage AI, common to all channels

Access Read/write

Data type Double (or a two-element vector of doubles)

Read-only when Yes

running

Values

The default value is zero.

Examples

Create the analog input object ai and add one channel to it.

```
ai = analoginput('winsound');
ch = addchannel(ai,1);
```

TriggerConditionValue

The trigger executes when a signal with a negative slope passing through 0.2 volts is detected on channel 1.

```
set(ai, 'TriggerChannel',ch)
set(ai, 'TriggerType', 'Software')
set(ai, 'TriggerCondition', 'Falling')
set(ai, 'TriggerConditionValue',0.2)
```

Create the analog input object ai for a National Instruments device and add four channels to it.

```
ai = analoginput('nidaq', 'Dev1');
ch = addchannel(ai,0:3);
```

The trigger executes when a signal with a positive slope passing through 4.5 volts is detected on PFI2.

```
set(ai, 'TriggerType', 'HwDigital')
set(ai, 'HwDigitalTriggerSource', 'PFI2')
set(ai, 'TriggerCondition', 'PositiveEdge')
set(ai, 'TriggerConditionValue', 4.5)
```

See Also Properties

TriggerCondition, TriggerType

TriggerDelay

Purpose

Specify delay value for data logging

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

You can define both pretriggers and postriggers. Pretriggers are specified with a negative TriggerDelay value while postriggers are specified with a positive TriggerDelay value. You can delay a trigger in units of time or samples with the TriggerDelayUnits property. Pretriggers are not defined for hardware triggers or when TriggerType is Immediate.

Pretrigger samples are included as part of the total samples acquired per trigger as specified by the SamplesPerTrigger property. If sample-time pairs are returned to the workspace with the getdata function, then the pretrigger samples are identified with negative time values.

Characteristics

Usage

AI, common to all channels

Access

Read/write

Data type

Double

Read-only when

Yes

running

Values

The default value is zero.

Examples

Create the analog input object ai and add one channel to it.

```
ai = analoginput('winsound');
ch = addchannel(ai,1);
```

Configure at to acquire 44,100 samples per trigger with 11,025 samples (0.25 seconds) acquired as pretrigger data.

```
set(ai, 'SampleRate',44100)
set(ai, 'TriggerType', 'Manual')
set(ai, 'SamplesPerTrigger',44100)
set(ai, 'TriggerDelay', -0.25)
```

See Also Properties

 ${\tt SamplesPerTrigger}, {\tt TriggerDelayUnits}$

TriggerDelayUnits

Purpose

Specify units in which trigger delay data is measured

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

TriggerDelayUnits can be Seconds or Samples. If TriggerDelayUnits is Seconds, then data logging is delayed by the specified time for each channel group member. If TriggerDelayUnits is Samples, then data logging is delayed by the specified number of samples for each channel group member.

The trigger delay value is given by the TriggerDelay property.

Characteristics

Usage AI, common to all channels

Access Read/write

Data type String

Read-only when Yes

running

Values

{Seconds} The trigger is delayed by the specified number of

seconds.

Samples The trigger is delayed by the specified number of

samples.

See Also Properties

TriggerDelay

Specify callback function to execute when trigger occurs

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

A trigger event is generated immediately after a trigger occurs. This event executes the callback function specified for TriggerFcn. Under most circumstances, the callback function is not guaranteed to complete execution until sometime after Logging is set to On for analog input (AI) objects, or Sending is set to On for analog output (AO) objects.

Trigger event information is stored in the Type and Data fields of the EventLog property. The Type field value is Trigger. The Data field values are given below.

	Data Field Value	Description
	AbsTime	The absolute time (as a clock vector) the event occurred.
	RelSample	The acquired (AI) or output (AO) sample number when the event occurred.
	Channel	The index number for each input channel serving as a trigger source (AI only).
	Trigger	The trigger number.
5	Usage	AI, AO, common to all channels

Characteristics

Usage	AI, AO, common to all channels
Access	Read/write
Data type	String
Read-only when running	No

TriggerFcn

Values The default value is an empty string.

See Also Functions

trigger

Properties

EventLog, Logging

Specify number of additional times trigger executes

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

You can configure a trigger to occur once (one-shot acquisition) or multiple times. If TriggerRepeat is set to its default value of zero, then the trigger executes once. If TriggerRepeat is set to a positive integer value, then the trigger executes once, and is repeated the specified number of times. For example, if the value is set to 2, you will get a total of 3 triggers. If TriggerRepeat is set to inf then the trigger executes continuously until a stop function is issued or an error occurs.

You can quickly evaluate how many triggers have executed by examining the TriggersExecuted property or by invoking the display summary for the device object. The display summary is invoked by typing the device object name at the MATLAB Command Window.

Note We have observed that National Instruments USB devices have a significant cycle time for the communications required to trigger the device. If you are using an NI USB device, we recommend that you set up longer acquisitions that use fewer triggers. That is, increase SamplesPerTrigger and decrease TriggerRepeat.

Characteristics

Usage AI, common to all channels

Access Read/write

Data type Double

Read-only when Yes

running

TriggerRepeat

Values The default value is zero.

See Also Functions

disp, stop

Properties

SamplesPerTrigger, TriggersExecuted, TriggerType

Indicate number of triggers that execute

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

You can find out how many triggers executed by returning the value of TriggersExecuted. The trigger number for each trigger executed is also recorded by the Data.Trigger field of the EventLog property.

Characteristics

Usage AI, AO, common to all channels

Access Read-only
Data type Double
Read-only when N/A

running

Values

The default value is zero.

Examples

Create the analog input object ai and add one channel to it.

```
ai = analoginput('winsound');
ch = addchannel(ai,1);
```

Configure at to acquire 40,000 samples with five triggers using the default sampling rate of 8000 Hz.

```
set(ai, 'TriggerRepeat',4)
start(ai)
```

TriggersExecuted returns the number of triggers executed.

```
ai.TriggersExecuted
ans =
5
```

TriggersExecuted

See Also Properties

EventLog

TriggersPerRun

Purpose Indicate the number of times the trigger executes in an operation

Description When working with the session-based interface, the TriggersPerRun

property indicates the number of times the specified trigger executes for

one acquisition or generation session.

See Also daq.Session.addTriggerConnection

Triggers Remaining

Purpose Indicates the number of trigger to execute in an operation

Description When working with the session-based interface, the TriggersRemaining

property indicates the number of trigger remaining for this acquisition or generation session. This value depends on the number of triggers

set using TriggersPerRun.

See Also daq.Session.addTriggerConnection

Specify type of trigger to execute

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

TriggerType can be Immediate, Manual, or Software. If TriggerType is Immediate, the trigger occurs immediately after the start function is issued. If TriggerType is Manual, the trigger occurs immediately after the trigger function is issued. If TriggerType is Software, the trigger occurs when the associated trigger condition is satisfied (AI only).

For a given hardware device, additional trigger types might be available. Some trigger types require trigger conditions and trigger condition values. Trigger conditions are specified with the TriggerCondition property, while trigger condition values are specified with the TriggerConditionValue property.

When a trigger occurs for an analog input object, data logging is initiated and the Logging property is automatically set to On. When a trigger occurs for an analog output object, data sending is initiated and the Sending property is automatically set to On.

Characteristics

Usage AI, AO, common to all channels

Access Read/write

Data type String
Read-only when Yes

running

TriggerType

Values All Supported Hardware

{Immediate} The trigger executes immediately after start is

issued. Pretrigger data cannot be captured.

Manual The trigger executes immediately after the

trigger function is issued.

Software The trigger executes when the associated

trigger condition is satisfied. Trigger conditions are given by the TriggerCondition property.

(AI only).

Measurement Computing

HwDigital The trigger source is an external digital signal

(AI only). Pretrigger data cannot be captured.

HwAnalog The trigger source is an external analog signal

(AI only).

National Instruments

HwDigital The trigger source is an external digital signal.

Pretrigger data cannot be captured. Control the trigger source with HwDigitalTriggerSource property. Specify the external digital

signal with the TriggerCondition and TriggerConditionValue properties.

HwAnalogChannel The trigger source is an external analog

signal (AI only). To set the trigger source, see

TriggerChannel property.

HwAnalogPin The trigger source is a low-range external

analog signal (AI only). Note that HwAnalogPin is supported only for Traditional NIDAQ devices. It is not supported for NIDAQmx

devices.

For 1200 Series hardware, HwDigital is the only device-specific TriggerType value for analog input subsystems. Analog output subsystems do not support any device-specific TriggerType values.

Note The Traditional NI-DAQ adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition ToolboxTM object for Traditional NI-DAQ adaptor beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release. See the supported hardware page at www.mathworks.com/products/daq/supportedio.html for more information.

See Also Functions

start, trigger

Properties

Logging, Sending, TriggerChannel, TriggerCondition, TriggerConditionValue

TriggerType

Purpose Type of trigger executed

Description When working with the session-based interface, use this read only

property displays the type of trigger that the source device executes to synchronize operations in the session. Currently all trigger types

re digital.

See Also Functions

daq.Session.addTriggerConnection

Properties

TriggerCondition

Indicate device object type, channel, or line

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Type is associated with device objects, channels, and lines. For device objects, Type can be Analog Input, Analog Output, or Digital I/O. Once a device object is created, the value of Type is automatically defined.

For channels, the only value of Type is Channel. For lines, the only value of Type is Line. The value is automatically defined when channels or lines are added to the device object.

Characteristics

Values

Usage AI, AO, common to all channels and per

channel; DIO, common to all lines and per

line

N/A

Access Read-only

Data type String

Read-only when

running

Device Objects

For device objects, Type has these possible values:

Analog Input The device object type is analog input.

Analog Output The device object type is analog output.

Digital IO The device object type is digital I/O.

The value is automatically defined after the device object is created.

Channels and Lines

For channels, the only value of Type is Channel. For lines, the only value of Type is Line. The value is automatically defined when channels or lines are added to the device object.

Purpose Display synchronization trigger type

Description When working with the session-based interface, this property displays

the trigger type

Characteristics Usage AI, AO, common to all channels and per

channel; DIO, common to all lines and per

line

Access Read-only

Data type String
Read-only when N/A

running

Values Device Objects

For device objects, Type has these possible values:

Analog Input The device object type is analog input.

Analog Output The device object type is analog output.

Digital IO The device object type is digital I/O.

The value is automatically defined after the device object is created.

Channels and Lines

For channels, the only value of Type is Channel. For lines, the only value of Type is Line. The value is automatically defined when channels or lines are added to the device object.

Units

Purpose

Specify unit of RTD measurement

Description

Use this property to specify the temperature unit of the analog input channel with RTD measurement type in the session-based interface.

You can specify temperature values as:

- Celsius (Default)
- Fahrenheit
- Kelvin
- Rankine

See Also

Class

daq.Session

Specify engineering units label

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Units is a string that specifies the engineering units label to associate with your data. You should use Units in conjunction with the UnitsRange property.

Characteristics

Usage AI, AO, per channel

Access Read/write

Data type String

Read-only when No

running

Values The default value is Volts.

See Also Properties

UnitsRange

UnitsRange

Purpose

Specify range of data as engineering units

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

You use UnitsRange to scale your data to reflect particular engineering units.

For analog input objects, the data is scaled while it is extracted from the engine with the getdata function according to the formula

scaled value = (A/D value)(units range)/(sensor range)

The A/D value is constrained by the InputRange property, which reflects the gain and polarity of your analog input channels. The sensor range is given by the SensorRange property, which reflects the range of data you expect from your sensor.

For analog output objects, the data is scaled when it is queued in the engine with the putdata function according to the formula

scaled value = (original value)(output range)/(units range)

The output range is constrained by the OutputRange property, which specifies the gain and polarity of your analog output channels.

For both objects, you can also use the Units property to associate a meaningful label with your data.

Characteristics

Usage AI, AO, per channel

Access Read/write

Data type Two-element vector of doubles

Read-only when

running

Values The default value is determined by the default value of the InputRange

or the OutputRange property.

See Also Functions

getdata, putdata

Properties

 ${\tt InputRange,\,OutputRange,\,SensorRange,\,Units}$

Store data to associate with device object

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

UserData stores data that you want to associate with the device object.

Note that if you return analog input object information to the MATLAB workspace using the dagread function, the UserData value is not restored.

Characteristics

Usage AI, AO, DIO, common to all channels and

lines

Access Read/write
Data type Any type

Read-only when No

running

Values

The default value is an empty vector.

Examples

Create the analog input object ai and add two channels to it.

```
ai = analoginput('nidaq','Dev1');
addchannel(ai,0:1);
```

Suppose you want to access filter coefficients during the acquisition. You can create a structure to store these coefficients, which can then be stored in UserData.

```
coeff.a = 1.0;
coeff.b = -1.25;
set(ai, 'UserData', coeff)
```

Vendor information associated with session object

Description

In the session-based interface, the Vendor property displays information about the vendor.

Examples

Use the dag.getVendors to get information about vendors.

```
s = daq.createSession('ni');
v = s.Vendor
v =
```

Data acquisition vendor 'National Instruments':

```
ID: 'ni'
```

FullName: 'National Instruments'

AdaptorVersion: '2.17 (R2010b)'
DriverVersion: '9.1 NI-DAQmx'

IsOperational: true

Values

a daq. Vendor object that represents the vendor associated with the session.

See Also

Properties

ScansQueued

Methods

daq.createSession

Class

daq.Session

ZResetCondition

Purpose

Reset condition for Z-indexing

Description

When working with the session-based interface, use the ZResetCondition property to specify reset conditions for Z-indexing of counter Input 'Position' channels. Accepted values are:

- 'BothHigh'
- 'BothLow'
- 'AHigh'
- 'BHigh'

See Also

Class

daq.Session, daq.Session.addCounterInputChannel

ZResetEnable

Purpose Enable reset for Z-indexing

Description When working with the session-based interface, use the ZResetEnable

property to specify if you will allow the Z-indexing to be reset on a

counter input 'Position' channel.

See Also Class

 ${\tt daq.Session.addCounterInputChannel}$

ZResetValue

Purpose Reset value for Z-indexing

Description When working with the session-based interface, use the ZResetValue

property to specify the reset value for Z-indexing on a counter input

'Position' channel.

See Also Class

daq.Session, daq.Session.addCounterInputChannel

Device-Specific Properties

Groups by Vendor

Device-specific properties apply only to hardware devices of a specific type or from a specific vendor. For example, the BitsPerSample property is supported only for sound cards, while the NumMuxBoards property is supported only for National Instruments devices.

The properties are grouped according to these supported vendors:

```
"Advantech" on page 3-3
```

You can display device-specific properties with the set function. The device-specific properties are displayed after the base properties.

Note Some device-specific property values are not available for all devices. Refer to your hardware documentation for detailed information about device-specific behavior.

[&]quot;Measurement Computing" on page 3-4

[&]quot;National Instruments" on page 3-5

[&]quot;Parallel Port" on page 3-7

[&]quot;Sound Card" on page 3-8

Advantech

Property Name	Description	Device Objects
TransferMode	Specify how data is transferred from data acquisition device to system memory	AI, AO

Measurement Computing

Property Name	Description	Device Objects
OutOfDataMode	Specify how value held by analog output subsystem is determined	AO
TransferMode	Specify how data is transferred from data acquisition device to system memory	AI, AO

National Instruments

Use this table to map to the National Instruments terminology.

Data Acquisition Toolbox	NI_DAQmx
Scan Clock	Sample Clock
Sample Clock	Convert Clock

Property Name	Description	Device Objects
Coupling	Specify input coupling mode	AI
ExternalClockDriveLine	Specify which signal is driven by the clock indicating that an analog output update has occurred	AO
ExternalClockSource	Specify which signal generates an analog output update across channels	AO
ExternalSampleClockDriveL	Specify which signal line is driven by the clock for sample conversions on each channel	AI
ExternalSampleClockSource	Specify which signal provides clock for sample conversions across channels	AI
ExternalScanClockDriveLine	Specify which signal is driven by the clock indicating the start of a series of conversions across channels	AI
ExternalScanClockSource	Specify which signal starts series of conversions across channels	AI
ExternalTriggerDriveLine	Specify which signal line is driven with a pulse when data acquisition or generation starts	AI

Property Name	Description	Device Objects
HwDigitalTriggerSource	Specify which signal initiates data acquisition	AI, AO
NumMuxBoards	Specify number of external multiplexer devices connected	AI
OutOfDataMode	Specify how value held by analog output subsystem is determined	AO
TransferMode	Specify how data is transferred from data acquisition device to system memory	AI, AO

Parallel Port

Property Name	Description	Device Objects
BiDirectionalBit	Specify BIOS control register bit that determines bidirectional operation	DIO
PortAddress	Indicate base address of parallel port	DIO

Sound Card

Property Name	Description	Device Objects
BitsPerSample	Specify number of bits sound card uses to represent each sample	AI, AO
StandardSampleRates	Specify whether valid sample rates snap to small set of standard values, or if you can set sample rate to any value within allowed bounds	AI, AO

Getting Command-Line Property Help

To get command-line property help, you should use the daqhelp function. For example, to get help for the sound card's BitsPerSample property:

daghelp BitsPerSample

Note You can use daghelp without creating a device object.

You can also get property characteristics, such as the default property value, using the propinfo function. For example, suppose you create the analog input object ai for a sound card and want to find the default value for the BitsPerSample property.

```
ai = analoginput('winsound');
out = propinfo(ai, 'BitsPerSample');
out.DefaultValue
ans =
    16
```

Device-Specific Properties — Alphabetical List

BiDirectionalBit

Purpose

Specify BIOS control register bit that determines bidirectional operation

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

BiDirectionalBit can be 5, 6, or 7. The default value is 5 because most parallel port hardware uses bit 5 of the BIOS control register to determine the direction (input or output) of port 0.

If port 0 is unable to input data, you need to configure the BiDirectionalBit value to 6 or 7. Typically, you will not know the bit value required by your port, and some experimentation is required.

Note The Parallel Port adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition Toolbox[™] object for 'parallel' beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release. See the supported hardware page at www.mathworks.com/products/dag/supportedio.html for more information.

Characteristics

Vendor Parallel port

Usage DIO, common to all lines

Access Read/write

Data type Double Read-only when Yes

running

BiDirectionalBit

Values

 $\{5\},\,6,\,\mathrm{or}\ 7$

The BIOS control register bit that determines bidirectional operation.

BitsPerSample

Purpose

Specify number of bits sound card uses to represent samples

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

BitsPerSample can be 8, 16, or any value between 17 and 32. The specified number of bits determines the number of unique values a sample can take on. For example, if BitsPerSample is 8, the sound card represents each sample with 8 bits. This means that each sample is represented by a number from 0 through 255. If BitsPerSample is 16, the sound card represents each sample with 16 bits. This means that each sample is represented by a number from 0 through 65,535.

For older Sound Blaster cards configured for full duplex operation, you might not be able to set BitsPerSample to 16 bits for both the analog input and analog output subsystems. Instead, you need to set one subsystem for 8 bits, and the other subsystem for 16 bits.

Note To use the high-resolution (greater than 16 bit) capabilities for some sound cards, you might need to configure BitsPerSample to either 24 or 32 even if your device does not use that number of bits.

Characteristics

Vendor Sound cards

Usage AI, AO, common to all channels

Access Read/write

Data type Double

Read-only when running Yes

BitsPerSample

Values

 $8, \{16\}, or 17-32$

Represent data with the specified number of bits.

Coupling

Purpose

Specify input coupling mode

Description

The Coupling property is visible only if the device you are using supports coupling and the value can be changed. Coupling can be DC or AC. If Coupling is DC, the input is connected directly to the amplifier. If Coupling is AC, a series capacitor is inserted between the input connector and the amplifier.

When AC coupling is selected, the DC bias component of the measured signal is filtered out of the waveform by the hardware. This is typically used with dynamic signals such as audio. When DC coupling is selected, the complete signal including the DC bias component is measured. This is typically used with slowly changing signals such as temperature or voltage readings.

Values

{AC} A series capacitor is inserted between the

input connector and the amplifier.

DC The input is connected directly to the amplifier.

The default is set to AC for

- National Instruments devices that use the NI-DAQmx interface and support AC coupling
- National Instruments DSA cards using the Traditional NI-DAQ interface

Note The Traditional NI-DAQ adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition Toolbox™ object for Traditional NI-DAQ adaptor beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release. See the supported hardware page at www.mathworks.com/products/daq/supportedio.html for more information.

In all other cards, the default is set to DC.

Examples

In the session-based interface, create a session and add an analog input channel.

```
s = daq.createSession('ni');
ch = s.addAnalogInputChannel('Dev4', 'ai1', 'Voltage')
Change the coupling type to DC:
ch.Coupling = 'DC';
```

In the legacy interface, create the analog input object ai for a National Instruments board, and add a hardware channel to it.

```
ai = analoginput('nidaq','Dev1');
addchannel(ai,0);
```

You can return the coupling modes supported by the board with the Coupling field of the daghwinfo function.

```
out = daqhwinfo(ai);
out.Coupling
ans =
   'AC,DC'
```

Configure the channel contained by ai to use dc-coupling:

```
ai.Channel.Coupling = 'DC';
ai.Channel.Coupling
ans=
DC
```

ExternalClockDriveLine

Purpose

Specify which signal is driven by the clock indicating that an analog output update has occurred

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

ExternalClockDriveLine defines which pin is pulsed when analog output channels are updated. You can use this property to synchronize the operations of multiple cards over the RTSI bus or via external PFI pins.

Note The National Instruments term for this clock is AO Sample Clock.

Characteristics

Vendor National Instruments

Usage AO

Access Read/write

Data type String
Read-only when Yes

running

Values

PFI0 to PFI15 Use specified pin from PFI0 through PFI15.

RTSIO to RTSI6 Use specified pin from RTSIO through RTSI6.

See Also Properties

ExternalClockSource

Purpose

Specify which signal generates an analog output update across channels

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

ExternalClockSource specifies the pin whose signal is used as the clock to update analog outputs across a group of channels. This property is in effect when the ClockSource property is set to External.

Note The National Instruments term for this clock is AO Sample Clock.

Characteristics

Vendor National Instruments

Usage AO

Access Read/write

Data type String Read-only when Yes

running

Values PEIO to PEI15 Use specified pin from PFI0 through PFI15.

> RTSIO to RTSI6 Use specified pin from RTSIO through RTSI6.

See Also **Properties**

ClockSource

ExternalSampleClockDriveLine

Purpose

Specify which signal line is driven by the clock for sample conversions on each channel

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

ExternalSampleClockDriveLine defines which pin is pulsed when conversions occur on each channel. Data acquisition cards with simultaneous sample and hold ignore this property. You can use this property to synchronize the operations of multiple cards over the RTSI bus or via external PFI pins.

Note The National Instruments term for this clock is AI Convert Clock.

Characteristics

Vendor National Instruments

Usage AI

Access Read/write

Data type String

Read-only when running Yes

Values

PFI0 to PFI15 Use specified pin from PFI0 through PFI15.

RTSIO to RTSI6 Use specified pin from RTSIO through

RTSI6.

See Also Properties

ExternalSampleClockSource

ExternalSampleClockSource

Purpose

Specify which signal provides clock for sample conversions across channels

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

ExternalSampleClockSource specifies the pin whose signal is used as the channel clock for conversions on each channel. This property is in effect when the ClockSource property is set to ExternalSampleCtrl or ExternalSampleAndScanCtrl.

Data acquisition cards with simultaneous sample and hold ignore this property.

Note The National Instruments term for this clock is AI Convert Clock.

Characteristics

Vendor National Instruments

Usage AI

Access Read/write
Data type String

Read-only when running Yes

Values

PFI0 to PFI15 Use specified pin from PFI0 through PFI15.

RTSIO to RTSI6 Use specified pin from RTSIO through

RTSI6.

${\bf External Sample Clock Source}$

See Also Properties

ClockSource, ExternalScanClockSource

ExternalScanClockDriveLine

Purpose

Specify which signal is driven by the clock indicating the start of a series of conversions across channels

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

ExternalScanClockDriveLine defines which pin is pulsed when a series of conversions across channels start. You can use this property to synchronize the operations of multiple cards over the RTSI bus or via external PFI pins.

Note The National Instruments term for this clock is AI Sample Clock.

Characteristics

Vendor National Instruments

Usage AI

Access Read/write

Data type String

Read-only when

running

Values PFI0 to PFI15 Use specified pin from PFI0 through PFI15.

Yes

RTSI0 to RTSI6 Use specified pin from RTSI0 through RTSI6.

See Also Properties

ExternalScanClockSource

ExternalScanClockSource

Purpose

Specify which signal starts series of conversions across channels

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

ExternalScanClockSource specifies the pin whose signal is used as the scan clock to initiate conversions across a group of channels. This property is in effect when the ClockSource property is set to ExternalScanCtrl or ExternalSampleAndScanCtrl.

Note The National Instruments term for this clock is AI Sample Clock.

Characteristics

Vendor National Instruments

Usage AI

Access Read/write
Data type String

Read-only when

running

Yes

Values

PFI0 to PFI15 Use specified pin from PFI0 through PFI15.

RTSIO to RTSI6 Use specified pin from RTSIO through RTSI6.

See Also Properties

 ${\tt ClockSource}, {\tt ExternalSampleClockSource}$

ExternalTriggerDriveLine

Purpose

Specify which signal line is driven with a pulse when data acquisition or generation starts

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

ExternalTriggerDriveLine defines which pin is pulsed when a data acquisition or generation starts. You can use this property to synchronize the operations of multiple cards over the RTSI bus or via external PFI pins.

Characteristics

Vendor National Instruments

Usage AI

Access Read/Write

Data type String
Read-only when Yes

running

Values PFI0 to PFI15 Use specified pin from PFI0 through PFI15.

RTSI0 to RTSI6 Use specified pin from RTSI0 through RTSI6.

See Also Properties

HwDigitalTriggerSource

HwDigitalTriggerSource

Purpose

Specify which signal initiates data acquisition

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

HwDigitalTriggerSource defines which pin is used to initiate a data acquisition when the TriggerType property is set to HwDigital.

Characteristics

Vendor National Instruments

Usage AI, AO

Access Read/write

Data type String

Read-only when running Yes

Values

PFI0 to PFI15

Use specified pin from PFI0 through PFI15.

RTSI0 to RTSI6

Use specified pin from RTSI0 through

RTSI6.

See Also

Properties

 ${\tt TriggerType}$

Purpose

Specify number of external multiplexer devices connected

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

NumMuxBoards specifies the number of AMUX-64T multiplexer devices connected to your hardware. NumMuxBoards can be 0, 1, 2, or 4. If you are using a 1200 Series board, then NumMuxBoards can only be 0.

Characteristics

Vendor National Instruments Traditional NI-DAQ

devices

Usage AI, common to all channels

No

Access Read/write
Data type Double

Read-only when

running

Note The Traditional NI-DAQ adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition ToolboxTM object for Traditional NI-DAQ adaptor beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release. See the supported hardware page at www.mathworks.com/products/daq/supportedio.html for more information.

Values

{0}, 1, 2, or 4 The number of AMUX-64T multiplexer

devices connected.

OutOfDataMode

Purpose

Specify how value held by analog output subsystem is determined

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

When queued data is output to the analog output (AO) subsystem, the hardware typically holds a value. For National Instruments and Measurement Computing devices, the value held is determined by OutOfDataMode.

OutOfDataMode can be Hold or DefaultValue. If OutOfDataMode is Hold, then the last value output is held by the AO subsystem. If OutOfDataMode is DefaultValue, then the value specified by the DefaultChannelValue property is held by the AO subsystem.

Characteristics

Vendor Measurement Computing, National

Instruments

Usage AO, common to all channels

Access Read/write
Data type String

Read-only when

running

Yes

Values

{Hold the last output value.

DefaultValue Hold the value specified by

DefaultChannelValue.

Examples

Create the analog output object ao and add two channels to it.

```
ao = analogoutput('nidaq','Dev1');
addchannel(ao,0:1);
```

You can configure **ao** so that when queued data is finished being output, a value of 1 volt is held for both channels.

```
ao.OutOfDataMode = 'DefaultValue';
ao.Channel.DefaultChannelValue = 1.0;
```

See Also Properties

DefaultChannelValue

PortAddress

Purpose

Indicate base address of parallel port

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

The PC supports up to three parallel ports that are assigned the labels LPT1, LPT2, and LPT3. You can use any of these standard ports as long as they use the usual base addresses, which are (in hex) 378, 278, and 3BC, respectively.

Additional ports, or standard ports not assigned the usual base addresses, are not accessible by the toolbox. Note that most PCs that support MATLAB will include a single parallel printer port with base address 378 (LPT1).

Note The Parallel Port adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition ToolboxTM object for 'parallel' beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release. See the supported hardware page at www.mathworks.com/products/daq/supportedio.html for more information.

Characteristics

Vendor Parallel port

Usage DIO, common to all lines

Access Read only

Data type String

Read-only when running Yes

Values

The value is automatically defined when the object is created.

Examples

Create a digital I/O object for parallel port LPT1 and return the ${\tt PortAddress}$ value.

```
dio = digitalio('parallel','LPT1');
get(dio,'PortAddress')
ans =
0x378
```

The returned value indicates that LPT1 uses the usual base address.

StandardSampleRates

Purpose

Specify whether valid sample rates snap to small set of standard values, or if you can set sample rate to any allowed value

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

StandardSampleRates can be On of Off. If StandardSampleRates is Off, then it is possible to set the sample rate to any value within the bounds supported by the hardware. For most sound cards, the lower bound is 8.000 kHz, while the upper bound is 44.1 kHz. For newer sound cards, an upper bound of 96.0 kHz might be supported. The specified sample rate is rounded up to the next integer value.

If StandardSampleRates is On, then the available sample rates snap to a small set of standard values. The standard values are 8.000 kHz, 11.025 kHz, 22.050 kHz, and 44.100 kHz. If you specify a sampling rate that is within one percent of a standard value, then the sampling rate snaps to that standard value. If you specify a sampling rate that is not within one percent of a standard value, then the sampling rate rounds up to the closest standard value.

Regardless of the StandardSampleRates value, if you specify a sampling rate that is outside the allowed limits, then an error is returned.

Characteristics

Vendor Sound cards

AI, AO, common to all channels Usage

Read/write Access

Data type String Yes

Read-only when

running

${\bf Standard Sample Rates}$

Values	On	The sample rate can be set only to a small set of standard values.
	{Off}	If supported by the hardware, the sample rate can be set to any value within the allowed bounds, up to a maximum of 96.0 kHz.

TransferMode

Purpose

Specify how data is transferred from data acquisition device to system memory

Description

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

For National Instruments NI-DAQmx hardware, this property is ignored. The device driver automatically selects the most efficient transfer mode available.

For National Instruments Traditional NI-DAQ hardware, TransferMode can be Interrupts or SingleDMA for both analog input and analog output subsystems. If TransferMode is Interrupts, then data is transferred from the hardware first-in, first-out memory buffer (FIFO) to system memory using interrupts. If TransferMode is SingleDMA, then data is transferred from the hardware FIFO to system memory using a single direct memory access (DMA) channel. Some boards also support a TransferMode of DualDMA for analog input subsystems. For example, the AT-MIO-16E-1 board supports this transfer mode. If TransferMode is DualDMA, then data is transferred from the hardware FIFO to system memory using two DMA channels. Depending on your system resources, data transfer via interrupts can significantly degrade system performance.

For Measurement Computing hardware, TransferMode can be Default, InterruptPerPoint, DMA, InterruptPerBlock, or InterruptPerScan. If TransferMode is Default, the transfer mode is automatically selected by the driver based on the board type and the sampling rate. If TransferMode is InterruptPerPoint, a single conversion is transferred for each interrupt. You should use this property value if your sampling rate is less the 5 kHz or you specify a small block size for memory buffering (as defined by the BufferingConfig property). If TransferMode is DMA, data is transferred using a single DMA channel. If TransferMode is InterruptPerBlock, a block of data is transferred for each interrupt. You should use this property value if your sampling rate is greater than 5 kHz and you are using a board that has a fast

maximum sampling rate. Note that a data block is defined by the board, and usually corresponds to half the FIFO size. If TransferMode is InterruptPerScan, data is not transferred until the entire scan is complete. This can only be used when the number of points acquired is less than or equal to the FIFO size. You should use this mode if your sampling rate is higher than the maximum continuous scan rate of the data acquisition device.

Note If your sampling rate is greater than $\sim 5~\mathrm{kHz}$, you should avoid using interrupts if possible. The recommended TransferMode setting for your application will be described in your hardware documentation, and depends on the specific board you are using and your platform configuration.

Characteristics

Vendor Measurement Computing, National

Instruments

Usage AI, AO, common to all channels

Access Read/write

Data type String

Read-only when Yes

running

Values Advantech

{InterruptPerPoint} Transfer single data points using interrupts.

InterruptPerBlock Transfer a block of data using interrupts (AI

only).

TransferMode

Measurement Computing

{Default} The transfer mode is automatically selected

by the driver based on the board type and the

sampling rate.

InterruptPerPoint Transfer single data points using interrupts.

DMA Transfer data using a single DMA channel

(AI only).

InterruptPerBlock Transfer a block of data using interrupts (AI

only).

InterruptPerScan Transfer all data when the acquisition is

complete (AI only).

National Instruments

Interrupts Transfer data using interrupts.

SingleDMA Transfer data using a single DMA channel.

DualDMA Transfer data using two DMA channels.

This default property value is supplied by the driver. For most devices that support data transfer via interrupts and DMA, SingleDMA is the default value.

Note The Traditional NI-DAQ adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition ToolboxTM object for Traditional NI-DAQ adaptor beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release. See the supported hardware page at www.mathworks.com/products/daq/supportedio.html for more information.

Examples

Set the $\mbox{TransferMode}$ property for a National Instruments board before acquiring data.

```
ai = analoginput('nidaq', 1);
set(ai, 'TransferMode', 'SingleDMA');
addchannel(ai, 1:2);
softscope(ai)
```

TransferMode

Block Reference

Analog Input
Analog Input (Single Sample)
Analog Output
Analog Output (Single Sample)
Digital Input
Digital Output

Analog Input

Purpose

Acquire data from multiple analog channels of data acquisition device

Library

Data Acquisition Toolbox

Note You cannot use certain devices with Data Acquisition Toolbox Simulink® blocks. Refer to the Supported Hardware page to see if your device supports Simulink use.

Description

The Analog Input block opens, initializes, configures, and controls an analog data acquisition device. The opening, initialization, and configuration of the device occur once at the start of the model's execution. During the model's run time, the block acquires data either synchronously (deliver the current block of data the device is providing) or asynchronously (buffer incoming data).

Note You need a license for both Data Acquisition Toolbox and Simulink software to use this block.

The block has no input ports. It has one or more output ports, depending on the configuration you choose in the Source Block Parameters dialog box. The following diagram shows the block configured with one port for both channels and with one port for each channel, in the case of a device that has two channels.

winsound 0 NVIDIA(R) n... 8000 samples/sec

Analog Input - one port

winsound 0 Left : NVIDIA(R) n... 8000 samples/sec Right :

Analog Input - two ports

Use the Analog Input block to incorporate live measured data into Simulink for:

- System characterization
- Algorithm verification
- System and algorithm modeling
- Model and design validation
- design control

Note You can use the Analog Input block only with devices that support clocked acquisition. The block will error out when the model is run with a device that does not support clocking. To acquire data using devices that do not support clocking, use the Analog Input (Single Sample) block.


You can use this block for signal applications by using it with basic Simulink and DSP System ToolboxTM.

You can use the Analog Input block either synchronously or asynchronously. Select the acquisition mode in the Source Block Parameters dialog box.

The following diagram shows the basic analog input usage scenario, in which you would:

- Acquire data at each time step or once per model execution.
- Analyze the data, or use it as input to a system in the model.
- Optionally display results.

Analog Input

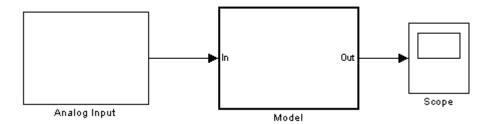

For an example of creating a model using the Analog Input block, see Example: Bringing Analog Data into a Model.

Other Supported Features

The Analog Input block supports the use of Simulink Accelerator™ mode. This feature speeds up the execution of Simulink models.


Note You need the C++ Compiler to use Simulink Accelerator mode.

The block supports the use of model referencing. This feature lets your model include other Simulink models as modular components.

For more information on these features, see the Simulink documentation.

Dialog Box

Use the Source Block Parameters dialog box to select your acquisition mode and to set other configuration options.

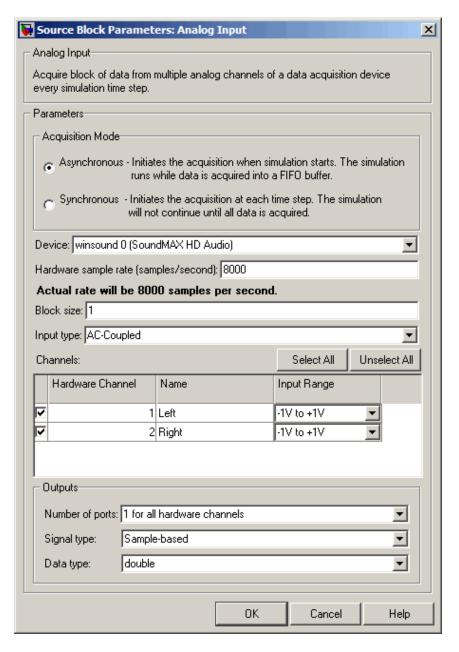

Acquisition Mode

Asynchronous


Initiates the acquisition when the simulation starts. The simulation runs while data is acquired into a FIFO (First in, First out) buffer. The acquisition is continuous; the block buffers data while outputting a scan/frame of data at each time step.

Synchronous

Initiates the acquisition at each time step. The simulation will not continue until the requested block of data is acquired. This is unbuffered input; the block will synchronously output the latest scan/frame of data at each time step.


The following diagrams show the difference between synchronous and asynchronous modes for the Analog Input block.

Synchronous Analog Input

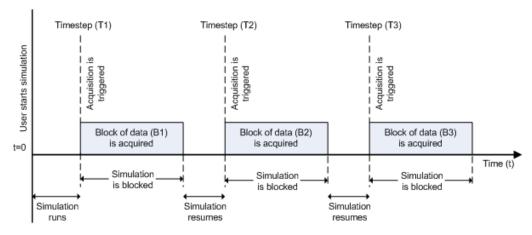


At the first time step (T1), the acquisition is initiated for the required block of data (B1). The simulation does not continue until B1 is completely acquired.

Asynchronous Analog Input - Scenario 1

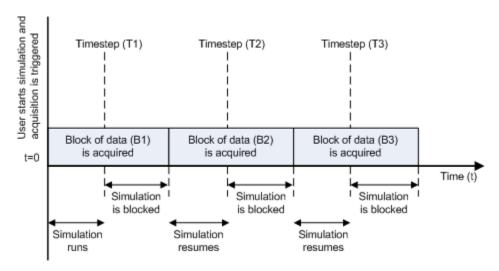

Scenario 1 shows the case when simulation speed outpaces data acquisition speed. At the first time step (T1), the required block of data (B1) is still being acquired. Therefore, the simulation does not continue until B1 is completely acquired.

Analog Input

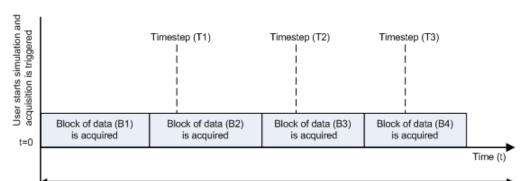

Simulation runs continuously

Asynchronous Analog Input - Scenario 2

Scenario 2 shows the case when data acquisition speed outpaces simulation speed. At the first time step (T1), the required block of data (B1) has been completely acquired. Therefore, the simulation runs continuously.

Note Several factors, including device hardware and model complexity, can affect the simulation speed, causing both scenarios 1 and 2 to occur within the same simulation.

Options

Device

The data acquisition device from which you want to acquire data. The items in the list vary, depending on which devices you have connected to your system. Devices in the list are specified by adaptor/vendor name and unique device ID, followed by the name of the device. The first available device is selected by default.

Hardware sample rate

The rate at which samples are acquired from the device, in samples per second. This is the sampling time for the hardware. The default is defined when a device is selected.

The sample rate must be a positive real number, and be within the range allowed for the selected hardware.

Block size

The desired number of data samples to output at each time step for each channel. Block size corresponds to the SamplesPerTrigger property for an analog input device. The default value for block size depends on the hardware selected. It must be a positive integer, and be within the range allowed for the selected hardware.

Input type

Specifies the hardware channel configuration, such as single-ended, differential, etc. The input type is defined by the capabilities of the selected device.

Channels

The channel configuration table lists your device's hardware channels and lets you configure them. Use the check boxes and selection buttons to specify which channels to acquire data from. These parameters are specified for each selected channel:

Hardware Channel — Displays the hardware channel ID specified by the device. The **Hardware Channel** column is read only and the parameters are defined when the device is selected.

The **Name** — Specifies the channel name. By default the table displays any names provided by the hardware, but you can edit the names. For example, if the device is a sound card with two channels, you can name them Left and Right.

Analog Input

Input Range — Specifies the input ranges available for each channel supported by the hardware, and is defined when a device is selected.

Outputs

Number of ports

Select 1 for all hardware channels (default) or 1 per hardware channel.

Using **1** for all hardware channels outputs data from a single port as a matrix, with a size of Block size x Number of Channels selected.

Using 1 per hardware channel outputs data from N ports, where N is equal to the number of selected channels. Each output port will be a column vector with a size of Block size x 1. For naming, each output port will use the channel name if one was specified, or otherwise use [HWChannel + channel ID], for example, HWChannel2.

Signal type

Select **Sample-based** or **Frame-based**. This option determines whether the signal type is sample-based or frame-based. **Sample-based** is the default.

Note The **Frame-based** option works only if you have the DSP System Toolbox software installed.

Data type

Select your data type to output from the block. The Analog Input block supports double and native data types, as supported by the hardware. double is the default. Native data types will be dynamically populated in this list based on the hardware that is selected. For example, if int16 is a native data type of a specific

Analog Input

hardware device, then one of the entries for ${\bf Data}\ {\bf type}$ will be int16 (native).

See Also

Analog Input (Single Sample), Analog Output, Analog Output (Single Sample), Digital Input, Digital Output

Purpose

Acquire single sample from multiple analog channels of data acquisition device

Library

Data Acquisition Toolbox

Note You cannot use certain devices with Data Acquisition Toolbox Simulink blocks. Refer to the Supported Hardware page to see if your device supports Simulink use.

Description


The Analog Input (Single Sample) block opens, initializes, configures, and controls an analog data acquisition device. The opening, initialization, and configuration of the device occur once at the start of the model's execution. The block acquires a single sample every sample time, synchronously from the device, during the model's run time.

Note You need a license for both Data Acquisition Toolbox and Simulink software to use this block.

The block has no input ports. It has one or more output ports, depending on the configuration you choose in the Source Block Parameters dialog box. The following diagram shows the block configured with one port for both channels and with one port for each channel, in the case of a device that has two channels.

mee 1 PMD-1208FS Analog Input

(Single Sample) - one port

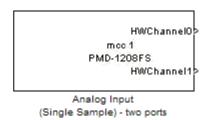

Use the Analog Input (Single Sample) block to incorporate live measured data into Simulink for:

- System characterization
- Algorithm verification
- System and algorithm modeling
- Model and design validation
- Controls design

Note You can use Analog Input (Single Sample) block only with devices that support single sample acquisition. The block will error out when the model is run with a device that does not support single sample acquisition. To acquire data from devices that do not support acquisition of a single sample (like devices designed for sound and vibration), use the Analog Input block.

You can use the Analog Input (Single Sample) block for signal applications by using it with basic Simulink and DSP System Toolbox.

Other Supported Features

The Analog Input (Single Sample) block supports the use of Simulink Accelerator mode. This feature speeds up the execution of Simulink models.


Note You need the C++ Compiler to use Simulink Accelerator mode.

This block supports the use of model referencing. This feature lets your model include other Simulink models as modular components.

For more information on these features, see the Simulink documentation.

Dialog Box

Use the Source Block Parameters dialog box to select your device and to set other configuration options.

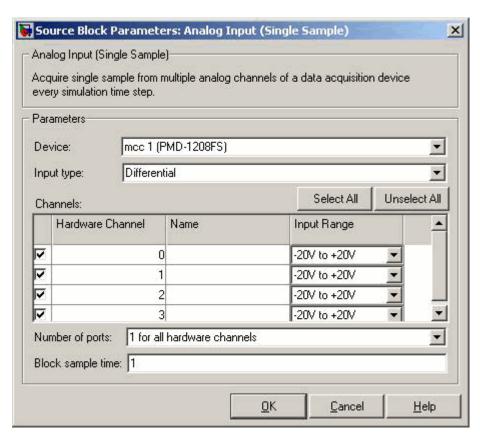

Device

The data acquisition device from which you want to acquire data. The items in the list vary, depending on which devices you have connected to your system. Devices in the list are specified by adaptor/vendor name and unique device ID, followed by the name of the device. The first available device is selected by default.

Input type

Specifies the hardware channel configuration, such as single-ended, differential, etc. When you select a device, the device capability defines the available values for input type.

Channels

The channel configuration table lists your device's hardware channels and lets you configure them. Use the check boxes and selection buttons to specify which channels to acquire data from. These parameters are specified for each selected channel:

Hardware Channel — Displays the hardware channel ID specified by the device. The **Hardware Channel** column is read-only and the parameters are defined when the device is selected.

Name — Specifies the channel name. By default the table will display any names provided by the hardware, but you can edit the names. For example, if you are using a device to acquire indoor and outdoor temperature from two channels, you can name them IndoorTemp and OutdoorTemp.

Input Range — Specifies the input ranges available for each channel supported by the hardware, and the available values are defined when a device is selected.

Number of ports

Select 1 for all hardware channels (default) or 1 per hardware channel.

Using **1** for all hardware channels, outputs data from a single port as a matrix, with a size of [1 x Number of Channels selected].

Using 1 per hardware channel, outputs data from N ports, where N is equal to the number of selected channels. Each output port will be a scalar value. For naming, each output port will use the channel name if one was specified, or otherwise use ["HWChannel" + channel ID], for example, HWChannel2.

Block sample time

Specifies the sample time of the block during the simulation. This is the rate at which the block is executed during simulation. The default value is 0.01 (seconds).

See Also

Analog Input, Analog Output, Analog Output (Single Sample), Digital Input, Digital Output

Purpose

Output data to multiple analog channels of data acquisition device

Library

Data Acquisition Toolbox

Note You cannot use certain devices with Data Acquisition Toolbox Simulink blocks. Refer to the Supported Hardware page to see if your device supports Simulink use.

Description


The Analog Output block opens, initializes, configures, and controls an analog data acquisition device. The opening, initialization, and configuration of the device occur once at the start of the model's execution. During the model's run time, the block outputs data to the hardware either synchronously (outputs the block of data as it is provided) or asynchronously (buffers output data).

Note You need a license for both Data Acquisition Toolbox and Simulink software to use this block.

The block has one or more input ports, depending on the option you choose in the Sink Block Parameters dialog box. It has no output ports. The following diagram shows the block configured with one port for both channels and with one port for each channel, in the case of a device that has two channels selected.

Analog Output

winsound 0 > NVIDIA(R) n... 8000 samples/sec

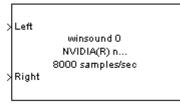

Analog Output - one port Analog Output - two ports

Note You can use the Analog Output block only with devices that support clocked generation. The block will error out when the model is run with a device that does not support clocking. To send data using devices that do not support clocking, use the Analog Output (Single Sample) block.

The Analog Output block inherits the sample time from the driving block connected to the input port. The valid data types of the signal at the input port are double or native data types supported by the hardware.

Other Supported Features

The Analog Output block supports the use of Simulink Accelerator mode. This feature speeds up the execution of Simulink models.


Note You need the C++ Compiler to use Simulink Accelerator mode.

The block supports the use of model referencing. This feature lets your model include other Simulink models as modular components.

For more information on these features, see the Simulink documentation.

Dialog Box

Use the Sink Block Parameters dialog box to select your acquisition mode and to set other configuration options.

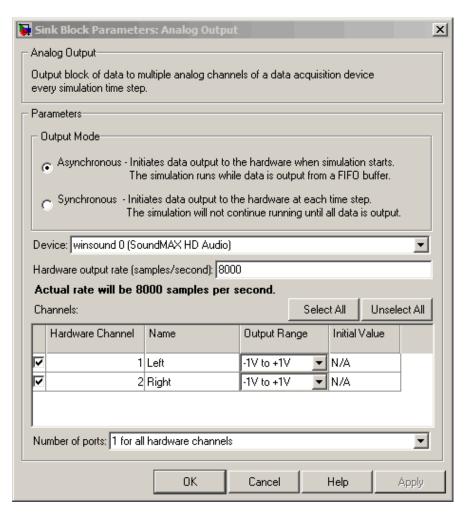

Analog Output

Output Mode

Asynchronous


Initiates data output to the hardware when simulation starts. The simulation runs while data is output from a FIFO (First In, First Out) buffer. This mode buffers and outputs data from the block, letting you perform a frame-based or sample-based output.

Synchronous

Initiates data output to the hardware at each time step. The simulation will not continue running until the current block of data is output. In synchronous mode, the block synchronously outputs a vector or frame of samples provided at each time step.


The following diagrams show the difference between synchronous and asynchronous analog output.

Synchronous Analog Output

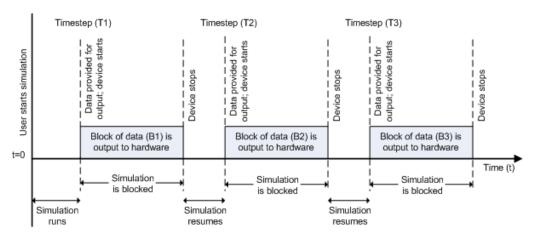


At the first time step (T1), data output is initiated and the corresponding block of data (B1) is output to the hardware. The simulation does not continue until B1 is output completely.

Asynchronous Analog Output - Scenario 1

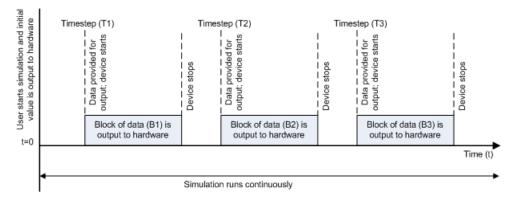

Scenario 1 shows the case when data output speed outpaces simulation speed. At the first time step (T1), data output is initiated and the corresponding block of data (B1) is output to the hardware. The simulation runs continuously in this mode.

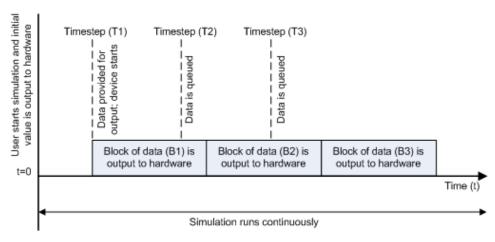

Asynchronous Analog Output - Scenario 2

Scenario 2 shows the case when simulation speed outpaces data acquisition speed. At the first time step (T1), data output is initiated and the corresponding block of data (B1) is output to the hardware. Data is queued at successive time steps and is output to the hardware once the previous block completes. The simulation runs continuously in this mode.

Note Several factors, including device hardware and model complexity, can affect the simulation speed, causing both scenarios 1 and 2 to occur within the same simulation.

Options

Device

The data acquisition device to which you want to output data. The items in the list vary, depending on which devices you have connected to your system. Devices in the list are specified by adaptor/vendor name and unique device ID, followed by the name of the device. The first available device is selected by default.

Hardware output rate

The rate at which samples are output to the device, in samples per second. This output rate for the hardware is defined when a device is selected. The output rate specified must be within the range supported by the selected device.

Channels

The channel configuration table lists your device's hardware channels and lets you configure them. Use the check boxes and selection buttons to specify which channels to send data to.

Hardware Channel — Displays the channel ID specified by the device, and is read only.

Name — specifies the channel name. By default the table displays any names provided by the hardware, but you can edit the names. For example, if the device is a sound card with two channels, you can name them Left and Right.

Output Range — Specifies the output ranges available for each channel supported by the hardware, and is defined by the selected device.

Initial Value — Specifies the initial value to be output at the start of the simulation, if you are using Asynchronous mode. The default value is 0. In Synchronous mode, the **Initial Value** column does not appear in the table.

Note For AC-coupled devices like a sound card, this column is not used and is read only.

Number of ports

Select 1 for all hardware channels (default) or 1 per hardware channel.

Analog Output

Using 1 for all hardware channels inputs data from a single port as a matrix, with a size of [S x Number of Channels selected], where S is number of samples provided as input.

Using 1 per hardware channel inputs data from N ports, where N is equal to the number of selected channels. Each input port will be a column vector with a size of [S x 1], where S is the number of samples provided as an input. For naming, each output port will use the channel name if one was specified, or otherwise use ["HWChannel" + channel ID], for example, HWChannel2.

See Also

Analog Input, Analog Input (Single Sample), Analog Output (Single Sample), Digital Input, Digital Output

Purpose

Output single sample to multiple analog channels of data acquisition device

Library

Data Acquisition Toolbox

Note You cannot use certain devices with Data Acquisition Toolbox Simulink blocks. Refer to the Supported Hardware page to see if your device supports Simulink use.

Description

The Analog Output (Single Sample) block opens, initializes, configures, and controls an analog data acquisition device. The opening, initialization, and configuration of the device occur once at the start of the model's execution. The block outputs a single sample every sample time, synchronously to the hardware, during the model's run time.

Note You need a license for both Data Acquisition Toolbox and Simulink software to use this block.

The block has one or more input ports, depending on the option you choose in the Sink Block Parameters dialog box. It has no output ports. The following diagram shows the block configured with one port for both channels and with one port for each channel, in the case of a device that has two channels selected.



> HWChannel0 nidaq Dev1 PCI-6251 > HWChannel1

(Single Sample) - two ports

Note You can use Analog Output (Single Sample) block only with devices that support single sample output. The block will error out when the model is run with a device that does not support single sample acquisition. To send data using devices that do not support acquisition of a single sample (like devices designed for sound and vibration), use the Analog Output block.

The Analog Output (Single Sample) block inherits the sample time from the driving block connected to the input port. The valid data type of the signal at the input port is double.

Other Supported Features

The Analog Output (Single Sample) block supports the use of Simulink Accelerator mode. This feature speeds up the execution of Simulink models.


Note You need the C++ Compiler to use Simulink Accelerator mode.

The Analog Output (Single Sample) block supports the use of model referencing. This feature lets your model include other Simulink models as modular components.

For more information on these features, see the Simulink documentation.

Dialog Box

Use the Sink Block Parameters dialog box to select your device and to set other configuration options.

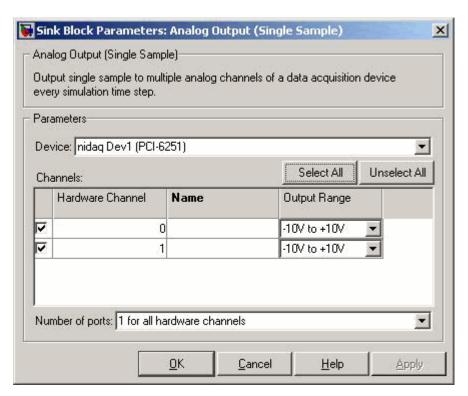

Device

The data acquisition device to which you want to output data. The items in the list vary, depending on which devices you have connected to your system. Devices in the list are specified by adaptor/vendor name and unique device ID, followed by the name of the device. The first available device is selected by default.

Channels

The channel configuration table lists your device's hardware channels and lets you configure them. Use the check boxes and

selection buttons to specify which channels to acquire data from. These parameters are specified for each selected channel:

Hardware Channel — Displays the hardware channel ID specified by the device. The **Hardware Channel** column is read-only and the parameters are defined when the device is selected.

Name — Specifies the channel name. By default the table will display any names provided by the hardware, but you can edit the names. For example, if you are sending data and trigger signals to an output device, you can name them Data and TriggerStatus.

Output Range — Specifies the output ranges available for each channel supported by the hardware, and the available values are defined when a device is selected.

Number of ports

Select 1 for all hardware channels (default) or 1 per hardware channel.

Using **1** for all hardware channels, receives data from a single port as a matrix, with a size of [Block size x Number of Channels selected].

Using 1 per hardware channel, receives data from N ports, where N is equal to the number of selected channels. Each input port will be a scalar. For naming, each output port will use the channel name if one was specified, or otherwise use ["HWChannel" + channel ID], for example, HWChannel2.

See Also

Analog Input, Analog Input (Single Sample), Analog Output, Digital Input, Digital Output

Purpose

Acquire latest set of values from multiple digital lines of data

acquisition device

Library

Data Acquisition Toolbox

Note You cannot use certain devices with Data Acquisition Toolbox Simulink blocks. Refer to the Supported Hardware page to see if your device supports Simulink use.

Description

The Digital Input block synchronously outputs the latest scan of data available from the digital lines selected at each simulation time step. It acquires unbuffered digital data, and the data delivered is a binary vector.


Note You need a license for both Data Acquisition Toolbox and Simulink software to use this block.

The block has no input ports. It has one or more output ports, depending on the option you choose in the Source Block Parameters dialog box. The following diagram shows the block configured with one port for all lines and with one port for each line, in the case of a device that has 17 lines selected.

Digital Input




Digital Input - one port for all lines

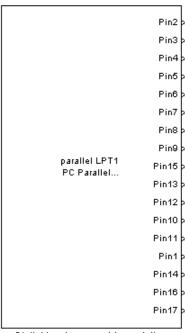

Digital Input - one port for each line

The block inherits the sample time of the model.

The output data is always a binary vector (binvec), i.e., a vector of logical values.

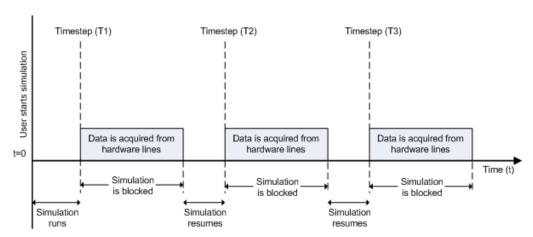

Digital input acquisition is done synchronously. The following diagram shows synchronous digital input.

At the first time step (T1), data is acquired from the selected hardware lines. The simulation does not continue until data is read from all lines.

Other Supported Features

The Digital Input block supports the use of Simulink Accelerator mode. This feature speeds up the execution of Simulink models.


Note You need the C++ Compiler to use Simulink Accelerator mode.

The block supports the use of model referencing. This feature lets your model include other Simulink models as modular components.

For more information on these features, see the Simulink documentation.

Dialog Box

Use the Source Block Parameters dialog box to set configuration options.

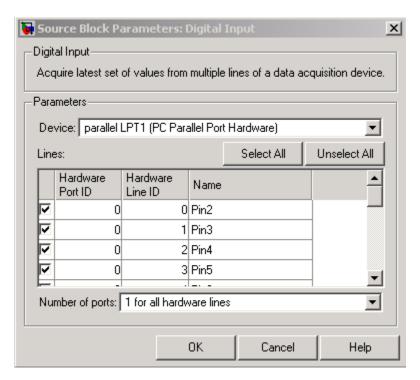

Device

The data acquisition device from which you want to acquire data. The items in the list vary, depending on which devices you have connected to your system. Devices in the list are specified by adaptor/vendor name and unique device ID, followed by the name of the device. The first available device is selected by default.

Lines

The line configuration table lists your device's lines and lets you configure them. The table lists all the lines that can be configured for input. Use the check boxes and selection buttons to specify which lines to acquire data from.

Hardware Port ID

Specifies the ID for each hardware port. This is automatically detected and filled in by the selected device, and is read only.

Hardware Line ID

Specifies the ID of the hardware line. This is automatically detected and filled in by the selected device, and is read only.

Name

Specifies the hardware line name. This is automatically detected and filled in from the hardware, though you can edit the name.

Number of ports

Select 1 for all hardware lines (default) or 1 per hardware line.

Using 1 for all hardware lines means that the block will have only one output port for all of the lines that are selected in the table. Data must be [S x number of lines], where S is the number of samples. Data will be a binary vector (binvec).

Using 1 per hardware line means the block will have one output port per selected line. The name of each output port is the name specified in the table for each line. If no name is provided, the name is "Port" + HwPort ID + "Line" + Line ID. For example, if line 2 of hardware port 3 is selected, and you did not specify a name in the line table, Port3Line2 appears in the block. Data will be [1 x 1].

See Also

Analog Input, Analog Input (Single Sample), Analog Output, Analog Output (Single Sample), Digital Output

Purpose

Output data to multiple digital lines of data acquisition device

Library

Data Acquisition Toolbox


Note You cannot use certain devices with Data Acquisition Toolbox Simulink blocks. Refer to the Supported Hardware page to see if your device supports Simulink use.

Description


The Digital Output block synchronously outputs the latest set of data to the hardware at each simulation time step. It outputs unbuffered digital data. The output data is always a binary vector (binvec).

Note You need a license for both Data Acquisition Toolbox and Simulink software to use this block.

The block has no output ports. It can have one or more input ports, depending on the option you choose in the Sink Block Parameters dialog box. The following diagram shows the block configured with one port for all lines and with one port for each line, in the case of a device that has 12 lines selected.

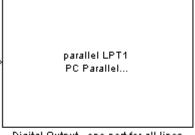


Digital Output - one port for all lines

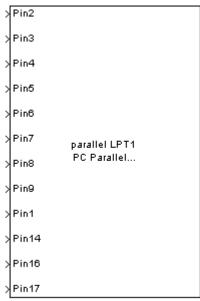

Digital Output - one port for each line

The Digital Output block inherits the sample time from the driving block connected to the input port. The data type of the signal at the input port must be a logical data type.

Digital output is done synchronously. The following diagram shows synchronous digital output.

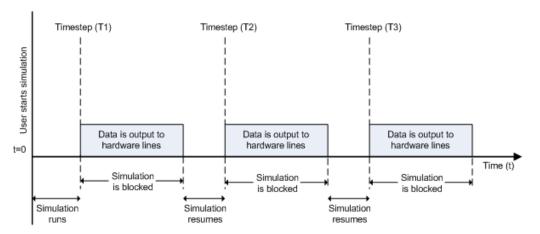

At the first time step (T1), data is output to the selected hardware lines. The simulation does not continue until data is output to all lines.

Other Supported Features

The Digital Output block supports the use of Simulink Accelerator mode. This feature speeds up the execution of Simulink models.


Note You need the C++ Compiler to use Simulink Accelerator mode.

The block supports the use of model referencing. This feature lets your model include other Simulink models as modular components.

For more information on these features, see the Simulink documentation.

Dialog Box

Use the Sink Block Parameters dialog box to set configuration options.

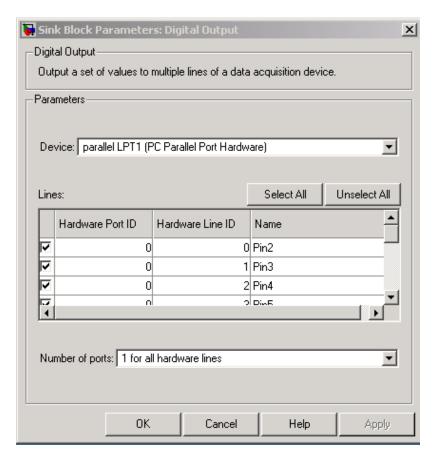

Device

The data acquisition device to which you want to output data. The items in the list vary, depending on which devices you have connected to your system. Devices in the list are specified by adaptor/vendor name and unique device ID, followed by the name of the device. The first available device is selected by default.

Lines

The line configuration table lists your device's lines and lets you configure them. Use the check boxes and selection buttons to specify which lines to send data to.

Hardware Port ID

Specifies the ID for each hardware port. This is automatically detected and filled in by the selected device, and is read only.

Hardware Line ID

Specifies the ID of the hardware line. This is automatically detected and filled in by the selected device, and is read only.

Name

Specifies the hardware line name. This is automatically detected and filled in by the selected device, though you can edit the name.

Number of ports

Select 1 for all hardware lines (default) or 1 per hardware line.

Using 1 for all hardware lines means that the block will have only one input port for all lines selected in the table. Data needs to be [S x number of lines], where S is the number of samples. Data at the input port needs to be a binary vector (binvec).

Using 1 per hardware line means the block will have one input port per selected line. The name of each input port is the name specified in the table for each line. If no name is provided, the name is "Port" + HwPort ID + Line + Line ID. For example, if line 2 of port 3 is selected, and you did not specify a name in the line table, Port3Line2 appears in the block. Data needs to be [1 x 1].

See Also

Analog Input, Analog Input (Single Sample), Analog Output, Analog Output (Single Sample), Digital Input

Function Reference

Device Object Creation (p. 6-2) Create device objects

Channels and Lines (p. 6-3) Add hardware channels and lines

Property Getting and Setting (p. 6-4) Display properties

Object Execution (p. 6-5) Execute device objects

Data Manipulation (p. 6-6) Work with acquired data

Information and Help (p. 6-7) Display help information

Utilities (p. 6-8) General Data Acquisition Toolbox

functionality

Session-Based Interface (p. 6-10) Session-based interface functions

and methods using the session-based

interface



Device Object Creation

analoginput an a log outputdigitalio

Create analog input object Create analog output object Create digital I/O object

Channels and Lines

addchannel Add hardware channels to analog

input or output object

addline Add hardware lines to digital I/O

object

addmuxchannel Add hardware channels to analog

input objects when using National Instruments multiplexer board



Property Getting and Setting

Device object properties get

Open Property Inspector inspect

Configure or display device object set

properties

setverify Configure and return specified

property

Object Execution

islogging Determine whether analog input

object is logging data

isrunning Determine whether device object is

running

issending Determine whether analog output

object is sending data

start Start device object

stop Stop device object

trigger Manually execute trigger for analog

input or output object

wait Wait until analog input or output

device object stops running



Data Manipulation

flushdata Remove analog input data from data

acquisition engine

getdata Extract analog input data, time,

and event information from data

acquisition engine

Immediately acquire one analog getsample

input sample

getvalue Read values from digital input and

output lines

peekdata Preview most recent acquired analog

input data

putdata Queue analog output data in engine

for eventual output

putsample Immediately output one analog

output sample

putvalue Write values to digital output lines

Information and Help

daqhelp Help for device objects, constructors,

adaptors, functions, and properties

daqhwinfo Data acquisition hardware

information

propinfo Property characteristics for device

objects, channels, or lines



Utilities

binvec2dec Convert digital input and output

binary vector to decimal value

clear Remove device objects from

MATLAB workspace

dagcallback Callback function that displays

event information for specified event

dagfind Return device objects, channels, or

lines from data acquisition engine to

MATLAB workspace

daqmem Allocate or display analog input and

output memory resources

dagread Read Data Acquisition Toolbox (.daq)

file for analog input

Register or unregister hardware dagregister

driver adaptor

dagreset Remove device objects, engine

MEX-file, and adaptor DLLs from

memory

dec2binvec Convert digital input and output

decimal value to binary vector

delete Remove device objects, channels, or

lines from data acquisition engine

disp Summary information for device

objects, channels, or lines

ischannel Check for channels

isdioline Check for lines

isvalid Determine whether device objects,

channels, or lines are valid

length Length of device object, channel

group, or line group

load Load device objects, channels, or

lines into MATLAB workspace

makenames List descriptive channel or line

names

muxchanidx Multiplexed scanned analog input

channel index

obj2mfile Convert device objects, channels, or

lines to MATLAB code

save Save device objects to MAT-file

showdaqevents Analog input and output event log

information

size Size of device object, channel group,

or line group

Session-Based Interface

addAnalogInputChannel

(daq.Session)

Add analog input channel

addAnalogOutputChannel

(daq.Session)

Add analog output channel

addClockConnection (daq.Session) Add clock connection

addCounterInputChannel

(dag.Session)

Add counter input channel

addCounterOutputChannel

(daq.Session)

Add counter output channel

addDigitalChannel (daq.Session) Add digital channel

Create event listener addlistener (dag.Session)

addTriggerConnection (daq.Session) Add trigger connection

binaryVectorToDecimal Convert binary vector value to

decimal

binaryVectorToHex Convert binary vector value to

hexadecimal

daq.createSession Create data acquisition session for

specific vendor hardware

daq.getDevices Display available National

Instruments devices

daq.getVendors Display available vendors

DataAvailable Event Notify when acquired data is

available to process

DataRequired Event Notify when additional data is

required for output on continuous

generation

decimalToBinaryVector Convert decimal value to binary

ErrorOccurred Event Notify when device-related errors

occur

hexToBinaryVector Convert hexadecimal value to binary

vector

inputSingleScan (daq.Session) Acquire single scan from all input

channels

outputSingleScan (daq.Session) Generate single scan on all output

channels

prepare (daq.Session) Prepare session for operation

queueOutputData (daq.Session) Queue data to be output release (daq.Session) Release session resources

removeChannel (daq.Session)

Remove channel from session object
removeConnection (daq.Session)

Remove clock or trigger connection

resetCounters (daq.Session) Reset counter channel to initial

count

startBackground (daq.Session) Start background operations

startForeground (daq.Session) Start foreground operations stop (daq.Session) Stop background operation

wait (daq.Session) Block MATLAB until background

operation completes

Class Reference

daq.Session

Represent data acquisition session using National Instruments devices

Purpose Represent data acquisition session using National Instruments devices

Description The session object configures and controls one or more devices including

devices plugged into a CompactDAQ chassis, using the session-based

interface. This class is not instantiated directly.

Construction

s = daq.createSession('vendor') creates the data acquisition session s to work with vendor devices. Currently the only supported vendor is National Instruments.

Input Arguments

vendor

Is the ID of the device vendor you want to use. Currently the only supported vendor is 'ni'.

Properties ActiveEdge Rising or falling edges of

EdgeCount signals

Active Pulse Active pulse measurement of

PulseWidth counter channel

ADCTimingMode Set channel timing mode

AutoSyncDSA Automatically Synchronize DSA

devices

BridgeMode Specify analog input device bridge

mode

Channels Array of channel objects

associated with session object

Connections Array of connections in a session

CountDirection Specify direction of counter

channel

Coupling Specify input coupling mode

Destination Indicates trigger destination

terminal

Device annel device

Direction Specify digital channel direction
DurationInSeconds Specify duration of acquisition
DutyCycle Duty cycle of counter output

channel

EncoderType Encoding type of counter channel

ExcitationCurrent Voltage of external source of

excitation

ExcitationSource External source of excitation

ExcitationVoltage Voltage of excitation source

External Trigger Timeout Indicate if external trigger timed

out

Frequency Frequency of generated pulses on

counter output channel

ID ID

IdleState Default state of counter output

channel

InitialCount Specify initial count point
InitialDelay Delay until output channel

generates pulses

IsContinuous Specify if operation continues

until manually stopped

IsDone Indicate if operation is complete
IsLogging Indicate if hardware is acquiring

or generating data

IsNotifyWhenDataAvailableExceedsCourtcrol if is set automatically

IsNotifyWhenScansQueuedBelowAutontrol if is set automatically

IsRunning Indicate if operation is still in

progress

IsSimulated Indicate if device is simulated

IsWaitingForExternalTrigger Indicates if synchronization is

waiting for an external trigger

MaxSoundPressureLevel Sound pressure level for

microphone channels

MeasurementType Type counter channel

measurement

Name Specify descriptive name for the

channel

Name Specify descriptive name for the

channel

NominalBridgeResistance Resistance of sensor

NotifyWhenDataAvailableExceeds Control firing of DataAvailable

event

NotifyWhenScansQueuedBelow Control firing of DataRequired

event

Number Of Scans for operation

when starting

R0 Specify resistance value

Range Specify channel measurement

range

Rate of operation in scans per

second

RateLimit Limit of rate of operation based

on hardware configuration

RTDConfiguration Specify wiring configuration of

RTD device

RTDType Specify sensor sensitivity

ScansAcquired Number of scans acquired during

operation

ScansOutputByHardware Indicate number of scans output

by hardware

ScansQueued Indicate number of scans queued

for output

Sensitivity Sensitivity of an channel

ShuntLocation Indicate location of channel's

shunt resistor

ShuntResistance Resistance value of channel's

shunt resistor

Source Indicates trigger source terminal

Terminal PFI terminal

TerminalConfig Specify terminal configuration

Terminals Terminals available on device or

CompactDAQ chassis

Thermocouple Type Select thermocouple type

TriggerCondition Specify condition that must be

satisfied before trigger executes

TriggersPerRun Indicate the number of times the

trigger executes in an operation

Triggers Remaining Indicates the number of trigger to

execute in an operation

TriggerType Type of trigger executed

Type Display synchronization trigger

type

Units Specify unit of RTD measurement

Vendor Vendor information associated

with session object

ZResetCondition Reset condition for Z-indexing
ZResetEnable Enable reset for Z-indexing
ZResetValue Reset value for Z-indexing

Methods

addAnalogInputChannel Add analog input channel addAnalogOutputChannel Add analog output channel

addClockConnection Add clock connection

addCounterInputChannel Add counter input channel addCounterOutputChannel Add counter output channel

addDigitalChannel Add digital channel
addlistener Create event listener
addTriggerConnection Add trigger connection

binaryVectorToDecimal Convert binary vector value to

decimal

binaryVectorToHex Convert binary vector value to

hexadecimal

DataAvailable Event Notify when acquired data is

available to process

DataRequired Event Notify when additional data is

required for output on continuous

generation

decimalToBinaryVector Convert decimal value to binary

vector

ErrorOccurred Event Notify when device-related errors

occur

hexToBinaryVector Convert hexadecimal value to

binary vector

inputSingleScan Acquire single scan from all input

channels

outputSingleScan Generate single scan on all output

channels

prepare Prepare session for operation

queueOutputData Queue data to be output
release Release session resources

removeChannel Remove channel from session

object

removeConnection Remove clock or trigger

connection

resetCounters Reset counter channel to initial

count

startBackground Start background operations
startForeground Start foreground operations
stop Stop background operation

wait Block MATLAB until background

operation completes

Events

DataAvailable Notify when acquired data is

available to process.

DataRequired Notify when additional data is

required for output on continuous

generation.

ErrorOccurred Notify when device-related errors

occur.

Examples s = daq.createSession('ni');

s.addAnalogInputChannel('cDAQ1Mod1', 'ai0', 'Voltage');

data = s.startForeground();

See Also daq.createSession | daq.getDevices | daq.getVendors

Functions — Alphabetical List

Purpose

Add hardware channels to analog input or output object

Syntax

```
chans = addchannel(obj,hwch)
chans = addchannel(obj,hwch,index)
chans = addchannel(obj,hwch,'names')
chans = addchannel(obj,hwch,index,'names')
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj An analog input or analog output object.

 $\label{eq:baseline} \mbox{hwch} \qquad \mbox{Specifies the numeric IDs of the hardware channels}$

added to the device object. Any MATLAB vector syntax $\,$

can be used.

index The MATLAB indices to associate with the hardware

channels. Any MATLAB vector syntax can be used provided the vector elements are monotonically

increasing.

'names' A descriptive channel name or cell array of descriptive

channel names.

chans A column vector of channels with the same length as hwch.

Description

chans = addchannel(obj, hwch) adds the hardware channels specified by hwch to the device object obj. The MATLAB indices associated with the added channels are assigned automatically. chans is a column vector of channels.

chans = addchannel(obj,hwch,index) adds the hardware channels specified by hwch to the device object obj. index specifies the MATLAB indices to associate with the added channels.

chans = addchannel(obj, hwch, 'names') adds the hardware channels specified by hwch to the device object obj. The MATLAB indices associated with the added channels are assigned automatically. names is a descriptive channel name or cell array of descriptive channel names.

chans = addchannel(obj,hwch,index,'names') adds the hardware channels specified by hwch to the device object obj. index specifies the MATLAB indices to associate with the added channels. names is a descriptive channel name or cell array of descriptive channel names.

Tips Rules for Adding Channels

- The numeric values you supply for hwch depend on the hardware you access. For National Instruments and Measurement Computing hardware, channels are "zero-based" (begin at zero). For sound cards, channels are "one-based" (begin at one).
- Hardware channel IDs are stored in the HwChannel property and the associated MATLAB indices are stored in the Index property.
- You can add individual hardware channels to multiple device objects.
- For sound cards, you cannot add a hardware channel multiple times to the same device object.
- You can configure sound cards in one of two ways: mono mode or stereo mode. For mono mode, hwch must be 1. For stereo mode, the first hwch value specified must be 1.

Note If you are using National Instruments AMUX-64T multiplexer boards, you must use the addmuxchannel function to add channels.

 When you use the sound card, and only one channel is added to an analog output object the card is put into mono mode. The same signal is output to both channels.

More About MATLAB Indices

Every hardware channel contained by a device object has an associated MATLAB index that is used to reference the channel. Index assignments are made either automatically by addchannel or explicitly with the index argument and follow these rules:

- If index is not specified and no hardware channels are contained by the device object, then the assigned indices automatically start at one and increase monotonically. If hardware channels have already been added to the device object, then the assigned indices automatically start at the next highest index value and increase monotonically.
- If index is specified but the indices are previously assigned, then the requested assignment takes precedence and the previous assignment is reindexed to the next available values. If the lengths of hwch and index are not equal, then an error is returned and no channels are added to the device object.
- The resulting indices begin at one and increase monotonically up to the size of the channel group.
- If you are using scanning hardware, then the indices define the scan order.
- Sound cards cannot be reindexed.

Rules for Adding Channels to National Instruments 1200 Series Boards

When using National Instruments 1200 Series hardware, you need to modify the above rules in these ways:

 Channel IDs are given in reverse order with addchannel. For example, to add eight single-ended channels to the analog input object ai:

```
addchannel(ai,7:-1:0);
```

- The scan order is from the highest ID to the lowest ID (which must be 0).
- There cannot be any gaps in the channel group.

• When channels are configured in differential mode, the hardware IDs are 0, 2, 4, and 6.

More About Descriptive Channel Names

You can assign hardware channels descriptive names, which are stored in the ChannelName property. Choosing a unique descriptive name can be a useful way to identify and reference channels. For a single call to addchannel, you can

- Specify one channel name that applies to all channels that are to be added
- Specify a different name for each channel to be added

If the number of names specified in a single addchannel call is more than one but not equal to the number of channels to be added, then an error is returned. If a channel is to be referenced by its name, then that name must not contain symbols. If you are naming a large number of channels, then the makenames function might be useful. If a channel is not assigned a descriptive name, then it must be referenced by index.

A sound card configured in mono mode is automatically assigned the name Mono, while a sound card configured in stereo mode is automatically assigned the names Left for the first channel and Right for the second channel. You can change these default channel names when the device object is created, or any time after the channel is added.

Examples National Instruments

Suppose you create the analog input object AI1 for a National Instruments board, and add the first four hardware channels (channels 0-3) to it.

```
AI1 = analoginput('nidaq','Dev1');
addchannel(AI1,0:3);
```

The channels are automatically assigned the indices 1-4. If you want to add the first four hardware channels to AI1 and assign descriptive names to the channels.

```
addchannel(AI1,0:3,{'chan1','chan2','chan3','chan4'});
```

Note that you can use the makenames function to create a cell array of channel names. If you add channels 4, 5, and 7 to the existing channel group,

```
addchannel(AI1,[4 5 7]);
```

the new channels are automatically assigned the indices 5-7. Suppose instead you add channels 4, 5, and 7 to the channel group and explicitly assign them indices 1-3.

```
addchannel(AI1,[4 5 7],1:3);
```

The new channels are assigned the indices 1-3, and the previously defined channels are reindexed as indices 4-7. However, if you assigned channels 4, 5, and 7 to indices 6-8, an error is returned because there is a gap in the indices (index 5 has no associated hardware channel).

Sound Card

Suppose you create the analog input object AI1 for a sound card. Most sound cards have only two channels that can be added to a device object. To configure the sound card to operate in mono mode, you must specify hwch as 1.

```
AI1 = analoginput('winsound');
addchannel(AI1,1);
```

The ChannelName property is automatically assigned the value Mono. You can now configure the sound card to operate in stereo mode by adding the second channel.

```
addchannel(AI1,2);
```

The ChannelName property is assigned the values Left and Right for the two hardware channels. Alternatively, you can configure the sound card to operate in stereo mode with one call to addchannel.

```
addchannel(AI1,1:2);
```

addchannel

See Also

delete | makenames | ChannelName | HwChannel | Index

Purpose

Add hardware lines to digital I/O object

Syntax

```
lines = addline(obj,hwline,'direction')
lines = addline(obj,hwline,port,'direction')
lines = addline(obj,hwline,'direction','names')
lines = addline(obj,hwline,port,'direction','names')
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj A digital I/O object.

hwline The numeric IDs of the hardware lines added to the

device object. Any MATLAB vector syntax can be used.

'direction' The line directions can be In or Out, and can be

specified as a single value or a cell array of values.

port The numeric IDs of the digital I/O port.

'names' A descriptive line name or cell array of descriptive

line names.

lines A column vector of lines with the same length as

hwline.

Description

lines = addline(obj,hwline,'direction') adds the hardware lines specified by hwline to the digital I/O object obj. direction configures the lines for either input or output. lines is a row vector of lines.

lines = addline(obj,hwline,port,'direction') adds the hardware lines specified by hwline from the port specified by port to the digital I/O object obj.

lines = addline(obj,hwline,'direction','names') adds the hardware lines specified by hwline to the digital I/O object obj. names is a descriptive line name or cell array of descriptive line names.

lines = addline(obj,hwline,port,'direction','names') adds the hardware lines specified by hwline from the port specified by port to the digital I/O object obj. direction configures the lines for either input or output. names is a descriptive line name or cell array of descriptive line names.

You cannot configure lines independently on devices that use the NI-DAQmx adaptor. Refer to "Line and Port Characteristics" for more information about line configurable devices.

Tips Rules for Adding Lines

- The numeric values you supply for hwline depend on the hardware you access. For National Instruments and Measurement Computing hardware, line IDs are "zero-based" (begin at zero).
- You can add a line only once to a given digital I/O object.
- Hardware line IDs are stored in the HwLine property and the associated MATLAB indices are stored in the Index property.
- For a single call to addline, you can add multiple lines from one port
 or the same line ID from multiple ports. You cannot add multiple
 lines from multiple ports.
- If a port ID is not explicitly referenced, lines are added first from port 0, then from port 1, and so on.
- You can specify the line directions as a single value or a cell array of values. If a single direction is specified, then all added lines have that direction. If supported by the hardware, you can configure individual lines by supplying a cell array of directions.

More About MATLAB Indices

Every hardware line contained by a device object has an associated MATLAB index that is used to reference the line. Index assignments are made automatically by addline and follow these rules:

• If no hardware lines are contained by the device object, then the assigned indices automatically start at one and increase monotonically. If hardware lines have already been added to the device object, then the assigned indices automatically start at the next highest index value and increase monotonically.

- The resulting indices begin at one and increase monotonically up to the size of the line group.
- The first indexed line represents the least significant bit (LSB) and the highest indexed line represents the most significant bit (MSB).

More About Descriptive Line Names

You can assign hardware lines descriptive names, which are stored in the LineName property. Choosing a unique descriptive name can be a useful way to identify and reference lines. For a single call to addline, you can

- Specify one line name that applies to all lines that are to be added
- Specify a different name for each line to be added

If the number of names specified in a single addline call is more than one but differs from the number of lines to be added, then an error is returned. If a line is to be referenced by its name, then that name must not contain symbols. If you are naming a large number of lines, then the makenames function might be useful. If a line is not assigned a descriptive name, then it must be referenced by index.

Examples

Create the digital I/O object dio and add the first four hardware lines (line IDs 0-3) from port 0.

```
dio = digitalio('nidaq','Dev1');
addline(dio,0:3,'in');
```

These lines are automatically assigned the indices 1-4. If you want to add the first four hardware lines to dio and assign descriptive names to the lines,

```
addline(dio,0:3,'in',{'line1','line2','line3','line4'});
```

Note that you can use the makenames function to create a cell array of line names. You can add the first four hardware lines (line IDs 0-3) from port 1 to the existing line group.

```
addline(dio,0:3,1,'out');
```

The new lines are automatically assigned the indices 5-8.

See Also

delete | makenames | HwLine | Index | LineName

addmuxchannel

Purpose

Add hardware channels to analog input objects when using National

Instruments multiplexer board

Syntax

```
addmuxchannel(obj)
addmuxchannel(obj,chanids)
chans = addmuxchannel(...)
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj An analog input object associated with a National

Instruments Traditional NI-DAQ board.

chanids The hardware channel IDs.

chans The channels that are added to obj.

Description

addmuxchannel(obj) adds as many channels to obj as is physically possible based on the number of National Instruments AMUX-64T multiplexer (mux) boards specified by the NumMuxBoards property. For one mux board, 64 channels are added. For two mux boards, 128 channels are added. For four mux boards, 256 channels are added.

addmuxchannel(obj,chanids) adds the channels specified by chanids to obj. chanids refers to the hardware channel IDs of the data acquisition board.

The actual number of channels added to obj depends on the number of mux boards used. For example, suppose you are using a data acquisition board with 16 channels connected to one mux board. If chanid is 0, then addmuxchannel adds four channels. Refer to the *AMUX-64T User Manual* for more information about adding mux channels based on hardware channel IDs and the number of mux boards used.

chans = addmuxchannel(...) returns the channels added to chans.

Tips

This function is not available for National Instruments NI-DAQmx boards.

Before using addmuxchannel, you must set the NumMuxBoards property to the appropriate value. You can use as many as four mux boards with one analog input object. addmuxchannel deletes all channels contained by obj before new channels are added.

Note The Traditional NI-DAQ adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition ToolboxTM object for Traditional NI-DAQ adaptor beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release. See the supported hardware page at www.mathworks.com/products/daq/supportedio.html for more information.

See Also

muxchanidx

Purpose

Create analog input object

Syntax

AI = analoginput('adaptor')
AI = analoginput('adaptor', ID)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Description

AI = analoginput('adaptor') creates the analog input object AI for a sound card having an ID of 0 (adaptor must be winsound). This is the only case where ID is not required.

AI = analoginput('adaptor', ID) creates the analog input object AI for the specified adaptor and for the hardware device with device identifier ID. ID can be specified as an integer or a string.

Note The Traditional NI-DAQ adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition Toolbox™ object for Traditional NI-DAQ adaptor beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release. See the supported hardware page at www.mathworks.com/products/daq/supportedio.html for more information.

Tips

Creating Analog Input Objects

- When an analog input object is created, it does not contain any
 hardware channels. To execute the device object, hardware channels
 must be added with the addchannel function.
- You can create multiple analog input objects that are associated with a particular analog input subsystem. However, you can typically execute only one object at a time.

- The analog input object exists in the data acquisition engine and in the MATLAB workspace. If you create a copy of the device object, it references the original device object in the engine.
- If ID is a numeric value, then you can specify it as an integer or a string. If ID contains any nonnumeric characters, then you must specify it as a string. (See the National Instruments example below.)
- The Name property is automatically assigned a descriptive name that is produced by concatenating *adaptor*, ID, and -AI. You can change this name at any time.

Notes When you create an analog input object, it consumes system resources. To avoid this issue, make sure that you do not create objects in a loop. If you must create objects in a loop, make sure you delete them within the loop.

Hardware Device Identifier

When data acquisition devices are installed, they are assigned a unique number which identifies the device in software. The device identifier is typically assigned automatically and can usually be manually changed using a vendor-supplied device configuration utility. National Instruments refers to this identifier as the device name.

For sound cards, the device identifier is typically not exposed to you through the Microsoft® Windows® environment. However, Data Acquisition Toolbox software automatically associates each sound card with an integer ID value. There are two cases to consider:

- If you have one sound card installed, then ID is 0. You are not required to specify ID when creating an analog input object associated with this device.
- If you have multiple sound cards installed, the first one installed has an ID of 0, the second one installed has an ID of 1, and so on. You must specify ID when creating analog input objects associated with devices not having an ID of 0.

There are two ways you can determine the ID for a particular device:

- Type daghwinfo('adaptor').
- Execute the vendor-supplied device configuration utility.

Input Arguments

adaptor

The hardware driver adaptor name. The supported adaptors are advantech, , mcc, nidag, and winsound.

ID

The hardware device identifier. ID is optional if the device object is associated with a sound card having an ID of 0.

Output Arguments

ΑI

The analog input object.

Properties

Basic Setup

SampleRate Specify per-channel rate at which

analog data is converted to digital

data, or vice versa

SamplesPerTrigger Specify number of samples to

acquire for each channel group member for each trigger that

occurs

TriggerType Specify type of trigger to execute

Channel Properties

Channel Name Specify descriptive channel name

HwChannel Specify hardware channel ID

HwLine Specify hardware line ID

Index MATLAB index of hardware

channel or line

InputRange Specify range of analog input

subsystem

NativeOffset Indicate offset to use when

converting between native data

format and doubles

NativeScaling Indicate scaling to use when

converting between native data

format and doubles

Parent Indicate parent (device object) of

channel or line

SensorRange Specify range of data expected

from sensor

Type Indicate device object type,

channel, or line

Units Specify engineering units label

UnitsRange Specify range of data as

engineering units

Trigger Properties

InitialTriggerTime Absolute time of first trigger

ManualTriggerHw0n Specify hardware device starts at

manual trigger

TriggerChannel Specify channel serving as trigger

source

TriggerCondition Specify condition that must be

satisfied before trigger executes

 ${\tt TriggerConditionValue} \qquad \qquad {\tt Specify\ voltage\ value(s)\ that}$

must be satisfied before trigger

executes

TriggerDelay Specify delay value for data

logging

TriggerDelayUnits Specify units in which trigger

delay data is measured

TriggerFcn Specify callback function to

execute when trigger occurs

TriggerRepeat Specify number of additional

times trigger executes

TriggersExecuted Indicate number of triggers that

execute

TriggerType Specify type of trigger to execute

Logging Properties

LogFileName Specify name of disk file

information is logged to

Logging Indicate whether data is being

logged to memory or disk file

LoggingMode Specify destination for acquired

data

LogToDiskMode Specify whether data, events, and

hardware information are saved

to one or more disk files

Status Properties

Logging Indicate whether data is being

logged to memory or disk file

Running Indicate whether device object is

running

SamplesAcquired Indicate number of samples

acquired per channel

Samples Available Indicate number of samples

available per channel in engine

Hardware Configuration Properties

ChannelSkew Specify time between consecutive

scanned hardware channels

ChannelSkewMode Specify how channel skew is

determined

ClockSource Specify clock that governs

hardware conversion rate

InputType Specify analog input hardware

channel configuration

SampleRate Specify per-channel rate at which

analog data is converted to digital

data, or vice versa

Callback Properties

DataMissedFcn Specify callback function to

execute when data is missed

InputOverRangeFcn Specify callback function to

execute when acquired data exceeds valid hardware range

RuntimeErrorFcn Specify callback function to

execute when run-time error

occurs

SamplesAcquired Indicate number of samples

acquired per channel

SamplesAcquiredFcn Specify callback function to

execute when predefined number of samples is acquired for each

channel group member

SamplesAcquiredFcnCount Specify number of samples to

acquire for each channel group member before samples acquired

event is generated

StartFcn Specify callback function to

execute before device object runs

StopFcn Specify callback function to

execute after device object runs

TimerFcn Specify callback function to

execute when predefined time

period passes

TimerPeriod Specify time period between timer

events

TriggerFcn Specify callback function to

execute when trigger occurs

General Purpose Properties

BufferingConfig Specify per-channel allocated

memory

BufferingMode Specify how memory is allocated

Channel Contain hardware channels

added to device object

EventLog Store information for specific

events

Name Specify descriptive name for the

channel

Tag Specify device object label

Timeout Specify additional waiting time to

extract or queue data

Type Indicate device object type,

channel, or line

UserData Store data to associate with

device object

Examples

To create an analog input object for a National Instruments device defined as 'Dev1':

AI = analoginput('nidaq','Dev1');

To create an analog input object for a Measurement Computing device defined as '1':

AI = analoginput('mcc','1');

Alternatives

See Also

addchannel | daghwinfo

analogoutput

Purpose

Create analog output object

Syntax

```
A0 = analogoutput('adaptor')
A0 = analogoutput('adaptor',ID)
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

'adaptor' The hardware driver adaptor name. The supported adaptors are advantech, mcc, nidag, and winsound.

ID The hardware device identifier. ID is optional if the

device object is associated with a sound card having

an ID of 0.

A0 The analog output object.

Description

AO = analogoutput('adaptor') creates the analog output object AO for a sound card having an ID of O (adaptor must be winsound). This is the only case where ID is not required.

AO = analogoutput('adaptor', ID) creates the analog output object AO for the specified adaptor and for the hardware device with device identifier ID. ID can be specified as an integer or a string.

Note The Traditional NI-DAQ adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition ToolboxTM object for Traditional NI-DAQ adaptor beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release. See the supported hardware page at www.mathworks.com/products/daq/supportedio.html for more information.

Tips More About Creating Analog Output Objects

- When an analog output object is created, it does not contain any hardware channels. To execute the device object, hardware channels must be added with the addchannel function.
- You can create multiple analog output objects that are associated with a particular analog output subsystem. However, you can typically execute only one object at a time.
- The analog output object exists in the data acquisition engine and in the MATLAB workspace. If you create a copy of the device object, it references the original device object in the engine.
- If ID is a numeric value, then you can specify it as an integer or a string. If ID contains any nonnumeric characters, then you must specify it as a string.
- The Name property is automatically assigned a descriptive name that is produced by concatenating *adaptor*, ID, and -AO. You can change this name at any time.

Notes When you create an analog output object, it consumes system resources. To avoid this issue, make sure that you do not create objects in a loop. If you must create objects in a loop, make sure you delete them within the loop.

More About the Hardware Device Identifier

When data acquisition devices are installed, they are assigned a unique number which identifies the device in software. The device identifier is typically assigned automatically and can usually be manually changed using a vendor-supplied device configuration utility. National Instruments refers to this number as the device number.

For sound cards, the device identifier is typically not exposed to you through the Microsoft Windows environment. However, Data Acquisition Toolbox software automatically associates each sound card with an integer ID value. There are two cases to consider:

analogoutput

- If you have one sound card installed, then ID is 0. You are not required to specify ID when creating an analog output object associated with this device.
- If you have multiple sound cards installed, the first one installed has an ID of 0, the second one installed has an ID of 1, and so on. You must specify ID when creating analog output objects associated with devices not having an ID of 0.

There are two ways you can determine the ID for a particular device:

- Type daqhwinfo('adaptor').
- Execute the vendor-supplied device configuration utility.

Properties

Basic Setup Properties

SampleRate Specify per-channel rate at which

analog data is converted to digital

data, or vice versa

TriggerType Specify type of trigger to execute

Channel Properties

Channel Name Specify descriptive channel name

DefaultChannelValue Specify value held by analog

output subsystem

HwChannel Specify hardware channel ID

Index MATLAB index of hardware

channel or line

NativeOffset Indicate offset to use when

converting between native data

format and doubles

NativeScaling Indicate scaling to use when

converting between native data

format and doubles

OutputRange Specify range of analog output

hardware subsystem

Parent Indicate parent (device object) of

channel or line

Type Indicate device object type,

channel, or line

Units Specify engineering units label

UnitsRange Specify range of data as

engineering units

Trigger Properties

InitialTriggerTime Absolute time of first trigger

TriggerFcn Specify callback function to

execute when trigger occurs

TriggersExecuted Indicate number of triggers that

execute

TriggerType Specify type of trigger to execute

Status Properties

Running Indicate whether device object is

running

SamplesAvailable Indicate number of samples

available per channel in engine

SamplesOutput Indicate number of samples

output per channel from engine

Sending Indicate whether data is being

sent to hardware device

Hardware Configuration Properties

ClockSource Specify clock that governs

hardware conversion rate

SampleRate Specify per-channel rate at which

analog data is converted to digital

data, or vice versa

Data Management Properties

MaxSamplesQueued Indicate maximum number of

samples that can be queued in

engine

RepeatOutput Specify number of additional

times queued data is output

Timeout Specify additional waiting time to

extract or queue data

Callback Properties

RuntimeErrorFcn Specify callback function to

execute when run-time error

occurs

SamplesOutputFon Specify callback function to

execute when predefined number of samples is output for each channel group member

SamplesOutputFcnCount Specify number of samples to

output for each channel group member before samples output

event is generated

StartFcn Specify callback function to

execute before device object runs

StopFcn Specify callback function to

execute after device object runs

TimerFcn Specify callback function to

execute when predefined time

period passes

TimerPeriod Specify time period between timer

events

TriggerFcn Specify callback function to

execute when trigger occurs

General Purpose Properties

BufferingConfig Specify per-channel allocated

memory

BufferingMode Specify how memory is allocated

Channel Contain hardware channels

added to device object

EventLog Store information for specific

events

Name Specify descriptive name for the

channel

OutOfDataMode Specify how value held by analog

output subsystem is determined

Tag Specify device object label

Type Indicate device object type,

channel, or line

UserData Store data to associate with

device object

Examples National Instruments

To create an analog output object for a National Instruments device defined as 'Dev1':

```
A0 = analogoutput('nidaq','Dev1');
```

To create an analog output object for a Measurement Computing device defined as '1':

```
A0 = analogoutput('mcc','1');
```

See Also

addchannel | daqhwinfo | Name

Purpose

Convert digital input and output binary vector to decimal value

Syntax

out = binvec2dec(bin)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

bin A binary vector.

out A double array.

Description

out = binvec2dec(bin) converts the binary vector bin to the equivalent decimal number and stores the result in out. All nonzero binary vector elements are interpreted as a 1.

Tips

A binary vector (binvec) is constructed with the least significant bit (LSB) in the first column and the most significant bit (MSB) in the last column. For example, the decimal number 23 is written as the binvec value [1 1 1 0 1].

Note The binary vector cannot exceed 52 values.

Examples

To convert the binvec value [1 1 1 0 1] to a decimal value:

```
binvec2dec([1 1 1 0 1])
ans =
```

23

See Also

dec2binvec

Purpose

Remove device objects from MATLAB workspace

Syntax

clear obj
clear obj.Channel(index)
clear obj.Line(index)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj A device object or array of device objects.

obj.Channel(index) One or more channels contained by obj.
obj.Line(index) One or more lines contained by obj.

Description

clear obj removes obj and all associated channels or lines from the MATLAB workspace, but not from the data acquisition engine.

clear obj.Channel(index) removes the specified channels contained by obj from the MATLAB workspace, but not from the data acquisition engine.

clear obj.Line(index) removes the specified lines contained by obj from the MATLAB workspace, but not from the data acquisition engine.

Tips

Clearing device objects, channels, and lines follows these rules:

- clear does not remove device objects, channels, or lines from the data acquisition engine. Use the delete function for this purpose.
- If multiple references to a device object exist in the workspace, clearing one reference will not invalidate the remaining references.
- You can restore cleared device objects to the MATLAB workspace with the dagfind function.

If you use the help command to display the file help for clear, then you must supply the pathname shown below.

help daq/private/clear

Examples

Create the analog input object ai, copy ai to a new variable aicopy, and then clear the original device object from the MATLAB workspace.

```
ai = analoginput('winsound');
ch = addchannel(ai,1:2);
aicopy = ai;
clear ai
```

Retrieve ai from the engine with daqfind, and show you that ai is identical to aicopy.

```
ainew = daqfind;
isequal(aicopy,ainew)
ans =
    1
```

See Also

daqfind | delete

daqcallback

Purpose

Callback function that displays event information for specified event

Syntax

daqcallback(obj,event)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj A device object.

event

A variable that captures the event information contained

by the EventLog property.

Description

daqcallback(obj,event) is an example callback function that displays information to the MATLAB Command Window. For all events, the information includes the event type and the name of the device object that caused the event to occur. For events that record the absolute time in EventLog, the event time is also displayed. For run-time error events, the error message is also displayed.

Tips

You specify daqcallback as the callback function to be executed for any event by specifying it as the value for the associated callback property. For analog input objects, daqcallback is the default value for the DataMissedFcn and RuntimeErrorFcn properties. For analog output objects, daqcallback is the default value for the RuntimeErrorFcn property.

You can use the showdaqevents function to easily display event information captured by the EventLog property.

Examples

Create the analog input object at and call daqcallback when a trigger event occurs.

```
ai = analoginput('winsound');
addchannel(ai,1);
```

daq callback

```
set(ai, 'TriggerRepeat',3)
set(ai, 'TriggerFcn',@daqcallback)
start(ai)
```

See Also

showdaqevents | DataMissedFcn | EventLog | RuntimeErrorFcn

daqfind

Purpose

Return device objects, channels, or lines from data acquisition engine to MATLAB workspace

Syntax

```
out = daqfind
out = daqfind('PropertyName',PropertyValue,...)
out = daqfind(S)
out = daqfind(obj,'PropertyName',PropertyValue,...)
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

'PropertyName'	A device object, channel, or line property name.
PropertyValue	A device object, channel, or line property value.
obj	A device object, array of device objects, channels, or lines.
S	A structure with field names that are property names and field values that are property values.
out	An array or cell array of device objects, channels, or lines.

Description

out = daqfind returns all device objects that exist in the data acquisition engine. The output out is an array.

out = daqfind('PropertyName', PropertyValue,...) returns all device objects, channels, or lines that exist in the data acquisition engine and have the specified property names and property values. The property name/property value pairs can be specified as a cell array.

out = daqfind(S) returns all device objects, channels, or lines that exist in the data acquisition and have the property names and property values specified by S. S is a structure with field names that are property names and field values that are property values.

out = daqfind(obj, 'PropertyName', PropertyValue,...) returns all device objects, channels, or lines listed by obj that have the specified property names and property values.

Tips More About Finding Device Objects, Channels, or Lines

daqfind is particularly useful in these circumstances:

- A device object is cleared from the MATLAB workspace, and it needs to be retrieved from the data acquisition engine.
- You need to locate device objects, channels, or lines that have particular property names and property values.

Rules for Specifying Property Names and Property Values

- You can use property name/property value string pairs, structures, and cell array pairs in the same call to daqfind. However, in a single call to daqfind, you can specify only device object properties or channel/line properties.
- You must use the same format as returned by get. For example, if get returns the ChannelName property value as Left, you must specify Left as the property value in daqfind (case matters). However, case does not matter when you specify enumerated property values. For example, daqfind will find a device object with a Running property value of On or on.

Examples

You can use dagfind to return a cleared device object.

```
ai = analoginput('winsound');
ch = addchannel(ai,1:2);
set(ch,{'ChannelName'},{'Joe';'Jack'})
clear ai
ainew = daqfind;
```

To return the channel associated with the descriptive name Jack:

```
ch2 = dagfind(ainew, 'ChannelName', 'Jack');
```

daqfind

To return the device object with a sampling rate of 8000 Hz and the descriptive name winsoundO-AI, you can pass a structure to daqfind.

```
S.Name = 'winsoundO-AI';
S.SampleRate = 8000;
daqobj = daqfind(S);
```

See Also

clear | get | propinfo

Purpose

Help for device objects, constructors, adaptors, functions, and properties

Syntax

```
daqhelp
out = daqhelp('name')
out = daqhelp(obj)
out = daqhelp(obj, 'name')
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

'name' A device object, constructor, adaptor, function, or property name.

obj A device object.

out Contains the specified help text.

Description

daqhelp displays a complete listing of Data Acquisition Toolbox constructors and functions along with a brief description of each.


out = daqhelp('name') returns help for the device object, constructor,
adaptor, function, or property specified by name. The help text is
returned to out.

out = daqhelp(obj) returns a complete listing of functions and properties for the device object obj to out. Help for obj's constructor is also displayed.

out = daqhelp(obj, 'name') returns help for name for the specified device object obj to out. name can be a constructor, adaptor, property, or function name.

Tips

As shown below, you can also display help via the Workspace browser by right-clicking a device object, and selecting **Explore > DAQ Help** from the context menu.

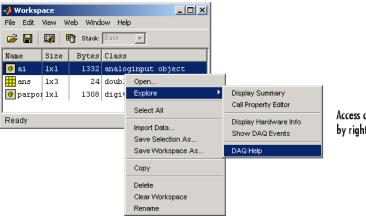

Access context (pop-up) menus by right-clicking a device object.

More About Displaying Help

- When displaying property help, the names in the "See Also" section that contain all uppercase letters are function names. The names that contain a mixture of upper- and lowercase letters are property names.
- When displaying function help, the "See Also" section contains only function names.

Rules for Specifying Names

For the daghelp('name') syntax:

- If name is the name of a constructor, a complete listing of the device object's functions and properties is displayed along with a brief description of each function and property. The constructor help is also displayed.
- You can display object-specific function information by specifying name as object/function. For example, to display the help for an analog input object's getdata function, name is analoginput/getdata.

• You can display object-specific property information by specifying name as obj.property. For example, to display the help for an analog input object's SampleRate property, name is analoginput.SampleRate.

For the daghelp(obj, 'name') syntax:

- If name is the name of a device object constructor and the .m extension is included, the constructor help is displayed.
- If name is the name of a function or property, the function or property help is displayed.

Examples

The following commands are some of the ways you can use daqhelp to obtain help on device objects, constructors, adaptors, functions, and properties.

```
daqhelp('analogoutput');
out = daqhelp('analogoutput.m');
daqhelp set
daqhelp analoginput/peekdata
daqhelp analoginput.TriggerDelayUnits
```

The following commands are some of the ways you can use daqhelp to obtain information about functions and properties for an existing device object.

```
ai = analoginput('winsound');
daqhelp(ai,'InitialTriggerTime')
out = daqhelp(ai,'getsample');
```

See Also

propinfo

daqhwinfo

Purpose

Data acquisition hardware information

Syntax

```
out = daqhwinfo
out = daqhwinfo('adaptor')
out = daqhwinfo(obj)
```

out = daqhwinfo(obj, 'FieldName')

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

'adaptor' The hardware driver adaptor name. The supported

adaptors are advantech, mcc, nidaq, parallel, and

winsound.

obj A device object or array of device objects.

'FieldName' A single field name or a cell array of field names.

Out A structure containing the requested hardware

information.

Description

out = daqhwinfo returns general hardware-related information as a structure to out. The returned information includes installed adaptors, the toolbox and the MATLAB software version, and the toolbox name.

out = daqhwinfo('adaptor') returns hardware-related information for the specified adaptor. The returned information includes the adaptor DLL name, the board names and IDs, and the device object constructor syntax.

Note If you are trying to discover National Instruments, including CompactDAQ or Counter/Timer subsystem devices, use the daq.getDevices method.


out = daqhwinfo('adaptor','FieldName') returns the hardware-related information specified by FieldName for adaptor. FieldName must be a single string. Out is a cell array. You can return a list of valid field names with the daghwinfo('adaptor') syntax.

out = daqhwinfo(obj) returns hardware-related information for the device object obj. If obj is an array of device objects, then out is a 1-by-n cell array of structures where n is the length of obj. The returned information depends on the device object type, and might include the maximum and minimum sampling rates, the channel gains, the hardware channel or line IDs, and the vendor driver version.

out = daqhwinfo(obj, 'FieldName') returns the hardware-related information specified by FieldName for the device object obj. FieldName can be a single field name or a cell array of field names. Out is an m-by-n cell array where m is the length of obj and n is the length of FieldName. You can return a list of valid field names with the daqhwinfo(obj) syntax.

Tips

As shown below, you can also return hardware information via the Workspace browser by right-clicking a device object, and selecting **Display Hardware Info** from the context menu.

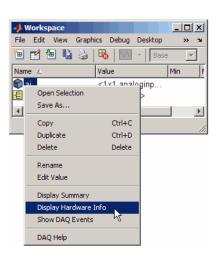

Examples

Display all installed adaptors. Note that this list might be different for your platform.

```
out = daqhwinfo;
out.InstalledAdaptors

ans =
    'advantech'
    'mcc'
    'nidaq'
    'parallel'
    'winsound'
```

To display the device object constructor names for all installed winsound devices:

```
out = daqhwinfo('winsound');
out.ObjectConstructorName
ans =
    'analoginput('winsound',0)'
    'analogoutput('winsound',0)'
```

Create the analog input object ai for a sound card. To display the input ranges for ai:

```
ai = analoginput('winsound');
out = daqhwinfo(ai);
out.InputRanges
ans =
    -1    1
```

To display the minimum and maximum sampling rates for ai:

```
out = daqhwinfo(ai,{'MinSampleRate','MaxSampleRate'})
out =
   [8000] [44100]
```

daqhwinfo

Notes The Traditional NI-DAQ adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition Toolbox[™] object for Traditional NI-DAQ adaptor beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release.

The Parallel adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition $Toolbox^{TM}$ object for 'parallel' beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release. See the supported hardware page at www.mathworks.com/products/daq/supportedio.html for more information.

daqmem

Purpose

Allocate or display analog input and output memory resources

Syntax

out = daqmem
out = daqmem(obj)
daqmem(obj,maxmem)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

Obj A device object or array of device objects.

maxmem The amount of memory to allocate.

Out A structure containing information about memory resources.

Description

out = dagmem returns the object out, which contains several properties describing the memory resources associated with your platform and Data Acquisition Toolbox. The fields are described below.

Field	Description
MemoryLoad	Specifies a number between 0 and 100 that gives a general idea of current memory utilization. 0 indicates no memory use and 100 indicates full memory use.
TotalPhys	Indicates the total number of bytes of physical memory.
AvailPhys	Indicates the number of bytes of physical memory available.

Field	Description
TotalPageFile	Indicates the total number of bytes that can be stored in the paging file. Note that this number does not represent the actual physical size of the paging file on disk.
AvailPageFile	Indicates the number of bytes available in the paging file.
TotalVirtual	Indicates the total number of bytes that can be described in the user mode portion of the virtual address space of the calling process.
AvailVirtual	Indicates the number of bytes of unreserved and uncommitted memory in the user mode portion of the virtual address space of the calling process.
UsedDaq	The total memory used by all device objects.

Note that all the above fields, except for UsedDaq, are identical to the fields returned by Windows' MemoryStatus function.

out = daqmem(obj) returns a 1-by-N structure out containing two fields: UsedBytes and MaxBytes for the device object obj. N is the number of device objects specified by obj. UsedBytes returns the number of bytes used by obj. MaxBytes returns the maximum number of bytes that can be used by obj.

daqmem(obj, maxmem) sets the maximum memory that can be allocated for obj to the value specified by maxmem.

Tips More About Allocating and Displaying Memory Resources

- For analog output objects, daqmem(obj,maxmem) controls the value of the MaxSamplesQueued property.
- If you manually configure the BufferingConfig property, then this
 value supersedes the values specified by daqmem(obj, maxmem) and
 the MaxSamplesQueued property.

Examples

Create the analog input object aiwin for a sound card and the analog input object aini for a National Instruments board, and add two channels to each device object.

```
aiwin = analoginput('winsound');
addchannel(aiwin,1:2);
aini = analoginput('nidaq','Dev1');
addchannel(aini,0:1);
```

To display the total memory used by all existing device objects:

```
out = daqmem;
out.UsedDaq
ans =
69120
```

To configure the maximum memory used by aiwin to 640 KB:

```
daqmem(aiwin,640000)
```

To configure the maximum memory used by each device object with one call to daqmem:

```
dagmem([aiwin aini],[640000 480000])
```

See Also

 ${\tt BufferingConfig} \mid {\tt MaxSamplesQueued}$

Purpose Read Data Acquisition Toolbox (.daq) file for analog input

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

See Also dagread

dagregister

Purpose

Register or unregister hardware driver adaptor

Syntax

```
daqregister('adaptor')
daqregister('adaptor', 'unload')
out = daqregister(...)
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

'adaptor' The hardware driver adaptor name. The supported

adaptors are advantech, mcc, nidaq, parallel, and

winsound.

'unload' Specifies that the hardware driver adaptor is to be

unloaded.

out Captures the message returned by daqregister.

Note The Traditional NI-DAQ adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition Toolbox[™] object for Traditional NI-DAQ adaptor beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release.

The Parallel adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition ToolboxTM object for 'parallel' beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release. See the supported hardware page at www.mathworks.com/products/daq/supportedio.html for more information.

Description

 ${\tt daqregister('adaptor')}$ registers the hardware driver adaptor specified by ${\tt adaptor}$.

Notes You must have administrative privileges to register or unregister hardware driver adaptors.

If you are using a Windows VistaTM machine, you must log on with Administrative privileges and run MATLAB. You should then execute daqregister with elevated permissions. This will allow the User Account Control feature on your computer to run correctly.

For third-party adaptors, *adaptor* must include the full pathname.

dagregister('adaptor', 'unload') unregisters the hardware driver adaptor specified by adaptor. For third-party adaptors, adaptor must include the full pathname.

out = dagregister(...) captures the resulting message in out.

Tips

A hardware driver adaptor must be registered so the data acquisition engine can make use of its services. Unless an adaptor is unloaded, registration is required only once.

For adaptors that are included with the toolbox, registration occurs automatically when you first create a device object. However, you might need to register third-party adaptors manually. In either case, you must install the associated hardware driver before registration can occur.

Examples

The following command registers the sound card adaptor provided with the toolbox.

```
dagregister('winsound');
```

The following command registers the third-party adaptor myadaptor.dll. Note that you must supply the full pathname to dagregister.

```
daqregister('D:/MATLABR12/toolbox/daq/myadaptors/
myadaptor.dll');
```

dagreset

Purpose

Remove device objects, engine MEX-file, and adaptor DLLs from

memory

Syntax

dagreset

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Description

daqreset removes all device objects existing in the engine, and unloads all data acquisition executables loaded by the engine (including the adaptor DLLs and the engine MEX-file).

You should use dagreset to return the MATLAB workspace to a known initial state of having no device objects and no data acquisition MEX-file or DLLs loaded in memory. When the MATLAB workspace returns to this state, the data acquisition hardware is reset.

Note dagreset only affects Data Acquisition Toolbox engine and its adaptors. It does not affect the hardware. To reset the hardware you must use the tools supplied by the hardware vendor. Refer to your hardware documentation for details.

See Also

clear | delete

Purpose

Convert digital input and output decimal value to binary vector

Syntax

out = dec2binvec(dec)
out = dec2binvec(dec,bits)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

dec A decimal value. dec must be nonnegative.

bits Number of bits used to represent the decimal number.

out A logical array containing the binary vector.

Description

out = dec2binvec(dec) converts the decimal value dec to an equivalent binary vector and stores the result as a logical array in out.

out = dec2binvec(dec,bits) converts the decimal value dec to an equivalent binary vector consisting of at least the number of bits specified by bits.

Tips More About Binary Vectors

A binary vector (binvec) is constructed with the least significant bit (LSB) in the first column and the most significant bit (MSB) in the last column. For example, the decimal number 23 is written as the binvec value [1 1 1 0 1].

More About Specifying the Number of Bits

- If bits is greater than the minimum number of bits required to represent the decimal value, then the result is padded with zeros.
- If bits is less than the minimum number of bits required to represent the decimal value, then the minimum number of required bits is used.

• If bits is not specified, then the minimum number of bits required to represent the number is used.

Examples

To convert the decimal value 23 to a binvec value:

```
dec2binvec(23)
ans =
     1     1     1     0     1
```

To convert the decimal value 23 to a binvec value using six bits:

```
dec2binvec(23,6)
ans =
    1    1    1    0    1    0
```

To convert the decimal value 23 to a binvec value using four bits, then the result uses five bits. This is the minimum number of bits required to represent the number.

See Also

binvec2dec

Purpose

Remove device objects, channels, or lines from data acquisition engine

Syntax

```
delete(obj)
delete(obj.Channel(index))
delete(obj.Line(index))
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj A device object or array of device objects.
obj.Channel(index) One or more channels contained by obj.
obj.Line(index) One or more lines contained by obj.

Description

delete(obj) removes the device object specified by obj from the engine. If obj contains channels or lines, they are removed as well. If obj is the last object accessing the driver, then the driver and associated adaptor are unloaded.

delete(obj.Channel(index)) removes the channels specified by index and contained by obj from the engine. As a result, the remaining channels might be reindexed.

delete(obj.Line(index)) removes the lines specified by index and contained by obj from the engine. As a result, the remaining lines might be reindexed.

Tips

Deleting device objects, channels, and lines follows these rules:

- delete removes device objects, channels, or lines from the data acquisition engine but not from the MATLAB workspace. To remove variables from the workspace, use the clear function.
- If multiple references to a device object exist in the workspace, then removing one device object from the engine invalidates the remaining

references. These remaining references should be cleared from the workspace with the clear function.

• If you delete a device object while it is running, then a warning is issued before it is deleted. You cannot delete a device object while it is logging or sending data.

You should use delete at the end of a data acquisition session. You can quickly delete all existing device objects with the command delete(dagfind).

If you use the help command to display the file help for delete, then you must supply the pathname shown below.

help dag/dagdevice/delete

Examples National Instruments

Create the analog input object ai for a National Instruments board, add hardware channels 0-7 to it, and make a copy of hardware channels 0 and 1.

```
ai = analoginput('nidaq','Dev1');
addchannel(ai,0:7);
ch = ai.Channel(1:2);
```

To delete hardware channels 0 and 1:

```
delete(ch)
```

These channels are deleted from the data acquisition engine and are no longer associated with ai. The remaining channels are reindexed such that the indices begin at 1 and increase monotonically to 6. To delete ai:

```
delete(ai)
```

Sound Card

Create the analog input object AI1 for a sound card, and configure it to operate in stereo mode.

```
AI1 = analoginput('winsound');
addchannel(AI1,1:2);
```

You can now configure the sound card for mono mode by deleting hardware channel 2.

```
delete(AI1.Channel(2))
```

If hardware channel 1 is deleted instead, an error is returned.

See Also

clear | dagreset

digitalio

Purpose

Create digital I/O object

Syntax

DIO = digitalio('adaptor', ID)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

'adaptor' The hardware driver adaptor name. The supported

adaptors are advantech, mcc, nidaq, and parallel.

ID The hardware device identifier.

DIO The digital I/O object.

Description

DIO = digitalio('adaptor', ID) creates the digital I/O object DIO for the specified adaptor and for the hardware device with device identifier ID. ID can be specified as an integer or a string.

Notes The Traditional NI-DAQ adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition Toolbox[™] object for Traditional NI-DAQ adaptor beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release.

The Parallel adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition $Toolbox^{TM}$ object for 'parallel' beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release. See the supported hardware page at www.mathworks.com/products/daq/supportedio.html for more information.

Tips More About Creating Digital I/O Objects

- When a digital I/O object is created, it does not contain any hardware lines. To execute the device object, hardware lines must be added with the addline function.
- You can create multiple digital I/O objects that are associated with a particular digital I/O subsystem. However, you can execute only one of these digital I/O objects at a time for the generation of timing events.
- The digital I/O object exists in the data acquisition engine and in the MATLAB workspace. If you create a copy of the device object, it references the original device object in the engine.
- The Name property is automatically assigned a descriptive name that is produced by concatenating *adaptor*, ID, and -DIO. You can change this name at any time.

Note When you create a digital input or output object, it consumes system resources. To avoid this issue, make sure that you do not create objects in a loop. If you must create objects in a loop, make sure you delete them within the loop.

The Parallel Port Adaptor

The toolbox provides basic DIO capabilities through the parallel port. The PC supports up to three parallel ports that are assigned the labels LPT1, LPT2, and LPT3. You can use only these ports. If you add additional ports to your system, or if the standard ports do not use the default memory resources, they will not be accessible by the toolbox. For more information about the parallel port, refer to "Parallel Port Characteristics".

More About the Hardware Device Identifier

When data acquisition devices are installed, they are assigned a unique number, which identifies the device in software. The device identifier is typically assigned automatically and can usually be manually

digitalio

changed using a vendor-supplied device configuration utility. National Instruments refers to this number as the device number.

There are two ways you can determine the ID for a particular device:

- Type daghwinfo('adaptor').
- Open the vendor-supplied device configuration utility.

Properties

Common Properties

Line Contain hardware lines added to

device object

Name Specify descriptive name for the

channel

Running Indicate whether device object is

running

Tag Specify device object label

TimerFcn Specify callback function to

execute when predefined time

period passes

TimerPeriod Specify time period between timer

events

Type Indicate device object type,

channel, or line

UserData Store data to associate with

device object

Line Properties

Direction Specify whether line is for input

or output

HwLine Specify hardware line ID

Index MATLAB index of hardware

channel or line

LineName Specify descriptive line name

Parent Indicate parent (device object) of

channel or line

Port Specify port ID

Type Indicate device object type,

channel, or line

Examples

Create a digital I/O object for a National Instruments device defined as 'Dev1'.

DIO = digitalio('nidaq','Dev1');

Create a digital I/O object for a Measurement Computing device defined as '1'.

DIO = digitalio('mcc','1');

Create a digital I/O object for parallel port LPT1.

DIO = digitalio('parallel','LPT1');

See Also

addline | daghwinfo | Name

Purpose

Summary information for device objects, channels, or lines

Syntax

```
disp(obj)
disp(obj.Channel(index))
disp(obj.Line(index))
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj A device object.

obj.Channel(index) One or more channels contained by obj.

One or more lines contained by obj.

Description

disp(obj) displays summary information for the specified device object obj, and any channels or lines contained by obj. Typing obj at the Command Window produces the same summary information.

disp(obj.Channel(index)) displays summary information for the specified channels contained by obj. Typing obj.Channel(index) at the Command Window produces the same summary information.


disp(obj.Line(index)) displays summary information for the specified lines contained by obj. Typing obj.Line(index) at the Command Window produces the same summary information.

Tips

You can invoke disp by typing the device object at the MATLAB Command Window or by excluding the semicolon when

- Creating a device object
- Adding channel or lines
- Configuring property values using the dot notation

As shown below, you can also display summary information via the Workspace browser by right-clicking a device object, a channel object, or a line object and selecting **Explore > Display Summary** from the context menu.

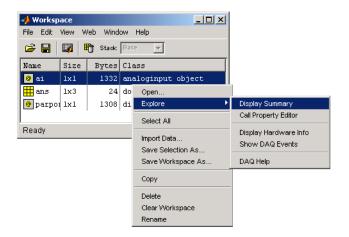

Access context (pop-up) menus by right-clicking a device object.

Examples

All the commands shown below produce summary information for the device object AI or the channels contained by AI.

```
AI = analoginput('winsound')
chans = addchannel(AI,1:2)
AI.SampleRate = 44100
AI.Channel(1).ChannelName = 'CH1'
chans
```

Remove analog input data from data acquisition engine

Syntax

```
flushdata(obj)
flushdata(obj, 'mode')
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj An analog input object or array of analog input objects.

'mode' Specifies how much data is removed from the engine.

Description

flushdata(obj) removes all data from the data acquisition engine and resets the SamplesAvailable property to zero.

flushdata(obj,'mode') removes data from the data acquisition engine depending on the value of mode:

- If mode is all, then flushdata removes all data from the engine and the SamplesAvailable property is set to 0. This is the same as flushdata(obj).
- If mode is triggers, then flushdata removes the data acquired during one trigger. triggers is a valid choice only when the TriggerRepeat property is greater than 0 and the SamplesPerTrigger property is not inf. The data associated with the oldest trigger is removed first.

Examples

Create the analog input object ai for a National Instruments board and add hardware channels 0-7 to it.

```
ai = analoginput('nidaq','Dev1');
addchannel(ai,0:7);
```

A 2-second acquisition is configured and the device object is executed.

```
set(ai, 'SampleRate',2000)
duration = 2;
ActualRate = get(ai, 'SampleRate');
set(ai, 'SamplesPerTrigger',ActualRate*duration)
start(ai)
wait(ai,duration+1)
```

Four thousand samples will be acquired for each channel group member. To extract 1000 samples from the data acquisition engine for each channel:

```
data = getdata(ai,1000);
```

You can use flushdata to remove the remaining 3000 samples from the data acquisition engine.

```
flushdata(ai)
ai.SamplesAvailable
ans =
    0
```

See Also

getdata | SamplesAvailable | SamplesPerTrigger | TriggerRepeat

Device object properties

Syntax

```
out = get(obj)
out = get(obj.Channel(index))
out = get(obj.Line(index))
out = get(obj,'PropertyName')
out = get(obj.Channel(index),'PropertyName')
out = get(obj.Line(index),'PropertyName')
get(...)
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

Obj A device object or array of device objects.

obj.Channel(index) One or more channels contained by obj.

obj.Line(index) One or more lines contained by obj.

'PropertyName' A property name or a cell array of property

names.

Description

out = get(obj) returns the structure out, where each field name is the name of a property of obj and each field contains the value of that property.

out = get(obj.Channel(index)) returns the structure out, where each
field name is the name of a channel property of obj and each field
contains the value of that property.

out = get(obj.Line(index)) returns the structure out, where each
field name is the name of a line property of obj and each field contains
the value of that property.

out = get(obj,'PropertyName') returns the value of the property specified by PropertyName to out. If PropertyName is replaced by a

1-by-n or n-by-1 cell array of strings containing property names, then get returns a 1-by-n cell array of values to out. If obj is an array of data acquisition objects, then out will be an m-by-n cell array of property values where m is equal to the length of obj and n is equal to the number of properties specified.

out = get(obj.Channel(index), 'PropertyName') returns the value of PropertyName to out for the specified channels contained by obj. If multiple channels and multiple property names are specified, then out is an m-by-n cell array where m is the number of channels and n is the number of properties.

out = get(obj.Line(index), 'PropertyName') returns the value of PropertyName to out for the specified lines contained by obj. If multiple lines and multiple property names are specified, then out is an m-by-n cell array where m is the number of lines and n is the number of properties.

get(...) displays all property names and their current values for the specified device object, channel, or line. Base properties are displayed first followed by device-specific properties.

Tips

If you use the help command to display the file help for get, then you must supply the pathname shown below.

help dag/dagdevice/get

Examples

Create the analog input object ai for a sound card and configure it to operate in stereo mode.

```
ai = analoginput('winsound');
addchannel(ai,1:2);
```

The commands shown below are some of the ways you can use get to return property values.

```
chan = get(ai, 'Channel');
out = get(ai, {'SampleRate', 'TriggerDelayUnits'});
out = get(ai);
```

```
get(chan(1), 'Units')
get(chan, {'Index', 'HwChannel', 'ChannelName'})
See Also
set | setverify
```

Extract analog input data, time, and event information from data acquisition engine

Syntax

```
data = getdata(obj)
data = getdata(obj,samples)
data = getdata(obj,samples,'type')
[data,time] = getdata(...)
[data,time,abstime] = getdata(...)
[data,time,abstime,events] = getdata(...)
[data,...] = getdata(obj, 'P1', V1, 'P2', V2,...)
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj	An analog input object.
samples	The number of samples to extract. If samples is not specified, the number of samples extracted is given by the SamplesPerTrigger property.
'type'	Specifies the format of the extracted data as double (the default) or as native.
data	An m-by-n array, where m is the number of samples extracted and n is the number of channels contained by obj.
time	An m-by-1 array of relative time values in seconds, where m is the number of samples extracted. time = 0 is defined as the point at which data logging begins, i.e., when the Logging property of obj is set to On. Measurement of time, with respect to 0, continues until the acquisition is stopped, i.e., when the Logging property of obj is set to Off.

abstime The absolute time of the first trigger returned as a

 ${\tt clock}$ vector. This value is identical to the value stored

by the InitialTriggerTime property.

events A structure containing a list of events that occurred

during the time period the samples were extracted.

Description

data = getdata(obj) extracts the number of samples specified by the SamplesPerTrigger property for each channel contained by obj. data is an m-by-n array, where m is the number of samples extracted and n is the number of channels.

data = getdata(obj, samples) extracts the number of samples specified by samples for each channel contained by obj.

data = getdata(obj, samples, 'type') extracts the number of samples specified by samples in the format specified by type for each channel contained by obj.

[data,time] = getdata(...) returns data as sample-time pairs. time is an m-by-1 array of relative time values, where m is the number of samples returned in data. Each element of time indicates the relative time, in seconds, of the corresponding sample in data, measured with respect to the first sample logged by the engine.

[data,time,abstime] = getdata(...) extracts data as sample-time pairs and returns the absolute time of the trigger. The absolute time is returned as a clock vector and is identical to the value stored by the InitialTriggerTime property.

[data,time,abstime,events] = getdata(...) extracts data as sample-time pairs, returns the absolute time of the trigger, and returns a structure containing a list of events that occurred during the time period the samples were extracted. The events returned are a sub set of those stored by the EventLog property.

[data,...] = getdata(obj, 'P1', V1, 'P2', V2,...) specifies the number of samples to be returned, the format of the data matrix, and whether to return a tscollection object.

The following	table shows a	summary	of properties.
			0- P-0P-0-0-0.

Property	Description
Samples	Specify the number of samples to return.
DataFormat	Specify the data format as double (default) or native.
OutputFormat	Specify the output format as matrix (default) or tscollection.

Note When the ClockSource property for this function is set to one of the External options, the timing will be controlled externally and the values returned in the time variable will not accurately reflect the actual relative time of each sample. It is however an approximation based on the SampleRate you have configured.

Tips More About getdata

- In most circumstances, getdata returns all requested data and does not miss any samples. In the unlikely event that the engine cannot keep pace with the hardware device, it is possible that data is missed. If data is missed, the DataMissedFcn property is called and the device object is stopped.
- getdata is a *blocking* function because it returns execution control to the MATLAB workspace only when the requested number of samples is extracted from the engine for each channel group member.
- You can issue ^C (Ctrl+C) while getdata is blocking. This will not stop the acquisition but will return control to the MATLAB software.
- The amount of data that you can extract from the engine is given by the SamplesAvailable property.
- It is a good practice to use a wait command before your getdata command if the getdata is going to get all data returned by the

getdata

analog input subsystem. For example, if your analog input object is ai and you have set duration to be the number of seconds for the acquisition, you could add the following line right before the getdata:

wait(ai,duration+1)

- Setting the OutputFormat property to tscollection causes getdata to return a tscollection object. In this case, only the data left-hand argument is used.
- For more information on using the Time Series functionality, see "Example: Time Series Objects and Methods" in the MATLAB documentation.

More About Extracting Data From the Engine

- After the requested data is extracted from the engine, the SamplesAvailable property value is automatically reduced by the number of samples returned.
- If the requested number of samples is greater than the samples to be acquired, then an error is returned.
- If the requested data is not returned in the expected amount of time, an error is returned. The expected time to return data is given by the time it takes the engine to fill one data block plus the time specified by the Timeout property.
- If multiple triggers are included in a single getdata call, a NaN is inserted into the returned data and time arrays and the absolute time returned is given by the first trigger.
- When you use multiple immediate triggers Data Acquisition Toolbox cannot determine the "dead" time between triggers. Because of this, the toolbox assumes the "dead" time = 1 sample. For example if the sample rate is 1000 samples per second the toolbox assumes the "dead" time between triggers is one millisecond. The time argument returned by getdata reflects this assumption.

Examples

Create the analog input object at for a National Instruments board and add hardware channels 0 to 3 to it.

```
ai = analoginput('nidaq','Dev1');
addchannel(ai,0:3);
```

Configure a 1-second acquisition with SampleRate set to 1000 samples per second and SamplesPerTrigger set to 1000 samples per trigger.

```
set(ai, 'SampleRate',1000)
set(ai, 'SamplesPerTrigger',1000)
start(ai)
```

The following getdata command blocks execution control until all sample-time pairs, the absolute time of the trigger, and any events that occurred during the getdata call are returned.

```
wait(ai,1)
[data,time,abstime,events] = getdata(ai);
```

data is returned as a 1000-by-4 array of doubles, time is returned as a 1000-by-1 vector of relative times, abstime is returned as a clock vector, and events is returned as a 3-by-1 structure array.

To extract the 1000 data samples from hardware channel 0 only, examine the first column of data.

```
chan0 data = data(:,1);
```

The three events returned are the start event, the trigger event, and the stop event. To return specific event information about the stop event, you must access the Type and Data fields.

```
EventType = events(3).Type;
EventData = events(3).Data;
```

See Also

```
flushdata | getsample | peekdata | timeseries | tscollection
| wait | DataMissedFcn | EventLog | SamplesAvailable |
SamplesPerTrigger | Timeout
```

getsample

Purpose

Immediately acquire one analog input sample

Syntax

sample = getsample(obj)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj An analog input object.

sample A row vector containing one sample for each channel

contained by obj.

Description

sample = getsample(obj) immediately returns a row vector containing one sample for each channel contained by obj.

Tips

Using getsample is a good way to test your analog input configuration. Additionally:

- getsample does not store samples in, or extract samples from, the data acquisition engine.
- You can execute **getsample** at any time after channels have been added to obj.
- getsample is not supported for sound cards and Dynamic Signal Acquisition and Generation (DSA) cards.

Note Refer to the "Hardware Limitations by Vendor" section before you access National Instruments devices with the NI-DAQmx adaptor simultaneously from multiple applications.

Examples

Create the analog input object ai and add eight channels to it.

```
ai = analoginput('nidaq','Dev1');
ch = addchannel(ai,0:7);
```

The following command returns one sample for each channel.

```
sample = getsample(ai);
```

See Also

getdata | peekdata

Read values from digital input and output lines

Syntax

```
out = getvalue(obj)
out = getvalue(obj.Line(index))
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj A digital I/O object.

obj.Line(index) One or more lines contained by obj.

out A binary vector.

Description

out = getvalue(obj) returns the current value from all lines contained by obj as a binary vector to out.

out = getvalue(obj.Line(index)) returns the current value from the lines specified by obj.Line(index).

Tips

More About Reading Values from Lines

- By default, out is returned as a binary vector (binvec). A binvec value is constructed with the least significant bit (LSB) in the first column and the most significant bit (MSB) in the last column. For example, the decimal number 23 is written as the binvec value [1 1 1 0 1].
- You can convert a binvec value to a decimal value with the binvec2dec function.
- If obj contains lines from a port-configurable device, the data acquisition engine will automatically read from all the lines even if they are not contained by the device object.
- When obj contains lines configured for output, getvalue returns the most recently output value set by putvalue. If you have not called

putvalue since you created the digitalio object, then getvalue returns a 0. getvalue cannot ascertain the current output value on the hardware.

Note Refer to the "Hardware Limitations by Vendor" section before you access National Instruments devices with the NI-DAQmx adaptor simultaneously from multiple applications.

Examples

Create the digital I/O object dio and add eight input lines to it.

```
dio = digitalio('nidaq','Dev1');
lines = addline(dio,0:7,'in');
```

To return the current values from all lines contained by dio as a binvec value:

```
out = getvalue(dio);
```

See Also

binvec2dec

inspect

Purpose

Open Property Inspector

Syntax

inspect(obj)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj

An object or an array of objects.

Description

inspect(obj) opens the Property Inspector and allows you to inspect and set properties for the object obj.

Tips

You can also open the Property Inspector via the Workspace browser by double-clicking an object in the Workspace list.

The Property Inspector does not automatically update its display. To refresh the Property Inspector, open it again.

Examples

Create the analog input object ai for a sound card and add two channels.

```
ai = analoginput('winsound');
addchannel(ai,1:2);
```


Open the Property Inspector for the object ai.

```
inspect(ai)
```

The Property Inspector is shown below.

You can expand the properties that are arrays of objects. In the following figure, the Channel property is expanded to enumerate the individual channel objects that make up this property.

You can also expand these individual channel objects to display their own properties, as shown for channel 1.

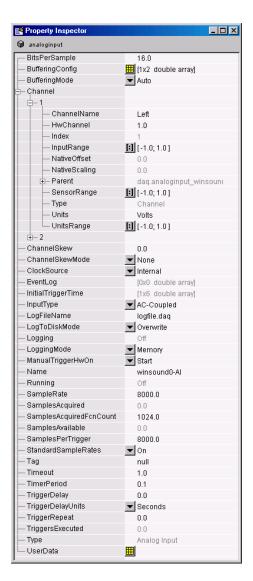

See Also dagfind | daghelp | get | propinfo | set

Check for channels

Syntax

out = ischannel(obj.Channel(index))

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj. Channel (index) One or more channels contained by obj.

out A logical value.

Description

out = ischannel(obj.Channel(index)) returns a logical 1 to out if obj.Channel(index) is a channel. Otherwise, a logical 0 is returned.

Tips

ischannel does not determine if channels are valid (associated with hardware). To check for valid channels, use the isvalid function.

Typically, you use ischannel directly only when you are creating your own files.

Examples

Suppose you create the function myfunc for use with Data Acquisition Toolbox software. If myfunc is passed one or more channels as an input argument, then the first thing you should do in the function is check if the argument is a channel.

```
function myfunc(chan)
% Determine if a channel was passed.
if ~ischannel(chan)
   error('The argument passed is not a channel.');
end
```

You can examine Data Acquisition Toolbox software files for examples that use ischannel.

ischannel

See Also isvalid

Check for lines

Syntax

out = isdioline(obj.Line(index))

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj.Line(index) One or more lines contained by obj.

out A logical value.

Description

out = isdioline(obj.Line(index)) returns a logical 1 to out if obj.Line(index) is a line. Otherwise, a logical 0 is returned.

Tips

isdioline does not determine if lines are valid (associated with hardware). To check for valid lines, use the isvalid function.

Typically, you use isdioline directly only when you are creating your own files.

Examples

Suppose you create the function myfunc for use with Data Acquisition Toolbox software. If myfunc is passed one or more lines as an input argument, then the first thing you should do in the function is check if the argument is a line.

```
function myfunc(line)
% Determine if a line was passed.
if ~isdioline(line)
    error('The argument passed is not a line.');
end
```

You can examine Data Acquisition Toolbox software files for examples that use isdioline.

isdioline

See Also isvalid

Determine whether analog input object is logging data

Syntax

```
bool = islogging(obj)
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Description

bool = islogging(obj) returns true if the analog input object obj is logging data, otherwise false. An analog input object is logging if the value of its Logging property is set to On.

If obj is an array of analog input objects, bool is a logical array where each element in bool represents the corresponding element in obj. If an object in obj is logging data, islogging sets the corresponding element in bool to true, otherwise false. If any of the analog input objects in obj is invalid, islogging returns an error.

Examples

Create an analog input object and add a channel.

```
ai = analoginput('winsound');
addchannel(ai, 1)
```

To put the analog input object in a logging state, start acquiring data. The example acquires 10 seconds of data to increase the amount of time that the object remains in the logging state.

```
set(ai, 'SamplesPerTrigger', 10*get(ai, 'SampleRate'))
start(ai)
```

When the call to the start function returns, and the object is still acquiring data, use islogging to check the state of the object.

```
bool = islogging(ai)
bool =
   1
```

Create a second analog input object.

```
ai2 = analoginput('winsound');
```

Start one of the analog input objects again, such as ai, and use islogging to determine which of the two objects is logging.

```
start(ai)
bool = islogging([ai ai2])
bool =
    1    0
```

See Also

isrunning | issending | start | stop | Logging | LoggingMode

Determine whether device object is running

Syntax

```
bool = isrunning(obj)
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Description

bool = isrunning(obj) returns true if the device object obj is running, otherwise false. A device object is running if the value of its Running property is set to On.

If obj is an array of device objects, bool is a logical array where each element in bool represents the corresponding element in obj. If an object in obj is running, the isrunning function sets the corresponding element in bool to true, otherwise false. If any of the device objects in obj is invalid, isrunning returns an error.

Examples

Create an analog input object and add a channel.

```
ai = analoginput('winsound');
addchannel(ai, 1)
```

To put the analog input object in a running state, configure a manual trigger and then start the object.

```
set(ai, 'TriggerType', 'Manual')
start(ai)
```

Use isrunning to check the state of the object.

```
bool = isrunning(ai)
bool =
    1
```

Create an analog output object.

```
ao = analogoutput('winsound');
```

Use isrunning to determine which of the two objects is running.

```
bool = isrunning([ai ao])
bool =
    1    0
```

See Also

islogging | issending | start | stop | Running

Determine whether analog output object is sending data

Syntax

bool = issending(obj)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Description

bool = issending(obj) returns true if the analog output object obj is sending data to the hardware device, otherwise false. An analog output object is sending if the value of its Sending property is set to On.

If obj is an array of analog output objects, bool is a logical array where each element in bool represents the corresponding element in obj. If an object in obj is sending, the issending function sets the corresponding element in bool to true, otherwise false. If any of the analog output objects in obj is invalid, issending returns an error.

Examples

Create an analog output object and add a channel.

```
ao = analogoutput('winsound');
addchannel(ao, 1);
```

To put the analog output object in a sending state, start acquiring data. The example sends 10 seconds of data to increase the amount of time that the object remains in the sending state.

```
putdata(ao, ones(10*get(ao, 'SampleRate'),1));
start(ao)
```

When the call to the start function returns, and the object is still sending data, use issending to check the state of the object.

```
bool = issending(ao)
bool =
    1
```

```
Create a second analog output object.
```

```
ao2 = analogoutput('winsound');
```

Start one of the analog output objects again, such as ao, and use issending to determine which of the two objects is sending.

```
putdata(ao, ones(10*get(ao, 'SampleRate'),1));
start(ao)
bool = issending([ao ao2])
bool =
    1    0
```

See Also

islogging | isrunning | start | stop | Sending

Determine whether device objects, channels, or lines are valid

Syntax

```
out = isvalid(obj)
out = isvalid(obj.Channel(index))
out = isvalid(obj.Line(index))
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj A device object or array of device objects.

obj.Channel(index) One or more channels contained by obj.

obj.Line(index) One or more lines contained by obj.

out A logical array.

Description

out = isvalid(obj) returns a logical 1 to out if obj is a valid device object. Otherwise, a logical 0 is returned.

out = isvalid(obj.Channel(index)) returns a logical 1 to out if the channels specified by obj.Channel(index) are valid. Otherwise, a logical 0 is returned.

out = isvalid(obj.Line(index)) returns a logical 1 to out if the lines specified by obj.Line(index) are valid. Otherwise, a logical 0 is returned.

Tips

Invalid device objects, channels, and lines are no longer associated with any hardware and should be cleared from the workspace with the clear function.

Typically, you use isvalid directly only when you are creating your own files.

Examples

Create the analog input object ai for a National Instruments board and add eight channels to it.

```
ai = analoginput('nidaq','Dev1');
ch = addchannel(ai,0:7);
```

To verify the device object is valid:

```
isvalid(ai)
ans =
    1
```

To verify the channels are valid:

If you delete a channel, then isvalid returns a logical 0 in the appropriate location:

Typically, you use isvalid directly only when you are creating your own files. Suppose you create the function myfunc for use with Data Acquisition Toolbox software. If myfunc is passed the previously defined device object ai as an input argument,

```
myfunc(ai)
```

the first thing you should do in the function is check if ai is a valid device object.

```
function myfunc(obj)
% Determine if an invalid handle was passed.
if ~isvalid(obj)
```

isvalid

```
error('Invalid data acquisition object passed.'); end
```

You can examine Data Acquisition Toolbox software files for examples that use isvalid.

See Also

clear | delete | ischannel | isdioline

Length of device object, channel group, or line group

Syntax

```
out = length(obj)
out = length(obj.Channel)
out = length(obj.Line)
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj A device object or array of device objects.

obj.Channel The channels contained by obj.
obj.Line The lines contained by obj.

out A double.

Description

out = length(obj) returns the length of the device object obj to out.

out = length(obj.Channel) returns the length of the channel group
contained by obj.

out = length(obj.Line) returns the length of the line group contained by obj.

Examples

Create the analog input object ai for a National Instruments board and add eight channels to it.

```
ai = analoginput('nidaq','Dev1');
aich = addchannel(ai,0:7);
```

Create the analog output object **ao** for a National Instruments board, add one channel to it, and create the device object array aiao.

```
ao = analogoutput('nidaq','Dev1');
aoch = addchannel(ao,0);
```

length

See Also

size

```
aiao = [ai ao]
Index:
          Subsystem:
                               Name:
              Analog Input
                                   nidaqmxDev1-AI
   1
   2
              Analog Output
                                   nidaqmxDev1-A0
To find the length of aiao:
length(aiao)
ans =
     2
To find the length of the analog input channel group:
length(aich)
ans =
     8
```

Load device objects, channels, or lines into MATLAB workspace

Syntax

```
load file
load file obj1 obj2. . .
out = load('file','obj1','obj2',. . .)
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

file

1110	The MAT-me name.
obj1 obj2	Device objects, an array of device objects, channels, or lines.
out	A structure containing the loaded device objects.

Description

load file returns all variables from the MAT-file file into the MATLAB workspace.

The MAT file name

load file obj1 obj2... returns the specified device objects from the MAT-file file into the MATLAB workspace.

out = load('file','obj1','obj2',...) returns the specified device objects from the MAT-file file as a structure to out instead of directly loading them into the workspace. The field names in out match the names of the loaded device objects. If no device objects are specified, then all variables existing in the MAT-file are loaded.

Tips

Loading device objects follows these rules:

- Unique device objects are loaded into the MATLAB workspace as well as the engine.
- If a loaded device object already exists in the engine but not the MATLAB workspace, the loaded device object automatically reconnects to the engine device object.

- If a loaded device object already exists in the workspace or the engine but has different properties than the loaded object, then these rules are followed:
 - The read-only properties are automatically reset to their default values.
 - All other property values are given by the loaded object and a warning is issued stating that property values of the workspace object have been updated.
- If the workspace device object is running, then it is stopped before loading occurs.
- If identical device objects are loaded, then they point to the same device object in the engine. For example, if you saved the array

x = [ai1 ai1 ai2]

only ai1 and ai2 are created in the engine, and x(1) will equal x(2).

- Values for read-only properties are restored to their default values upon loading. For example, the EventLog property is restored to an empty vector. Use the propinfo function to determine if a property is read only.
- Values for the BufferingConfig property when the BufferingMode property is set to Auto, and the MaxSamplesQueued property might not be restored to the same value because both these property values are based on available memory.

Note load is not used to read in acquired data that has been saved to a log file. You should use the dagread function for this purpose.

If you use the help command to display the help for load, then you must supply the pathname shown below.

help dag/private/load

Examples

This example illustrates the behavior of load when the loaded device object has properties that differ from the workspace object.

```
ai = analoginput('winsound');
addchannel(ai,1:2);
save ai
ai.SampleRate = 10000;
load ai
Warning: Loaded object has updated property values.
```

See Also

dagread | propinfo | save

makenames

Purpose

List descriptive channel or line names

Syntax

names = makenames('prefix',index)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

'prefix' A string that constitutes the first part of the name.

index Numbers appended to the end of prefix — any

MATLAB vector syntax can be used to specify index as

long as the numbers are positive.

names An m-by-1 cell array of channel names where m is the

length of index.

Description

names = makenames('prefix',index) generates a cell array of descriptive channel or line names by concatenating prefix and index.

Tips

You can pass names as an input argument to the addchannel or addline function.

If names contains more than one descriptive name, then the size of names must agree with the number of hardware channels specified in addchannel, or the number of hardware lines specified in addline.

If the channels or lines are to be referenced by name, then prefix must begin with a letter and contain only letters, numbers, and underscores. Otherwise the names can contain any character.

Examples

Create the analog input object AI. You can use makenames to define descriptive names for each channel that is to be added to AI.

AI = analoginput('nidaq','Dev1');

muxchanidx

Purpose

Multiplexed scanned analog input channel index

Syntax

scanidx = muxchanidx(obj,muxboard,muxidx)

scanidx = muxchanidx(obj,absmuxidx)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj An analog input object associated with a National

Instruments Traditional NI-DAQ board.

muxboard The multiplexer board.

muxidx The index number of the multiplexed channel.

absmuxidx The absolute index number of the multiplexed

channel.

scanidx The scanning index number of the multiplexed

channel.

Description

scanidx = muxchanidx(obj,muxboard,muxidx) returns the scanning index number of the multiplexed channel specified by muxidx. The multiplexer (mux) board is specified by muxboard. For each mux board, muxidx can range from 0-31 for differential inputs and 0-63 for single-ended inputs. muxboard and muxidx are vectors of equal length.

scanidx = muxchanidx(obj,absmuxidx) returns the scanning index
number of the multiplexed channel specified by absmuxidx. absmuxidx
is the absolute index of the channel independent of the mux board.

For single-ended inputs, the first mux board has absolute index values that range between 0 and 63, the second mux board has absolute index values that range between 64 and 127, the third mux board has absolute index values that range between 128 and 191, the fourth mux board has absolute index values that range between 192 and 255. For

example, the absolute index value of the second single-ended channel on the fourth mux board (muxboard is 4 and muxidx is 1) is 193.

Note The Traditional NI-DAQ adaptor will be deprecated in a future version of the toolbox. If you create a Data Acquisition ToolboxTM object for Traditional NI-DAQ adaptor beginning in R2008b, you will receive a warning stating that this adaptor will be removed in a future release. See the supported hardware page at www.mathworks.com/products/daq/supportedio.html for more information.

Tips

scanidx identifies the column number of the data returned by getdata and peekdata.

Refer to the *AMUX-64T User Manual* for more information about adding mux channels based on hardware channel IDs and the number of mux boards used.

Examples

Create the analog input object ai for a National Instruments board that is connected to four AMUX-64T multiplexers, and add 256 channels to ai using addmuxchannel.

```
ai = analoginput('nidaq',1);
ai.InputType = 'SingleEnded';
ai.NumMuxBoards = 4;
addmuxchannel(ai);
```

The following two commands return a scanned index value of 14.

```
scanidx = muxchanidx(ai,4,1);
scanidx = muxchanidx(ai,193);
```

See Also

addmuxchannel

Convert device objects, channels, or lines to MATLAB code

Syntax

```
obj2mfile(obj,'file')
obj2mfile(obj,'file','syntax')
obj2mfile(obj,'file','all')
obj2mfile(obj,'file','syntax','all')
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

Obj A device object, array of device objects, channels, or lines.

The file that the MATLAB code is written to. The full pathname can be specified. If an extension is not

specified, the .m extension is used.

'syntax' Syntax of the converted the MATLAB code. By default,

the set syntax is used. If dot is specified, then the subscripted referencing syntax is used. If named is specified, then named referencing is used (if defined).

'all' If all is specified, all properties are written to file. If

all is not specified, only properties that are not set to

their default values are written to file.

Description

obj2mfile(obj,'file') converts obj to the equivalent MATLAB code using the set syntax and saves the code to file. By default, only those properties that are not set to their default values are written to file.

obj2mfile(obj,'file','syntax') converts obj to the equivalent MATLAB code using syntax and saves the code to file. The values for syntax can be set, dot, or named. set uses the set syntax, dot uses subscripted assignment (dot notation), and named uses named referencing (if defined).

obj2mfile(obj, 'file', 'all') converts obj to the equivalent MATLAB code using the set syntax and saves the code to file. all specifies that all properties are written to file.

obj2mfile(obj,'file','syntax','all') converts obj including all of obj's properties to the equivalent MATLAB code using syntax and saves the code to file.

Tips

If the UserData property is not empty or if any of the callback properties are set to a cell array of values or a function handle, then the data stored in those properties is written to a MAT-file when the object is converted and saved. The MAT-file has the same name as the file containing the object code (see the example below).

You can recreate the saved device objects by typing the name of the file at the Command Window. You can also recreate channels or lines, by typing the name of the file with a device object as the only input.

Examples

Create the analog input object ai for a sound card, add two channels, and set values for several properties.

```
ai = analoginput('winsound');
addchannel(ai,1:2);
set(ai,'Tag','myai','TriggerRepeat',4)
set(ai,'StartFcn',{@mycallback,2,magic(10)})
```

The following command writes MATLAB code to the files myai.m and myai.mat.

```
obj2mfile(ai, 'myai.m', 'dot')
```

myai.m contains code that recreates the analog input code shown above using the dot notation for all properties that have their default values changed. Because StartFcn is set to a cell array of values, this property appears in myai.m as

```
ai.StartFcn = startfcn1;
and is saved in myai.mat as
```

obj2mfile

```
startfcn1 = {@mycallback,2,magic(10)};
To recreate ai and assign the device object to a new variable ainew:
ainew = myai;
The associated MAT-file, myai.mat, is automatically loaded.
```

Preview most recent acquired analog input data

Syntax

```
data = peekdata(obj,samples)
data = peekdata(obj,samples,'type')
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj	An analog input object.
-----	-------------------------

samples The number of samples to preview for each channel

contained by obj.

'type' Specifies the format of the extracted data as double

(the default) or as native.

data An m-by-n matrix where m is the number of samples

and n is the number of channels.

Description

data = peekdata(obj,samples) returns the latest number of samples specified by samples to data.

data = peekdata(obj,samples,'type') returns the number of samples specified by samples in the format specified by type for each channel contained by obj. If type is specified as native, the data is returned in the native data format of the device. If type is specified as double (the default), the data is returned as doubles.

Tips More About Using peekdata

- Unlike getdata, peekdata is a *nonblocking* function that immediately returns control to the MATLAB workspace. Because peekdata does not block execution control, data might be missed or repeated.
- peekdata takes a "snapshot" of the most recent acquired data and does not remove samples from the data acquisition engine. Therefore,

the SamplesAvailable property value is not affected when peekdata is called.

Rules for Using peekdata

- You can call peekdata before a trigger executes. Therefore, peekdata
 is useful for previewing data before it is logged to the engine or to
 a disk file.
- In most cases, you will call peekdata while the device object is running. However, you can call peekdata once after the device object stops running.
- If samples is greater than the number of samples currently acquired, all available samples are returned with a warning message stating that the requested number of samples were not available.
- If you start an analog input object and LoggingMode is Memory or Disk&Memory, extract the acquired data from the engine, using getdata. You can also flush it out using flushdata. If you do not extract or flush data, you receive a DataMissed event when the amount of acquired data reaches the MaxBytes limit for the object as seen by dagmem. The acquisition then stops.

Examples

Create the analog input object ai for a National Instruments board, add eight input channels, and configure ai for a two-second acquisition.

```
ai = analoginput('nidaq','Dev1');
addchannel(ai,0:7);
set(ai,'SampleRate',2000)
set(ai,'SamplesPerTrigger',4000)
```

After issuing the start function, you can preview the data.

```
start(ai)
data = peekdata(ai,100);
```

peekdata returns 100 samples of data for each of the eight channels added to the object. If 100 samples are not available, then whatever

peekdata

samples are available will be returned and a warning message is issued. The data is not removed from the data acquisition engine.

See Also

daqmem | flushdata | getdata | getsample | SamplesAvailable

propinfo

Purpose

Property characteristics for device objects, channels, or lines

Syntax

```
out = propinfo(obj)
out = propinfo(obj, 'PropertyName')
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj A device object, channels, or lines.

'PropertyName' A valid obj property name.

out A structure whose field names are the property

names for obj (if PropertyName is not specified).

Description

out = propinfo(obj) returns the structure out whose field names are the property names for obj. Each property name in out contains the fields shown below.

Field Name	Description
Туре	The property data type. Possible values are any, callback, double, and string.
Constraint	The type of constraint on the property value. Possible values are bounded, callback, enum, and none.
ConstraintValue	The property value constraint. The constraint can be a range of valid values or a list of valid string values.
DefaultValue	The property default value.

Field Name	Description
ReadOnly	Indicates when the property is read-only. Possible values are always, never, and whileRunning.
DeviceSpecific	If the property is device-specific, a 1 is returned. If a 0 is returned, the property is supported for all device objects of a given type.

out = propinfo(obj, 'PropertyName') returns the structure out for the property specified by PropertyName. If PropertyName is a cell array of strings, a cell array of structures is returned for each property.

Examples

Create the analog input object ai for a sound card and configure it to operate in stereo mode.

```
ai = analoginput('winsound');
addchannel(ai,1:2);
```

To capture all property information for all common ai properties:

```
out = propinfo(ai);
```

To display the default value for the SampleRate property:

```
out.SampleRate.DefaultValue
ans =
     8000
```

To display all the property information for the InputRange property:

propinfo

See Also

daqhelp

Queue analog output data in engine for eventual output

Syntax

putdata(obj,data)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj An analog output object.

data The data to be queued in the engine.

Description

putdata(obj,data) queues the data specified by data in the engine for eventual output to the analog output subsystem. data must consist of a column of data for each channel contained by obj. That is, data must be an m-by-n matrix, where m rows correspond to the number of samples and n columns correspond to the number of channels in obj.

data can consist of doubles or native data types but cannot contain NaNs. data must contain a column of data for each channel contained in obj. If data contains any data points that are not within the UnitsRange of the channel it pertains to, the data points will be clipped to the bounds of the UnitsRange property.

data can be a tscollection object or timeseries object. If data is a tscollection object, there must be one timeseries per channel in obj. If data is a timeseries object, there must be only one channel in obj. If the tscollection or timeseries object contains gaps, or is sampled at a different rate than the SampleRate of obj, the data will be resampled at the rate of obj using a zero order hold.

For more information on using the Time Series functionality, see "Example: Time Series Objects and Methods" in the MATLAB documentation.

Tips More About Queuing Data

- Data must be queued in the engine before obj is executed.
- putdata is a *blocking* function because it returns execution control to the MATLAB workspace only when the requested number of samples is queued in the engine for each channel group member.
- If the value of the RepeatOutput property is greater than 0, then all queued data is automatically requeued until the RepeatOutput value is reached. RepeatOutput must be configured before start is issued.
- After obj executes, you can continue to queue data unless RepeatOutput is greater than 0.
- Due to buffering constraints on certain devices, additional data queued close to the termination of the previous data may not be output to the device. To insure that all data is output, queue additional data well before the device has output all data.
- You can queue data in the engine until the value specified by the MaxSamplesQueued property is reached, or the limitations of your hardware or computer are reached.
- You should not modify the BitsPerSample, InputRange, SensorRange, and UnitsRange properties after calling putdata.
 If these properties are modified, all data is deleted from the data acquisition engine. If you add a channel after calling putdata, all data will be deleted from the buffer.
- The timeseries object must contain a single column of data.

More About Outputting Data

- Data is output as soon as a trigger occurs.
- An error is returned if a NaN is included in the data stream.
- You can specify data as the native data type of the hardware.
- If the output data is not within the range specified by the OutputRange property, then the data is clipped.

- The SamplesOutput property keeps a running count of the total number of samples that have been output per channel.
- The SamplesAvailable property tells you how many samples are ready to be output from the engine per channel. After data is output, SamplesAvailable is automatically reduced by the number of samples sent to the hardware.

Examples

Create the analog output object **ao** for a National Instruments board, add two output channels to it, and generate 10 seconds of data to be output.

```
ao = analogoutput('nidaq','Dev1');
ch = addchannel(ao,0:1);
set(ao,'SampleRate',1000)
data = linspace(0,1,10000)';
```

Before you can output data, it must be queued in the engine using putdata.

```
putdata(ao,[data data])
start(ao)
```

See Also

putsample | timeseries | tscollection | MaxSamplesQueued |
OutputRange | RepeatOutput | SamplesAvailable | SamplesOutput
| Timeout | UnitsRange

putsample

Purpose

Immediately output one analog output sample

Syntax

putsample(obj,data)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj An analog output object.

data The data to be queued in the engine.

Description

putsample(obj,data) immediately outputs the row vector data, which consists of one sample for each channel contained by obj.

Tips

Using putsample is a good way to test your analog output configuration. Additionally:

- putsample does not store samples in the data acquisition engine.
- putsample can be executed at any time after channels have been added to obj.
- putsample is not supported for sound cards and Dynamic Signal Acquisition and Generation (DSA) cards.

Note Refer to the "Hardware Limitations by Vendor" section before you access National Instruments devices with the NI-DAQmx adaptor simultaneously from multiple applications.

Examples

Create the analog output object **ao** for a National Instruments board and add two hardware channels to it.

putsample

```
ao = analogoutput('nidaq','Dev1');
ch = addchannel(ao,0:1);

To call putsample for ao:
putsample(ao,[1 1])

See Also
putdata
```

Write values to digital output lines

Syntax

putvalue(obj,data)
putvalue(obj.Line(index),data)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj A digital I/O object.

obj.Line(index) One or more lines contained by obj.

data A decimal value or binary vector.

Description

putvalue(obj,data) writes data to the hardware lines contained by the digital I/O object obj.

putvalue(obj.Line(index),data) writes data to the hardware lines specified by obj.Line(index).

Tips

More About Writing Values to Lines

- You can specify data as either a decimal value or a binary vector. A binary vector (or *binvec*) is constructed with the least significant bit (LSB) in the first column and the most significant bit (MSB) in the last column. For example, the decimal number 23 is written as the binary vector [1 1 1 0 1].
- If obj contains lines from a port-configurable device, then all lines will be written to even if they are not contained by the device object.
- An error will be returned if data is written to an input line.
- An error is returned if you attempt to write a negative value.
- If a decimal value is written to a digital I/O object and the value is too large to be represented by the hardware, then an error is returned.

Note Refer to the "Hardware Limitations by Vendor" section before you access National Instruments devices with the NI-DAQmx adaptor simultaneously from multiple applications.

Examples

Create the digital I/O object dio and add four output lines to it.

```
dio = digitalio('nidaq','Dev1');
lines = addline(dio,0:3,'out');
```

Write the value 8 as a decimal value and as a binary vector.

```
putvalue(dio,8)
putvalue(dio,[0 0 0 1])
```

Save device objects to MAT-file

Syntax

save file obj1 obj2...

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

file The MAT-file name.

obj1 obj2... One or more device objects or an array of device

objects.

Description

save file saves all the MATLAB variables to the MAT-file file. If an extension is not specified for file, then a .MAT extension is used.

save file obj1 obj2... saves the specified device objects to file.

Tips

Saving device objects follows these rules:

- You can use save in the functional form as well as the command form shown above. When using the functional form, you must specify the filename and device objects as strings.
- Samples associated with a device object are not stored in the MAT-file. You can bring these samples into the MATLAB workspace with the getdata function, and then save them to the MAT-file using a separate variable name. You can also log samples to disk by configuring the LoggingMode property to Disk or Disk&Memory.
- Values for read-only properties are restored to their default values upon loading. For example, the EventLog property is restored to an empty vector. Use the propinfo function to determine if a property is read only.

 Values for the BufferingConfig property (if the BufferingMode property is set to Auto) and the MaxSamplesQueued property might not be restored because both these property values are based on available memory.

If you use the help command to display the help for save, then you must supply this pathname:

help daq/private/save

See Also getdata | load | propinfo

Configure or display device object properties

Syntax

```
set(obj)
props = set(obj)
set(obj,'PropertyName')
props = set(obj,'PropertyName')
set(obj,'PropertyName',PropertyValue,...)
set(obj,PN,PV)
set(obj,S)
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj A device object, array of device objects, channels,

or lines.

'PropertyName' A property name.

PropertyValue A property value.

PN A cell array of property names.
PV A cell array of property values.

S A structure whose field names are device object,

channel, or line properties.

props A structure array whose field names are the

property names for obj, or a cell array of possible

values.

Description

set(obj) displays all configurable properties for obj. If a property has a finite list of possible string values, then these values are also displayed.

props = set(obj) returns all configurable properties to props. props is
a structure array with fields given by the property names, and possible

property values contained in cell arrays. if the property does not have a finite set of possible values, then the cell array is empty.

set(obj,'PropertyName') displays the valid values for the property specified by PropertyName. PropertyName must have a finite set of possible values.

props = set(obj,'PropertyName') returns the valid values for PropertyName to props. props is a cell array of possible values or an empty cell array if the property does not have a finite set of possible values.

set(obj,'PropertyName',PropertyValue,...) sets multiple property values with a single statement. Note that you can use structures, property name/property value string pairs, and property name/property value cell array pairs in the same call to set.

set(obj,PN,PV) sets the properties specified in the cell array of strings PN to the corresponding values in the cell array PV. PN must be a vector. PV can be m-by-n where m is equal to the specified number of device objects, channels, or lines and n is equal to the length of PN.

set(obj,S) where S is a structure whose field names are device object properties, sets the properties named in each field name with the values contained in the structure.

Tips

If you use the help command to display the help for set, then you must supply the pathname shown below.

help dag/dagdevice/set

Examples

Create the analog input object ai for a sound card and configure it to operate in stereo mode.

```
ai = analoginput('winsound');
addchannel(ai,1:2);
```

To display all of ai's configurable properties and their valid values:

set(ai)

To set the value for the SampleRate property to 10000:

```
set(ai, 'SampleRate', 10000)
```

The following two commands set the value for the SampleRate and InputType properties using one call to set.

```
set(ai, 'SampleRate', 10000, 'TriggerType', 'Manual')
set(ai, {'SampleRate', 'TriggerType'}, {10000, 'Manual'})
```

You can also set different channel property values for multiple channels.

```
ch = ai.Channel(1:2);
set(ch,{'UnitsRange','ChannelName'},{[-1 1] 'Name1'; [-2 2]
'Name2'})
```

See Also

get | setverify

Configure and return specified property

Syntax

Actual = setverify(obj,'PropertyName',PropertyValue)

Actual = setverify(obj.Channel(index), 'PropertyName', PropertyValue)
Actual = setverify(obj.Line(index), 'PropertyName', PropertyValue)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj A device object or array of device objects.

'PropertyName' A property name.

PropertyValue A property value.

obj.Channel(index) One or more channels contained by obj.
obj.Line(index) One or more lines contained by obj.

Actual The actual value for the specified property.

Description

Actual = setverify(obj,'PropertyName',PropertyValue) sets PropertyName to PropertyValue for obj, and returns the actual property value to Actual.

Actual =

setverify(obj.Channel(index),'PropertyName',PropertyValue) sets PropertyName to PropertyValue for the channels specified by index, and returns the actual property value to Actual.

Actual =

setverify(obj.Line(index),'*PropertyName*',*PropertyValue*) sets *PropertyName* to *PropertyValue* for the lines specified by index, and returns the actual property value to Actual.

Tips

setverify is equivalent to the commands

```
set(obj,'PropertyName',PropertyValue)
Actual = get(obj,'PropertyName')
```

Using setverify is not required for setting property values, but it does provide a convenient way to verify the actual property value set by the data acquisition engine.

setverify is particularly useful when setting the SampleRate, InputRange, and OutputRange properties because these properties can only be set to specific values accepted by the hardware. You can use the propinfo function to obtain information about the valid values for these properties.

If a property value is specified but does not match a valid value, then

- If the specified value is within the range of supported values,
 - For the SampleRate and InputRange properties, the value is automatically rounded up to the next highest supported value.
 - For all other properties, the value is automatically selected to be the nearest supported value.
- If the value is not within the range of supported values, an error is returned and the current property value remains unchanged.

Examples

Create the analog input object ai for a National Instruments AT-MIO-16DE-10 board, add eight hardware channels to it, and set the sample rate to 10,000 Hz using setverify.

```
ai = analoginput('nidaq','Dev1');
ch = addchannel(ai,0:7);
ActualRate = setverify(ai,'SampleRate',10000);
```

Suppose you use setverify to set the input range for all channels contained by ai to -8 to 8 volts.

```
ActualInputRange = setverify(ai.Channel, 'InputRange', [-8 8]);
```

The InputRange value was actually rounded up to -10 to 10 volts.

```
ActualInputRange{1}
ans =
-10 10
```

See Also

get | propinfo | set | InputRange | OutputRange | SampleRate

showdaqevents

Purpose

Analog input and output event log information

Syntax

showdaqevents(obj)
showdaqevents(obj,index)
showdaqevents(struct)
showdaqevents(struct,index)
out = showdaqevents(...)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj An analog input or analog output object.

index The event index.

struct An event structure.

out A one column cell array of event information.

Description

showdaqevents(obj) displays a summary of the event log for obj.

showdagevents(obj,index) displays a summary of the events specified by index for obj.

showdagevents(struct) displays a summary of the events stored in the structure struct.


showdaqevents(struct,index) displays a summary of the events specified by index stored in the structure struct.

out = showdaqevents(...) outputs the event information to a one column cell array out. Each element of out is a string that contains the event information associated with that index value.

Tips

You can pass a structure of event information to showdaqevents. This structure can be obtained from the getdata function, the daqread function, or the EventLog property.

As shown below, you can also display event information via the Workspace browser by right-clicking a device object and selecting **Explore > Show DAQ Events** from the context menu.

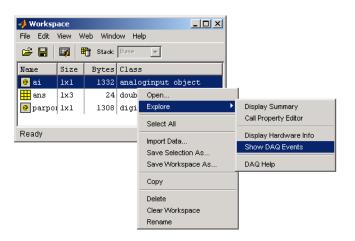

Access context (pop-up) menus by right-clicking a device object.

Examples

Create the analog input object ai for a sound card, add two channels, and configure ai to execute three triggers.

```
ai = analoginput('winsound');
ch = addchannel(ai,1:2);
set(ai,'TriggerRepeat',2)
```

Start ai and display the trigger event information with showdagevents.

```
start(ai)
showdaqevents(ai,2:4)
```

```
2 Trigger#1 ( 17:07:06, 0 ) Channel: N/A
3 Trigger#2 ( 17:07:07, 8000 ) Channel: N/A
4 Trigger#3 ( 17:07:08, 16000 ) Channel: N/A
```

showdaqevents

See Also

daqread | getdata | EventLog

Size of device object, channel group, or line group

Syntax

```
d = size(obj)
[m1,m2,m3,...,mn] = size(obj)
m = size(obj,dim)
d = size(obj.Channel)
[m1,m2,m3,...,mn] = size(obj.Channel)
m = size(obj.Channel,dim)
d = size(obj.Line)
[m1,m2,m3,...,mn] = size(obj.Line)
m = size(obj.Line,dim)
```

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

ohi

00)	if device object of array of device objects.
dim	The dimension.
obj.Channel	The channels contained by obj.
obj.Line	The lines contained by obj.
d	A two-element row vector containing the number of rows and columns in obj.
m1,m2,m3,,mn	Each dimension of obj is captured in a separate variable.
m	The length of the dimension specified by dim.

A device object or array of device objects

Description

d = size(obj) returns the two-element row vector d = [m,n]
containing the number of rows and columns in obj.

[m1,m2,m3,...,mn] = size(obj) returns the length of the first n
dimensions of obj to separate output variables. For example, [m,n] =

size(obj) returns the number of rows to m and the number of columns to n.

m = size(obj,dim) returns the length of the dimension specified by the scalar dim. For example, size(obj,1) returns the number of rows.

d = size(obj.Channel) returns the two-element row vector d = [m,n] containing the number of rows and columns in the channel group obj.Channel.

[m1,m2,m3,...,mn] = size(obj.Channel) returns the length of the first
n dimensions of the channel group obj.Channel to separate output
variables. For example, [m,n] = size(obj.Channel) returns the
number of rows to m and the number of columns to n.

m = size(obj.Channel,dim) returns the length of the dimension specified by the scalar dim. For example, size(obj.Channel,1) returns the number of rows.

d = size(obj.Line) returns the two-element row vector d = [m,n]
containing the number of rows and columns in the line group obj.Line.

[m1,m2,m3,...,mn] = size(obj.Line) returns the length of the first n dimensions of the line group obj.Line to separate output variables. For example, [m,n] = size(obj.Line) returns the number of rows to m and the number of columns to n.

m = size(obj.Line,dim) returns the length of the dimension specified by the scalar dim. For example, size(obj.Line,1) returns the number of rows.

Examples

Create the analog input object ai for a National Instruments board and add eight channels to it.

```
ai = analoginput('nidaq','Dev1');
ch = addchannel(ai,0:7);
```

To find the size of the device object:

```
size(ai)
ans =
```

1 1

To find the size of the channel group:

size(ch) ans = 8

See Also length

softscope

Purpose

Open data acquisition oscilloscope

Syntax

softscope
softscope(obj)
softscope('fname.si')

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj An analog input object.

fname.si Name of the file containing Oscilloscope settings.

Description

softscope opens the Hardware Configuration graphical user interface (GUI), which allows you to configure the hardware device to be used with the Oscilloscope. The Oscilloscope opens when you click the **OK** button, and at least one hardware channel is selected.

softscope(obj) opens the Oscilloscope configured to display the data acquired from the analog input object, obj. obj must contain at least one hardware channel.

softscope('fname.si') pens the Oscilloscope using the settings saved in the softscope file specified by fname. fname is generated from the Oscilloscope's File > Save or File > Save As menu item.

Tips


The Oscilloscope is a graphical user interface (GUI) that allows you to

- Stream acquired data into a display.
- Scale displayed data, and configure triggers and measurements.
- Configure analog input hardware settings.
- Export measurements and acquired data.

To support these tasks, the Oscilloscope includes several helper GUIs, which are described below.

Hardware Configuration

The Hardware Configuration GUI allows you to add channels from a particular hardware device to the Oscilloscope GUI. You can configure the device's sample rate and input type, as well as the input range for each added channel. The GUI shown below is configured to add both sound card channels using the default sample rate.

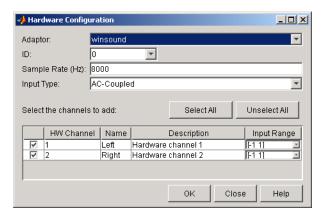

Oscilloscope

The Oscilloscope GUI consists of these panes:

- **Display** pane The display pane contains the hardware channel data (a trace) and the measurements, if defined. The display area also contains labels for each channel's horizontal and vertical units, and indicators for
 - Each trace
 - The trigger level (if defined)
 - The location of the start of the trigger (used for pretriggers)
- **Channel** pane The channel pane lists the hardware channels, math channels, and reference channels that are currently being

softscope


viewed in a display. The Channel Panel also contains knobs for configuring

- The display's horizontal offset and horizontal scale
- The selected channel's vertical offset and vertical scale
- **Trigger** pane The trigger pane allows you to define how data acquisition is initiated. There are three trigger types:
 - One-shot Acquire the specified number of samples once.
 - Continuous Continuously acquire the specified number of samples.
 - Sequence Continuously acquire the specified number of samples, and use the dependent trigger type each time.

For each trigger type, the Oscilloscope begins to acquire data after you press the **Trigger** button.

- **Measurement** pane The measurement pane lists all measurements that are currently being taken. When defining a measurement, you must specify
 - The hardware, math, or reference channel
 - The measurement type
 - Whether the measurement result is drawn as a cursor in the display


The Oscilloscope GUI shown below is configured to display the sound card channels in separate displays.

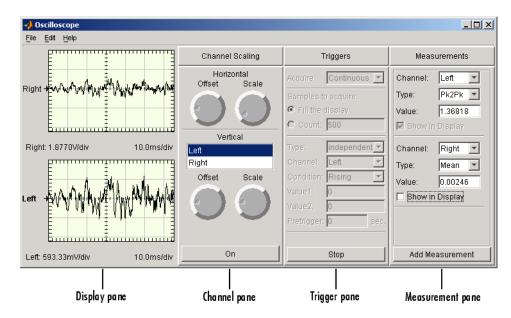

Channel Exporter

The Channel Exporter allows you to export the data associated with a hardware channel, a math channel, or a reference channel. You can export the channel data to one of four destinations:

- The MATLAB workspace as an array
- The MATLAB workspace as a structure
- A MATLAB figure window
- A MAT-file

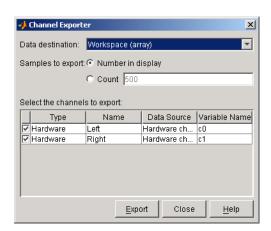


All channels added to the oscilloscope are listed in the GUI.

Measurement Exporter

The Measurement Exporter allows you to export the data associated with a measurement. You can export the measurement to one of three destinations:

- The MATLAB workspace
- A MATLAB figure window
- A MAT-file

The number of measurements exported depends on the BufferSize property value. By default, BufferSize is 1 indicating that the last measurement value calculated is available to export.

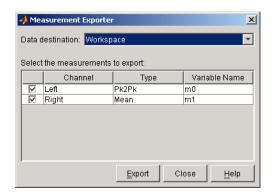


Scope Editor

The Scope Editor consists of two panes:

• Scope — Add and remove displays, the channel pane, the measurement pane, and the trigger pane. Note that you can define as many displays as you want, but there can only be only one channel pane, measurement pane, and trigger pane in the Oscilloscope at a time.

softscope

• **Scope Properties** — Configure properties for the displays, the channel pane, the measurement pane, and the trigger pane.

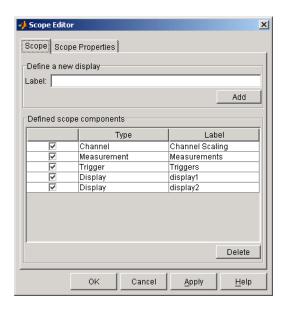


Channel Editor

The Channel Editor consists of three panes:

- **Channel** Add or delete math channels and reference channels, and select which defined channels are available to the Oscilloscope.
- Channel Properties Configure properties for defined hardware channels, math channels, and reference channels.

• Channel Display — Select the Oscilloscope display for each defined channel, or choose to not display a channel.

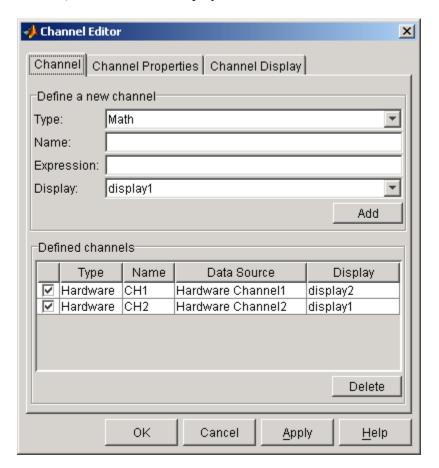


Measurement Editor

The Measurement Editor consists of three panes:

• **Measurement** — Add or delete measurements, and select which defined measurements are available to the Oscilloscope.

softscope

- **Measurement Properties** Configure properties for the defined measurements.
- **Measurement Type** Add or delete measurement types, and select which defined measurement types are available to the Oscilloscope.

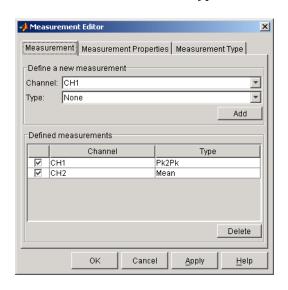

Purpose

Start device object

Syntax

start(obj)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj

A device object or an array of device objects.

Description

start(obj) initiates the execution of the device object obj.

Tips

When start is issued for an analog input or analog output object,

- The callback function specified for StartFcn is executed.
- The Running property is set to On.
- The start event is recorded in the EventLog property.
- Data existing in the engine is flushed.

Although an analog input or analog output object might be executing, data logging or sending is not necessarily initiated. Data logging or sending requires a trigger event to occur, and depends on the TriggerType property value.

For any device object, you can specify start as the value for a callback property.

ai.StopFcn = @start;

Note You typically execute a digital I/O object to periodically update and display its state. Refer to the diopanel example to see this behavior.

start

If you want to synchronize the input and output of data, or you require more control over when your hardware starts, you should use the ManualTriggerHwOn property.

See Also

stop | trigger | EventLog | ManualTriggerHwOn | Running | Sending | TriggerType

Purpose

Stop device object

Syntax

stop(obj)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj

A device object or an array of device objects.

Description

stop(obj) terminates the execution of the device object obj.

Tips

An analog input object automatically stops when the requested samples are acquired or data is missed. An analog output object automatically stops when the queued data is output. These two device objects can also stop executing under one of these conditions:

- The Timeout property value is reached.
- A run-time error occurs.

For analog input objects, stop must be used when the TriggerRepeat property or SamplesPerTrigger property is set to inf. For analog output objects, stop must be used when the RepeatOutput property is set to inf. When stop is issued for either of these device objects,

- The Running property is set to Off.
- The Logging property or Sending property is set to Off.
- The callback function specified for StopFcn is executed.
- The stop event is recorded in the EventLog property.
- All pending callbacks for this object are discarded.

For any device object, you can specify stop as the value for a callback property.

ao.TimerFcn = @stop;

Note Issuing stop is the only way to stop an executing digital I/O object. You typically execute a digital I/O object to periodically update and display its state. Refer to the diopanel example.

See Also

start | trigger | EventLog | Logging | RepeatOutput | Running |
SamplesPerTrigger | Sending | Timeout | TriggerRepeat

Purpose

Manually execute trigger for analog input or output object

Syntax

trigger(obj)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj

An analog input or analog output object or an array of these device objects.

Description

trigger(obj) manually executes a trigger.

Tips

After trigger is issued,

- The absolute time of the trigger event is recorded by the InitialTriggerTime property.
- The Logging property or Sending property is set to On.
- The callback function specified by TriggerFcn is executed.
- The trigger event is recorded in the EventLog property.

You can issue trigger only if TriggerType is set to Manual, Running is On, and Logging is Off.

See Also

start | stop | InitialTriggerTime | Logging | Running | Sending |
TriggerFcn | TriggerType

Purpose

Wait until analog input or output device object stops running

Syntax

wait(obj,waittime)

Note You cannot use the legacy interface on 64-bit MATLAB. See "Using the Session-Based Interface" to acquire and generate data.

Arguments

obj A device object or an array of device objects.

waittime The maximum time to wait for obj to stop running.

Description

wait(obj,waittime) blocks the MATLAB Command Window, and waits for obj to stop running. You specify the maximum waiting time, in seconds, with waittime. waittime overrides the value specified for the Timeout property. If obj is an array of device objects, then wait might wait up to the specified time for each device object in the array.

wait is particularly useful if you want to guarantee that the specified data is acquired before another task is performed.

Tips

If obj is not running when wait is issued, or if an error occurs while obj is running, then wait immediately relinquishes control of the Command Window.

When obj stops running, its Running property is automatically set to Off. obj can stop running under one of these conditions:

- The requested number of samples is acquired (analog input) or sent out (analog output).
- The stop function is issued on that object.
- A run-time error occurs.
- The Timeout property value is reached (waittime supersedes this value).

All callbacks, including the StopFcn, are executed before wait returns.

Examples

Create the analog input object ai for a National Instruments board, add eight channels to it, and configure a 25-second acquisition.

```
ai = analoginput('nidaq','Dev1');
ch = addchannel(ai,0:7);
ai.SampleRate = 2000;
ai.TriggerRepeat = 4;
ai.SamplesPerTrigger = 10000;
```

You can use wait to block the MATLAB Command Window until all the requested data is acquired. Because the expected acquisition time is 25 seconds, the waittime argument is 26. If the acquisition does not complete within this time, then a timeout occurs.

```
start(ai)
wait(ai,26)
```

See Also

EventLog | Running | StopFcn | Timeout

daq.createSession

Purpose Create data acquisition session for specific vendor hardware

Syntax session = daq.createSession ('vendor')

Description session = daq.createSession ('vendor') creates a session object

that you can configure to perform operations using a CompactDAQ

device.

Input Arguments vendor

Enter the vendor name for the device you want to create a session object for. The session-based interface currently supports National Instruments devices only, represented with the abbreviation ni.

Output Arguments

session

The data acquisition session object that represents a session with

hardware from a specific vendor.

Properties

Session acquisition and generation properties:

Channels Array of channel objects

associated with session object

Connections Array of connections in a session

DurationInSeconds Specify duration of acquisition

IsContinuous Specify if operation continues

until manually stopped

IsDone Indicate if operation is complete

IsLogging Indicate if hardware is acquiring

or generating data

IsNotifyWhenDataAvailableExcee(SaAutrol if is set automatically IsNotifyWhenScansQueuedBelowAutontrol if is set automatically

daq.createSession

 ${\tt NotifyWhenDataAvailableExceedsControl\ firing\ of\ DataAvailable}$

event

NotifyWhenScansQueuedBelow Control firing of DataRequired

event

Number of scans for operation

when starting

Range Specify channel measurement

range

Rate of operation in scans per

second

RateLimit Limit of rate of operation based

on hardware configuration

ScansAcquired Number of scans acquired during

operation

ScansOutputByHardware Indicate number of scans output

by hardware

ScansQueued Indicate number of scans queued

for output

Vendor information associated

with session object

Examples

Create a session object s:

s = daq.createSession ('ni')

s =

Data acquisition session using National Instruments hardware:

Will run for 1 second (1000 scans) at 1000 scans/second.

No channels have been added.

daq.createSession

See Also daq.Session | daq.Session.addAnalogInputChannel |

daq.Session.addAnalogOutputChannel | daq.getDevices |

daq.getVendors

How To

• "Using the Session-Based Interface"

Purpose Display available National Instruments devices

Syntax daq.getDevices

device = daq.getDevices

Description daq.getDevices lists National Instruments, including CompactDAQ

devices available to your system. Use

device = daq.getDevices stores this list in the variable device.

Tips Devices not supported by the toolbox are denoted with an *. For a

complete list of supported CompactDAQ devices, see the Supported Hardware page in the Data Acquisition Toolbox area of the MathWorks

Web site.

Output device Arguments The ve

The variable that you want to store a list of National Instruments

devices available to your system.

Examples Get a list of all devices available to your system and store it in the variable d:

d = daq.getDevices

d =

Data acquisition devices:

index	Vendor	Device ID		Description		
1	ni	cDAQ1Mod1	National	Instruments	NI	9205
2	ni	cDAQ1Mod2	National	Instruments	NI	9263
3	ni	cDAQ1Mod3	National	Instruments	ΝI	9234
4	ni	cDAQ1Mod4	National	Instruments	ΝI	9201
5	ni	cDAQ1Mod5	National	Instruments	ΝI	9402
6	ni	cDAQ1Mod6	National	Instruments	ΝI	9213

```
7
             cDAQ1Mod7 National Instruments NI 9219
      ni
8
      ni
             cDAQ1Mod8 National Instruments NI 9265
9
      ni
             cDAQ2Mod1 National Instruments NI 9201
10
      ni
             Dev1
                        National Instruments USB-6211
11
                        National Instruments USB-6218
             Dev2
      ni
12
             Dev3
                        National Instruments USB-6255
      ni
13
      ni
             Dev4
                        National Instruments USB-6363
14
      ni
             Dev5
                        National Instruments PCIe-6363
             PXI1Slot2 National Instruments PXI-4461
15
      ni
16
             PXI1Slot3 National Instruments PXI-4461
      ni
```

To get detailed information about a module on the chassis, type d(index). For example, to get information about NI 9265, which has the index 5, type:

```
d(5)
ans =

ni: National Instruments NI 9402 (Device ID: 'cDAQ1Mod5')
  Counter input subsystem supports:
    Rates from 0.1 to 80000000.0 scans/sec
    4 channels ('ctr0','ctr1','ctr2','ctr3')
    'EdgeCount','PulseWidth','Frequency','Position' measurement types

Counter output subsystem supports:
    Rates from 0.1 to 80000000.0 scans/sec
    4 channels ('ctr0','ctr1','ctr2','ctr3')
    'PulseGeneration' measurement type

This module is in slot 5 of the 'cDAQ-9178' chassis with the name 'cDAQ1'.
```

You can also click on the name of the device in the list. You can now access detailed device information which includes:

subsystem type

daq.getDevices

- rate
- number of available channels
- measurement type

See Also

daq.Session | daq.getVendors | daq.createSession

How To

• "Using the Session-Based Interface"

daq.getVendors

Purpose Display available vendors

Syntax daq.getVendors

vendor = daq.getVendors

Description daq.getVendors lists vendors available to your machine and MATLAB.

vendor = daq.getVendors stores this list in the variable vendor.

Note The current release of Data Acquisition Toolbox only supports National Instruments devices.

Output Arguments

vendor

The variable that stores a list of vendors available to your device.

Examples

Get a list of all vendors available to your machine and MATLAB and store it in the variable v:

v = daq.getVendors

v =

Data acquisition vendor 'National Instruments':

ID: 'ni'

FullName: 'National Instruments'

AdaptorVersion: '3.2 (R2012b)'
DriverVersion: '9.2.3 NI-DAQmx'

IsOperational: true

Data Acquisition Toolbox currently only supports National Instruments, including CompactDAQ devices, denoted with the abbreviation 'ni'.

See Also

daq.Session | daq.getDevices | daq.createSession

daq.getVendors

How To

• "Using the Session-Based Interface"

Purpose

Add analog input channel

Syntax

 $\verb|s.addAnalogInputChannel| (\textit{deviceID}, \textit{channelID},$

measurementType)

ch = s.addAnalogInputChannel(deviceID, channelID,

measurementType)

[ch, idx] = s.addAnalogInputChannel(deviceID, channelID,

measurementType)

Description

s.addAnalogInputChannel(deviceID, channelID, measurementType) adds a channel on the device represented by deviceID, with the specified channelID, and channel measurement type, represented by measurementType, to the session s. Measurement types are vendor specific.

ch = s.addAnalogInputChannel(deviceID, channelID,
measurementType) creates and displays the object ch.

[ch, idx] = s.addAnalogInputChannel(deviceID, channelID, measurementType) creates and displays the object ch, representing the channel that was added and the index, idx, which is an index into the array of the session object's Channels property.

Tips

- Use daq.createSession to create a session object before you use this method.
- To use counter channels, see daq.Session.addCounterInputChannel.

Input Arguments

deviceID

Specify the vendor-defined ID of the device. The specified channel is created for this device. The device ID is the ID of the device that you obtain by calling daq.getDevices.

channelID

Specify the ID of the channel added to the device. You can also add a range of channels.

For NI devices, use either a terminal name, like 'ai2', or a numeric equivalent like 2 for the channel ID.

Note Channel ID is the physical location of the channel on the device. The index for this channel displayed in the session indicates this channels position in the session. If you add a channel with channel ID `ai2' as the first channel in your session, the session index will be 1.

measurementType

Specify a string that represents a vendor-defined measurement type. Measurement types include:

- 'Voltage'
- 'Thermocouple'
- 'Current'
- 'Accelerometer'
- 'RTD'
- 'Bridge'
- 'Microphone'
- 'IEPE'

Output Arguments

ch

Object representing the channel you added.

idx

An index into the array of the session object's Channels property.

Properties ADCTimingMode Set channel timing mode

BridgeMode Specify analog input device bridge

mode

Coupling Specify input coupling mode

Device annel device

ExcitationCurrent Voltage of external source of

excitation

ExcitationSource External source of excitation

ExcitationVoltage Voltage of excitation source

ExternalTriggerTimeout Indicate if external trigger timed

out

ID ID

MaxSoundPressureLevel Sound pressure level for

microphone channels

MeasurementType Type counter channel

measurement

Name Specify descriptive name for the

channel

NominalBridgeResistance Resistance of sensor

R0 Specify resistance value

Range Specify channel measurement

range

RTDConfiguration Specify wiring configuration of

RTD device

RTDType Specify sensor sensitivity

Scans Acquired Number of scans acquired during

operation

Sensitivity Sensitivity of an channel

ShuntLocation Indicate location of channel's

shunt resistor

ShuntResistance Resistance value of channel's

shunt resistor

TerminalConfig Specify terminal configuration

Thermocouple Type Select thermocouple type

Units Specify unit of RTD measurement

Examples

Add an analog input current channel:

```
s = daq.createSession ('ni')
s.addAnalogInputChannel('cDAQ1Mod3','ai0', 'Current');
```

Add an analog input thermocouple channel. Specify output arguments to represent the channel object and the index:

```
s = daq.createSession ('ni')
[ch, idx] = s.addAnalogInputChannel('cDAQ2Mod6', 'ai0', 'Thermocouple')
```

Add analog input voltage channels 0, 2, and 4:

```
s = daq.createSession ('ni')
s.addAnalogInputChannel('cDAQ1Mod1',[0 2 4], 'Voltage');
```

See Also

daq.Session.startBackground | daq.Session.startForeground | daq.Session.addAnalogOutputChannel | daq.Session.removeChannel | daq.Session

How To

· "Using the Session-Based Interface"

Purpose

Add analog output channel

Syntax

- s.addAnalogOutputChannel(deviceID, channelID,
 measurementType)
- ch = s.addAnalogOutputChannel(deviceID, channelID,
 measurementType)
- [ch, idx] = s.addAnalogOutputChannel(deviceID, channelID, measurementType)

Description

- s.addAnalogOutputChannel(deviceID, channelID, measurementType) adds an analog output channel on the device represented by deviceID, with the specified channelID, and channel measurement type, defined by measurementType, on the session object, s. Measurement types are vendor specific.
- ch = s.addAnalogOutputChannel(deviceID, channelID,
 measurementType) creates and displays the object ch, representing the
 channel that was added.
- [ch, idx] = s.addAnalogOutputChannel(deviceID, channelID, measurementType) creates and displays the object ch, representing the channel that was added and the index, idx, which is an index into the array of the session object's Channels property.

Tips

- Use daq.createSession to create a session object before you use this method.
- To use counter channels, see daq.Session.addCounterInputChannel.

Input Arguments

deviceID

Specify the vendor-defined ID of the device on which you will create the analog output channel.

channelID

Specify the ID of the channel added to the device. You can also add a range of channels.

For NI devices, use either a terminal name, like 'ao2', or a numeric equivalent like 2, for the channel ID.

measurementType

Specify a string that represents a vendor-defined measurement type. Measurement types include: .

- 'Voltage'
- 'Current'

Output Arguments

ch

Object representing the channel you added.

idx

An index into the array of the session object's Channels property.

Properties

Device annel device

ExcitationCurrent Voltage of external source of

excitation

ExcitationSource External source of excitation

ExternalTriggerTimeout Indicate if external trigger timed

out

ID ID

MaxSoundPressureLevel Sound pressure level for

microphone channels

MeasurementType Type counter channel

measurement

Name Specify descriptive name for the

channel

Range Specify channel measurement

range

ScansOutputByHardware Indicate number of scans output

by hardware

ScansQueued Indicate number of scans queued

for output

Sensitivity Sensitivity of an channel

TerminalConfig Specify terminal configuration

Examples

Add an analog output voltage channel:

```
s = daq.createSession ('ni')
s.addAnalogOutputChannel('cDAQ1Mod2','ao0', 'Voltage');
```

Add four output current channels:

```
s = daq.createSession ('ni')
s.addAnalogOutputChannel('cDAQ1Mod8',0:3, 'Current');
```

See Also

daq.session.startBackground | daq.session.startForeground |
daq.Session.addAnalogInputChannel | daq.Session.removeChannel
| daq.Session

How To

· "Using the Session-Based Interface"

daq.Session.removeChannel

Purpose

Remove channel from session object

Syntax

s.removeChannel(idx)

Description

s.removeChannel(idx) removes the channel specified by idx from the session objects.

Input Arguments

idx

Index of the channel in the session.

Examples

Start with a session s, with two analog input and two analog output voltage channels and display channel information.

3

s =

Data acquisition session using National Instruments hardware:

No data queued. Will run at 1000 scans/second.

Operation starts immediately.

Number of channels: 4
index Type Device Channel

index	Туре	Device	Channel	${\tt InputType}$	Range	Name
1	ai	cDAQ1Mod4	ai0	SingleEnd	-10 to +10 Volts	
2	ai	cDAQ1Mod4	ai1	SingleEnd	-10 to +10 Volts	
3	ao	cDAQ1Mod2	ao0	n/a	-10 to +10 Volts	
4	ao	cDAQ1Mod2	ao1	n/a	-10 to +10 Volts	

Remove channel 'ai0' currently with the index 1 from the session:

s.removeChannel(1)

To see how the indexes shift after you remove a channel, type:

S

s =

daq.Session.removeChannel

```
Data acquisition session using National Instruments hardware:

No data queued. Will run at 1000 scans/second.

All devices synchronized using cDAQ1 CompactDAQ chassis backplane. (Details)

Number of channels: 3

index Type Device Channel InputType Range Name

1 ai cDAQ1Mod4 ai1 SingleEnd -10 to +10 Volts
2 ao cDAQ1Mod2 ao0 n/a -10 to +10 Volts
3 ao cDAQ1Mod2 ao1 n/a -10 to +10 Volts
```

Remove the first output channel 'aoo' currently at index 2:

```
s.removeChannel(2)
```

The session displays one input and one output channel:

```
s.Channels

ans =

Number of channels: 2

index Type Device Channel InputType Range Name

1 ai cDAQ1Mod4 ai1 SingleEnd -10 to +10 Volts
2 ao cDAQ1Mod2 ao1 n/a -10 to +10 Volts
```

To remove multiple channels together, type the channel indexes as an array. Create a session with multiple channels:

```
s = daq.createSession('ni');
s.addAnalogInputChannel('cDAQ1Mod7',1:3, 'Voltage');
s.addAnalogInputChannel('cDAQ1Mod7',0, 'Current');
s.addCounterInputChannel('Dev2',0:1,'EdgeCount');
s
```

daq.Session.removeChannel

s =

Data acquisition session using National Instruments hardware: Will run for 1 second (2 scans) at 2 scans/second.

Number of channels: 6

index	Туре	Device	Channel	${\tt MeasurementType}$	Range	Name
1	ai	cDAQ1Mod7	ai1	Voltage (Diff)	-60 to +60 Volts	
2	ai	cDAQ1Mod7	ai2	Voltage (Diff)	-60 to +60 Volts	
3	ai	cDAQ1Mod7	ai3	Voltage (Diff)	-60 to +60 Volts	
4	ai	cDAQ1Mod7	ai0	Current	-0.025 to +0.025 A	
5	ci	Dev2	ctr0	EdgeCount	n/a	
6	ci	Dev2	ctr1	EdgeCount	n/a	

Remove channel 2 to channel 5 with one daq.Session.removeChannel command:

See Also

daq.Session.addAnalogInputChannel |
daq.Session.addAnalogOutputChannel

daq.Session.startBackground

Purpose

Start background operations

Syntax

s.startBackground()

Description

s.startBackground() starts the operation of the session object, s, without blocking MATLAB command line and other code. To block MATLAB execution, use daq.Session.startForeground.

When you use startBackground() with analog input channels, the operation uses the DataAvailable event to deliver the acquired data. This event is fired periodically while an acquisition is in progress. For more information, see "Events and Listeners — Concepts".

When you add analog output channels to the session, you must call queueOutputData() before calling startBackground().

During a continuous generation, the DataRequired event is fired periodically to request additional data to be queued to the session. See DataRequired for more information.

By default, the IsContinuous property is set to false and the operation stops automatically. If you have set it to true, use daq.Session.stop to stop background operations explicitly.

Use daq.Session.wait to block MATLAB execution until a background operation is complete.

Tips

- If your session has analog input channels, you must use a DataAvailable event to receive the acquired data in a background acquisition.
- If your session has analog output channels and is continuous, you can use a DataRequired event to queue additional data during background generations.
- Create an acquisition session and add a channel before you use this method. See dag.createSession for more information.
- Call daq.Session.prepare to reduce the latency associated with startup and to preallocate resources.

• Use an ErrorOccurred event to display errors during an operation.

Examples

Acquire data in the background by creating a session and adding a listener to access the acquired data using an anonymous function:

```
s = daq.createSession ('ni');
s.addAnalogInputChannel('cDAQ1Mod1', 'ai0', 'Voltage');
lh = s.addlistener('DataAvailable', @plotData);

function plotData(src,event)
        plot(event.TimeStamps, event.Data)
end

Start the session and perform other MATLAB operations:
s.startBackground();
% perform other MATLAB operations.

Since this is not a continuous operation, the operation stops automatically.

Delete the listener:
delete (lh)
```

For a continuous background generation, add a listener event to queue additional data to be output:

daq.Session.startBackground

```
% perform other MATLAB operations during the generation.
The operation is continuous,
s.stop();
delete(lh);
```

See Also

daq.createSession | daq.Session.startForeground | daq.Session.addAnalogOutputChannel | daq.Session.addAnalogInputChannel | daq.Session.prepare | daq.Session.addlistener | daq.Session | DataAvailable | DataRequired | ErrorOccurred

daq.Session.startForeground

Purpose

Start foreground operations

Syntax

s.startForeground

data = s.startForeground

[data,timeStamps,triggerTime] = s.startForeground

Description

s.startForeground starts operations of the session object, s, and blocks MATLAB command line and other code until the session operation is complete.

To perform other MATLAB operations while the session operation executes, use dag.Session.startBackground.

data = s.startForeground returns the data acquired in the output
parameter, data.

[data,timeStamps,triggerTime] = s.startForeground returns the data acquired, timestamps relative to the time the operation is triggered, and a trigger time indicating the absolute time the operation was triggered.

If a session includes output channels, call queueOutputData before calling startForeground.

You cannot perform continuous operations using startForeground. To perform continuous operations use daq.Session.startBackground.

Output Arguments

data

An mxn matrix of doubles where m is the number of scans acquired, and n is the number of input channels in the session.

timeStamps

The timestamps relative to the time the operation is triggered in an mx1 array where m is the number of scans.

triggerTime

A MATLAB serial date time stamp representing the absolute time when timeStamps = 0.

daq.Session.startForeground

Examples

Acquire data by creating a session with an analog input channel:

```
s = daq.createSession('ni');
s.addAnalogInputChannel('cDAQ1Mod1','ai0','Voltage');
```

Start the acquisition and save the acquired data into the variable data:

```
data = s.startForeground;
```

Generate a signal by creating a session with an analog output channel:

```
s = daq.createSession('ni');
s.addAnalogOutputChannel('cDAQ1Mod2', 'ao0', 'Voltage')
```

Create and queue an output signal and start the generation:

```
outputSignal = linspace(-1, 1, 1000);
s.queueOutputData(outputSignal);
s.startForeground;
```

See Also

```
daq.createSession | daq.Session.startBackground
| daq.Session.addAnalogInputChannel |
daq.Session.addAnalogOutputChannel | daq.Session.prepare |
daq.Session
```

Purpose

Create event listener

Syntax

```
1h = addlistener('eventName',@callback)
1h = addlistener('eventName', @(src, event) expr)
```

Description

1h = addlistener('eventName', @callback) creates a listener for the
specified event, eventName, and fires the callback function, callback.
1h is the variable in which the listener handle is stored. Create a
callback function that executes when the listener detects the specified
event. The callback can be any MATLAB function.

1h = addlistener('eventName', @(src, event) expr) creates a
listener for the specified event, eventName, and fires an anonymous
callback function. The anonymous function uses the specified input
arguments and executes the operation specified in the expression
expr. Anonymous functions provide a quick means of creating simple
functions without storing them to a file. For more information, see
Anonymous Functions.

Tips

• You must delete the listener once the operation is complete.

```
delete (1h)
```

Input Arguments

'eventName'

Name of the event to listen for. Available events include:

- 'DataAvailable'
- 'DataRequired'
- 'ErrorOccurred'

callback

Name of the function to execute when the specified event occurs.

src

The session object, where the event occurred.

event

Specified event object. For more information, see Session Events.

expr

Expression that represents the body of the function.

Output Arguments

lh

Handle to the event listener returned by addlistener. Delete the listener once the operation completes.

Examples

Add a listener to an acquisition session by first:

- Creating a session
- Adding an analog input channel

```
s = daq.createSession('ni');
s.addAnalogInputChannel('cDAQ1Mod1', 'ai0', 'Voltage');
```

Add a listener for the DataAvailable event:

```
lh = s.addlistener('DataAvailable', @plotData);
```

Create the plotData callback function and save it as plotData.m:

Acquire data in the background:

```
s.startBackground();
```

Wait for the operation to complete and delete the listener:

```
delete (lh)
```

Add a listener using an anonymous function to a signal generation session by first:

• Creating a session.

```
• Setting IsContinuous to true.

    Adding two analog output channels.

s = daq.createSession('ni');
s.IsContinuous = true;
s.addAnalogOutputChannel('cDAQ1Mod2', 0:1, 'Voltage');
Create output data for the two channels:
outputData0 = linspace(-1, 1, 1000)';
outputData1 = linspace(-2, 2, 1000)';
Queue the output data:
s.queueOutputData([outputData0 outputData1]);
Add an anonymous listener and generate the signal in the background:
lh = s.addlistener('DataRequired', @(src,event)...
     src.queueOutputData([outputData0 outputData1]));
Generate signals in the background:
s.startBackground();
Perform other MATLAB operations, and then stop the session.
s.stop ()
Delete the listener:
delete (lh)
```

See Also daq.createSession | daq.Session.addAnalogInputChannel

| daq.Session.addAnalogOutputChannel |

daq.Session.startBackground | DataAvailable | DataRequired

| ErrorOccurred

How To• "Working with the Session-Based Interface"

daq.Session.prepare

Purpose Prepare session for operation

Syntax s.prepare

Description s.prepare configures and allocates hardware resources for the session

and reduces the latency of daq.Session.startBackground and

 ${\tt daq.Session.startForeground}$ methods. This method is optional and

is automatically called as needed.

See Also daq.Session.addAnalogInputChannel |

daq.Session.addAnalogInputChannel |

daq.Session.addAnalogOutputChannel | daq.Session.release

Purpose

Block MATLAB until background operation completes

Syntax

```
s.wait ()
s.wait (timeout)
```

Description

s.wait () blocks MATLAB until the background operation completes.

Press Ctrl+C to abort the wait.

 ${\tt s.wait\ (\it timeout)\ blocks\ MATLAB\ until\ the\ operation\ completes\ or}$

the specified time-out occurs.

Tips

- You cannot call wait if you have set the session's IsContinuous property to true.
- To terminate the operation, use daq.Session.stop

Input Arguments

timeout

Maximum time in seconds before the wait throws an error.

Examples

Create a session and add an analog output channel:

```
s = daq.createSession('ni');
s.addAnalogOutputChannel('cDAQ1Mod2', 'ao0', 'Voltage');
```

Queue some output data:

```
s.queueOutputData(zeros(10000,1));
```

Start the session and issue a wait. This blocks MATLAB for all data is output.

```
s.startBackground();
% perform other MATLAB operations.
s.wait()
```

Queue more data and wait for up to 15 seconds:

```
s.queueOutputData(zeros(10000,1));
```

daq.Session.wait

```
s.startBackground();
% perform other MATLAB operations.
s.wait(15);

See Also
daq.Session.startBackground | daq.Session.stop
```

daq.Session.stop

Purpose Stop background operation

Syntax s.stop()

Description s.stop() stops the session and all associated hardware operations in

progress. If your operation has acquired data and the DataAvailable event has not yet fired, the stop command will fire the event and

deliver the data acquired up to that point.

See Also daq.Session.startBackground | daq.Session.wait | daq.Session

Purpose

Release session resources

Syntax

s.release ()

Description

s.release () releases all session hardware resources.

In most cases this method is optional. Data Acquisition Toolbox has exclusive access to the data acquisition hardware that it is using in an operation. This improves performance and is handled automatically when you associate a device to a session object. In certain cases, like when you use multiple session objects with the same hardware, or you use applications other than MATLAB to access the hardware, you may choose to explicitly release the hardware.

Hardware resources associated with a session are automatically released when you delete the session object, or you assign a different value to the variable containing your session object.

Examples

Create a session and add an analog input voltage channel and acquire data in the foreground:

```
s1 = daq.createSession('ni');
s1.addAnalogInputChannel('cDAQ3Mod1', 'ai0', 'Voltage');
s1.startForeground;
```

Release the session hardware and create another session object with an analog input voltage channel on the same device as the previous session. Acquire in the foreground:

```
s1.release()
s2 = daq.createSession('ni');
s2.addAnalogInputChannel('cDAQ3Mod1', 'ai2', 'Voltage');
s2.startForeground;
```

See Also

```
daq.Session.prepare | daq.Session.startForeground |
daq.Session.startBackground | daq.Session
```

daq.Session.inputSingleScan

Purpose

Acquire single scan from all input channels

Syntax

```
data = s.inputSingleScan()
[data, triggerTime] = s.inputSingleScan()
```

Description

data = s.inputSingleScan() immediately acquires a single scan from
each input channel in the session and returns the data acquired in a
1xn array of doubles, data, where n is the number of input channels
in the session.

[data, triggerTime] = s.inputSingleScan() immediately acquires a single scan from each input channel in the session and returns the data acquired in a 1xn array of doubles, data, where n is the number of input channels in the session. triggerTime is a MATLAB serial date time stamp representing the time the data is acquired.

Tips

To acquire more than a single input, use daq.Session.startForeground.

Output Arguments

data

A 1xn array of doubles, where n is the number of input channels in the session.

triggerTime

A MATLAB serial date time stamp representing the time the data is acquired.

Examples

Analog Input

Acquire a single input from an analog channel.

1

Create a session and add two analog input channels:

```
s = daq.createSession('ni');
s.addAnalogInputChannel('cDAQ1Mod1', 1:2, 'Voltage');
```

```
2
  Input a single scan:
   data = s.inputSingleScan
  data =
     -0.1495
              0.8643
  Digital Input
  Acquire a single input from a digital channel.
1
  Create a session and add two digital channels with InputOnly
  measurement type:
  s = daq.createSession('ni');
  s.addDigitalChannel('dev1', 'Port0/Line0:1', 'InputOnly');
2
  Input a single scan:
   data = s.inputSingleScan
  Counter Input
  Acquire a single input from a counter channel.
1
  Create a session and add a counter input channel with EdgeCount
  measurement type:
  s = dag.createSession('ni');
  s.addCounterInputChannel('Dev1',0,'EdgeCount');
2
```

daq.Session.inputSingleScan

Input a single edge count:

data = s.inputSingleScan

See Also

daq.createSession | daq.Session.addAnalogInputChannel
| daq.Session.addCounterInputChannel |
daq.Session.addDigitalChannel | daq.Session.outputSingleScan
| daq.Session.startForeground

daq.Session.queueOutputData

Purpose

Queue data to be output

Syntax

s.queueOutputData (data)

Description

s.queueOutputData (data) queues data to be output. When using analog output channels, you must queue data before you call daq.Session.startForeground or daq.Session.startBackground.

Input Arguments

data

An mxn matrix of doubles where m is the number of scans to generate, and n is the number of output channels in the session.

Examples

Create a session, add an analog output channel, and queue some data to output:

```
s = daq.createSession('ni');
s.addAnalogOutputChannel('cDAQ1Mod2', 'ao0', 'Voltage');
s.queueOutputData (linspace(-1, 1, 1000)');
s.startForeground;
```

Queue output data for multiple channels:

```
s = daq.createSession('ni');
s.addAnalogOutputChannel('cDAQ1Mod2', 0:1, 'Voltage');
data0 = linspace(-1, 1, 1000)';
data1 = linspace(-2, 2, 1000)';
s.queueOutputData ([data0 data1]);
s.startBackground
```

See Also

daq.createSession | | daq.Session.addAnalogOutputChannel |
daq.Session.startBackground | daq.Session.startForeground

daq.Session.outputSingleScan

Purpose Generate single scan on all output channels

Syntax s.outputSingleScan(data)

Description s.outputSingleScan(data) outputs a single scan of data on one or

more analog output channels.

Input Arguments

data

A 1xn matrix of doubles where n is the number of output channels in the session.

Examples

Analog Output

Output a single scan on two analog output voltage channels

1

Create a session and add two analog output channels.

```
s = daq.createSession('ni');
s.addAnalogOutputChannel('cDAQ1Mod2', 0:1, 'Voltage');
```

2

Create an output value and output a single scan for each channel added.

```
s.outputSingleScan([1.5 4]);
```

Digital Output

Output one value each on 2 lines on a digital channel

1

Create a session and add two digital channels from port 0 that measures output only:

```
s = daq.createSession('ni');
s.addDigitalChannel('dev1', 'Port0/Line0:1', 'OutputOnly')
```

daq.Session.outputSingleScan

2

Output one value each on the two lines:

s.outputSingleScan([0 1])

See Also

daq.createSession | daq.Session.addAnalogOutputChannel | |
daq.Session.addDigitalChannel | daq.Session.outputSingleScan
| daq.Session.inputSingleScan

DataAvailable Event

Purpose

Notify when acquired data is available to process

Syntax

```
lh = session.addlistener('DataAvailable',callback);
lh = session.addlistener('DataAvailable',
@(src, event) expr)
```

Description

1h = session.addlistener('DataAvailable',callback); creates
a listener for the DataAvailable event. When data is available to
process, the callback is executed. The callback can be any MATLAB
function with the (src, event) signature.

1h = session.addlistener('DataAvailable', @(src, event)
expr) creates a listener for the DataAvailable event and fires an
anonymous callback function. The anonymous function requires the
specified input arguments and executes the operation specified in
the expression expr. Anonymous functions provide a quick means of
creating simple functions without storing your function to a file. For
more information see Anonymous Functions.

The callback has two required parameters: src and event. src is the session object for the listener and event is a daq.DataAvailableInfo object containing the data associated and timing information. Properties of daq.DataAvailableInfo are:

Data

An mxn matrix of doubles where m is the number of scans acquired, and n is the number of input channels in the session.

TimeStamps

The timestamps relative to TriggerTime in an mx1 array where m is the number of scans acquired.

TriggerTime

A MATLAB serial date time stamp representing the absolute time the acquisition trigger occurs.

Tips

 Frequency with which the DataAvaialble event is fired, is controlled by NotifyWhenDataAvailableExceeds

Examples

Create a session, add an analog input channel, and change the duration of the acquisition:

See Also

daq.Session.addlistener | daq.Session |
daq.Session.startBackground | NotifyWhenDataAvailableExceeds
| IsNotifyWhenDataAvailableExceedsAuto

DataRequired Event

Purpose

Notify when additional data is required for output on continuous generation

Syntax

```
1h = session.addlistener('DataRequired',callback);
1h = session.addlistener('DataRequired',
@(src,event) expr);
```

Description

1h = session.addlistener('DataRequired',callback); creates a
listener for the DataRequired event. When more data is required, the
callback is executed. The callback can be any MATLAB function with
the (src, event) signature.

1h = session.addlistener('DataRequired', @(src,event)
expr); creates a listener for the DataRequired event and fires an
anonymous function. The anonymous function requires the specified
input arguments and executes the operation specified in the expression
expr. Anonymous functions provide a quick means of creating simple
functions without storing your function to a file. For more information
see Anonymous Functions.

The callback has two required parameters:. src and event. src is the session object for the listener and event is a daq.DataRequiredInfo object.

Tips

- The callback is typically used to queue more data to the device.
- Frequency is controlled by NotifyWhenScansQueuedBelow.

Examples

Add an anonymous listener to a signal generation session by first:

- Creating a session
- Adding two analog output channels

```
s = daq.createSession('ni');
s.IsContinuous = true
s.addAnalogOutputChannel('cDAQ1Mod2', 0:1, 'Voltage');
```

Create output data for the two channels:

DataRequired Event

```
outputData0 = (linspace(-1, 1, 1000))';
outputData1 = (linspace(-2, 2, 1000))';

Queue the output data and add an anonymous listener and generate the signal in the background:
s.queueOutputData([outputData0, outputData1]);
lh = s.addlistener('DataRequired', ...
    @(src,event) src.queueOutputData([outputData0, outputData1]));

Generate data and pause for up to 15 seconds:
s.startBackground();
pause (15)

Delete the listener:
delete (lh)

daq.Session.addlistener | IsContinuous | daq.Session |
NotifyWhenScansQueuedBelow | IsNotifyWhenScansQueuedBelowAuto
```

See Also

Purpose

Notify when device-related errors occur

Syntax

```
1h = session.addlistener('ErrorOccurred', callback);
1h = session.addlistener('ErrorOccurred' @(src,event)
expr);
```

Description

1h = session.addlistener('ErrorOccurred', callback); creates a
listener for the ErrorOccurred event. When an error occurs, the call
back is executed. The callback can be any MATLAB function with the
(src, event) signature.

1h = session.addlistener('ErrorOccurred' @(src,event)
expr); creates a listener for the ErrorOccurred event and fires an
anonymous function. The anonymous function requires the specified
input arguments and executes the operation specified in the expression
expr. Anonymous functions provide a quick means of creating simple
functions without storing your function to a file. For more information,
see Anonymous Functions

The callback has two required parameters: src and event. src is the session object for the listener and event is a daq.ErrorOccurredInfo object. The daq.ErrorOccurredInfo object contains the Error property, which is the MException associated with the error. You could use the MException.getReport method to return a formatted message string that uses the same format as errors thrown by internal MATLAB code.

Examples

Create a session, and add an analog input channel:

```
s = daq.createSession('ni');
s.addAnalogInputChannel('cDAQ1Mod1', 'ai0', 'Voltage');
To get a formatted report of the error, type:

lh = s.addlistener('ErrorOccurred' @(src,event) disp(event.Error.getReport()));
Acquire data, wait and delete the listener:
s.startBackground();
```

ErrorOccurred Event

```
s.wait()
delete(lh)
```

See Also

daq.Session.addlistener | daq.Session.startBackground |
daq.Session | MException

daq.Session.addCounterInputChannel

Purpose

Add counter input channel

Syntax

- s.addCounterInputChannel(deviceID, channelID,
 measurementType)
- ch = s.addCounterInputChannel(deviceID, channelID,
 measurementType)
- [ch, idx] = s.addCounterInputChannel(deviceID, channelID,
 measurementType)

Description

- s.addCounterInputChannel(deviceID, channelID, measurementType) adds a counter channel on the device represented by deviceID, with the specified channelID, and channel measurement type, represented by measurementType, to the session s. Measurement types are vendor specific.
- ch = s.addCounterInputChannel(deviceID, channelID,
 measurementType) returns the object ch.
- [ch, idx] = s.addCounterInputChannel(deviceID, channelID, measurementType) returns the object ch, representing the channel that was added and the index, idx, which is an index into the array of the session object's Channels property.

Tips

Use daq.createSession to create a session object before you use this method.

Input Arguments

deviceID

Specify the vendor-defined ID of the device. The specified channel is created for this device. Obtain by calling daq.getDevices.

channelID

Specify the ID of the channel added to the device. You can also add a range of channels.

For NI devices, use either a terminal name, like 'ctr2', or a numeric equivalent like 2 for the channel ID.

measurementType

daq.Session.addCounterInputChannel

Specify a string that represents a vendor-defined measurement type. Measurement types include:

- 'EdgeCount'
- 'PulseWidth'
- 'Frequency'
- 'Position'

To see a list of all measurement types supported by a device, type daq.getDevices and then click on the device in the list.

Output Arguments

ch

Object representing the channel you added.

idx

An index into the array of the session object's Channels property.

Properties

ActiveEdge Rising or falling edges of

EdgeCount signals

Active pulse measurement of

PulseWidth counter channel

CountDirection Specify direction of counter

channel

Device annel device

EncoderType Encoding type of counter channel

ID ID

InitialCount Specify initial count point

MeasurementType Type counter channel

measurement

daq.Session.addCounterInputChannel

Name Specify descriptive name for the

channel

Terminal PFI terminal

ZResetCondition Reset condition for Z-indexing
ZResetEnable Enable reset for Z-indexing
ZResetValue Reset value for Z-indexing

Examples

Add a counter input EdgeCount channel:

```
s = daq.createSession ('ni')
s.addCounterInputChannel('cDAQ1Mod5','ctr0','EdgeCount')
```

Add a counter input Frequency channel. Specify output arguments to represent the channel object and the index:

```
s = daq.createSession ('ni')
[ch, idx] = s.addCounterInputChannel('cDAQ1Mod5',1,'Frequency')
```

Add counter input channels 0, 2, and 4:

```
s = daq.createSession ('ni')
s.addCounterInputChannel('cDAQ1Mod5',[0 2 4], 'EdgeCount');
```

See Also

daq.Session.inputSingleScan | daq.Session.startForeground | daq.Session.startBackground | daq.Session.removeChannel | daq.Session

How To

daq.Session.addCounterOutputChannel

Purpose

Add counter output channel

Syntax

- s.addCounterOutputChannel(deviceID, channelID,
 measurementType)
- ch = s.addCounterOutputChannel(deviceID, channelID,
 measurementType)
- [ch, idx] = s.addCounterOutputChannel(deviceID, channelID,
 measurementType)

Description

- s.addCounterOutputChannel(deviceID, channelID, measurementType) adds a counter output channel on the device represented by deviceID, with the specified channelID, and channel measurement type, defined by measurementType, on the session object, s. Measurement types are vendor specific.
- ch = s.addCounterOutputChannel(deviceID, channelID,
 measurementType) returns the object ch, representing the channel that
 was added.
- [ch, idx] = s.addCounterOutputChannel(deviceID, channelID, measurementType) returns the object ch, representing the channel that was added and the index, idx, which is an index into the array of the session object's Channels property.

Tips

Use daq.createSession to create a session object before you use this method.

Input Arguments

deviceID

Specify the vendor-defined ID of the device. The specified channel is created for this device. The device ID is the ID of the device that you obtain by calling daq.getDevices.

channelID

Specify the ID of the channel added to the device. You can also add a range of channels.

daq.Session.addCounterOutputChannel

For NI devices, use either a terminal name, like 'ctr2', or a numeric equivalent like 2, for the channel ID.

measurementType

Specify a string that represents a vendor-defined measurement type. A valid output measurement type is 'PulseGeneration'.

Output Arguments

ch

Object representing the channel you added.

idx

An index into the array of the session object's Channels property.

Properties

Device annel device

DutyCycle Duty cycle of counter output

channel

Frequency Frequency of generated pulses on

counter output channel

ID ID

IdleState Default state of counter output

channel

InitialDelay Delay until output channel

generates pulses

MeasurementType Type counter channel

measurement

Name Specify descriptive name for the

channel

Examples

Add an counter output PulseGeneration channel:

```
s = daq.createSession ('ni')
```

s.addCounterOutputChannel('cDAQ1Mod3','ctr0','PulseGeneration')

daq.Session.addCounterOutputChannel

Add two counter output PulseGeneration channels:

```
s = daq.createSession ('ni')
s.addCounterOutputChannel('cDAQ1Mod3',0:1,'PulseGeneration')
```

See Also

```
daq.session.startBackground |
daq.Session.addCounterInputChannel |
daq.Session.removeChannel | daq.Session
```

How To

daq.Session.resetCounters

Purpose

Reset counter channel to initial count

Syntax

s.resetCounters()

Description

s.resetCounters() restarts the current value of counter channels configured in the session object, s to the specified InitialCount property on each channel.

Tips

- Reset counters only if you are performing on-demand operations using daq.Session.inputSingleScan or daq.Session.outputSingleScan
- Create an acquisition session and add a channel before you use this method. See daq.createSession for more information.
- Use daq.Session.startCounters to start the counter channels.
- Use daq.Session.stopCounters to stop the counter channels.

Examples

1 Create a session with a counter channel with an 'EdgeCount' measurement type:

```
s = daq.createSession ('ni');
s.addCounterInputChannel('cDAQ1Mod5', 0, 'EdgeCount');
```

2 Acquire data:

```
s.inputSingleScan
ans =
    756
```

3 Reset the counter to the default value, 0, and acquire again:

```
s.resetCounters
s.inputSingleScan
ans =
```

daq.Session.resetCounters

303

See Also daq.createSession | | | daq.Session

Tutorials • "Counter Channels"

How To ·

daq.Session.startCounters

Purpose Start counter channels **Syntax** s.startCounters() **Description** s.startCounters() starts counter channels configured in the session object, s. **Tips** • Create an acquisition session and add a channel before you use this method. See dag.createSession for more information. **Examples** 1 Create a session with a counter channel with an 'EdgeCount' measurement type: s = daq.createSession ('ni'); s.addCounterInputChannel('cDAQ1Mod5', 0, 'EdgeCount'); **2** Start the counter and acquire data: s.startCounters; s.inputSingleScan ans = 756 See Also daq.createSession | daq.Session.startForeground | dag.Session.addCounterInputChannel | dag.Session.addCounterOutputChannel | dag.Session How To "Using the Session-Based Interface"

daq.Session.stopCounters

Purpose Stop counter channels

Syntax

s.stopCounters()

Description

s.stopCounters() stops counter channels configured in the session object, s.

Tips

- Create an acquisition session and add a channel before you use this method. See daq.createSession for more information.
- Use daq.Session.resetCounter to restart the count from the channels initial count.

Examples

Add an counter input EdgeCount channel:

```
s = daq.createSession ('ni')
s.addCounterInputChannel('cDAQ1Mod5','ctr0','EdgeCount')
```

Start the counters and input a single scan:

s.startCounters

s.inputSingleScan

ans =

133

Stop the counter:

s.stopCounters;

See Also

```
daq.createSession | daq.Session.startCounters
| daq.Session.addCounterOutputChannel |
daq.Session.addCounterInputChannel | daq.Session
```

daq.Session.stopCounters

How To• "Using the Session-Based Interface"

•

Purpose

Add trigger connection

Syntax

s.addTriggerConnection(source,destination,type)
tc = s.addTriggerConnection(source,destination,type)
[tc,idx] = s.addTriggerConnection(source,destination,type)

Description

s.addTriggerConnection(source, destination, type) establishes a trigger connection from the specified source device and terminal to the specified destination device and terminal, of the specified connection type.

tc = s.addTriggerConnection(source, destination, type) establishes a trigger connection from the specified source and terminal to the specified destination device and terminal, of the specified connection type and displays it in the variable tc.

[tc,idx] = s.addTriggerConnection(source,destination,type) establishes a trigger connection from the specified source device and terminal to the specified destination device and terminal of the specified connection type and displays the connection in the variable tc and the connection index. idx.

Tips

 Before adding trigger connections, create a session using daq.createSession, and add channels to the session.

Input Arguments

source

Specify a source for the trigger connection. Valid values are:

'external'

When your trigger is based on an external event.

'deviceID/terminal'

When your trigger source is on a specific terminal on a device in your session. For example, 'Dev1/PFI1', for more information on device ID see Device. For more information on terminal see Terminals.

'chassisId/terminal'

When your trigger source is on a specific terminal on a chassis in your session, for example, 'cDAQ1/PFI1'. For more information on terminal see Terminals.

You can have only one trigger source in a session.

destination

Specify a destination for the trigger connection. Valid values are:

'external'

When your trigger source is connected to an external device.

'deviceID/terminal'

When your trigger source is connected to another device in your session, for example, 'Dev1/PFI1'. For more information on device ID see Device. For more information on terminal see Terminals.

'chassisId/terminal'

When your trigger source is connected to a chassis in your session, for example, 'cDAQ1/PFI1'. For more information on terminal see Terminals.

You can also specify multiple destination devices as an array, for example, {'Dev1/PFI1', 'Dev2/PFI1'}.

type

Specify the trigger connection type. StartTrigger is the only connection type available for trigger connections at this time.

Output Arguments

tc

Object representing the trigger connection you added.

idx

An index into the array of the session object's Connections property.

ndicates	trigger	destination
]	ndicates	ndicates trigger

terminal

ExternalTriggerTimeout Indicate if external trigger timed

out

IsWaitingForExternalTrigger Indicates if synchronization is

waiting for an external trigger

Source Indicates trigger source terminal

Terminals available on device or

CompactDAQ chassis

TriggerCondition Specify condition that must be

satisfied before trigger executes

TriggersPerRun Indicate the number of times the

trigger executes in an operation

TriggersRemaining Indicates the number of trigger to

execute in an operation

TriggerType Type of trigger executed

Examples Add External Start Trigger Connection

Create a session and add an analog input channel from Dev1 to the session.

```
s=daq.createSession('ni')
s.addAnalogInputChannel('Dev1','ai0', 'Voltage');
```

Add a trigger connection from an external device to terminal PFI1 on Dev1 using the 'StartTrigger' connection type:

```
s.addTriggerConnection('external','Dev1/PFI1','StartTrigger')
```

Export Trigger to External Device

To Add trigger connection going to an external destination, create a session and add an analog input channel from Dev1 to the session.

```
s=daq.createSession('ni')
s.addAnalogInputChannel('Dev1','ai0','Voltage');
```

Add a trigger from terminal PFIO on Dev1 to an external device using the 'StartTrigger' connection type:

```
s.addTriggerConnection('Dev1/PFI1','external','StartTrigger')
```

Save Trigger Connection in Variables

To display a trigger connection in a variable, create a session and add an analog input channel from Dev1 and Dev2 to the session.

```
s=daq.createSession('ni')
s.addAnalogInputChannel('Dev1','ai0','Voltage');
s.addAnalogInputChannel('Dev2','ai1','Voltage');
```

Add a trigger connection from terminal PFI1 on Dev1 to terminal PFI0 on Dev2 using the 'StartTrigger' connection type and store it in to

```
tc = s.addTriggerConnection('Dev1/PFI1','Dev2/PFI0','StartTrigger');
```

See Also

daq.createSessiondaq.Session.addClockConnection |
daq.Session.removeConnection |

Related Examples

- "Multiple-Device Synchronization"
- "Multiple-Chassis Synchronization"

${\bf daq. Session. add Trigger Connection}$

Concepts

Purpose

Add clock connection

Syntax

s.addClockConnection(source,destination,type)
cc = s.addClockConnection(source,destination,type)
[cc,idx] = s.addClockConnection(source,destination,type)

Description

s.addClockConnection(source, destination, type) adds a clock connection from the specified source device and terminal to the specified destination device and terminal, of the specified connection type.

cc = s.addClockConnection(source,destination,type) adds a clock connection from the specified source device and terminal to the specified destination device and terminal, of the specified connection type and displays it in the variable cc.

[cc,idx] = s.addClockConnection(source,destination,type) adds a clock connection from the specified source device and terminal to the specified destination device and terminal, of the specified connection type and displays the connection in the variable cc and the connection index. idx.

Tips

• Before adding clock connections, create a session using daq.createSession, and add channels to the session.

Input Arguments

source

Specify a source for the clock connection. Valid values are:

'external'

When your clock is based on an external event.

'deviceID/terminal'

When your clock source is on a specific terminal on a device in your session, for example, 'Dev1/PFI1'. For more information on device ID see Device. For more information on terminal see Terminals.

'chassisId/terminal'

When your clock source is on a specific terminal on a chassis in your session, for example, 'cDAQ1/PFI1'. For more information on terminal see Terminals.

You can have only one clock source in a session.

destination

Specify a destination for the clock connection. Valid values are:

'external'

When your clock source is connected to an external device.

'deviceID/terminal'

When your clock source is connected to another device in your session, for example, 'Dev1/PFI1'. For more information on device ID see Device. For more information on terminal see Terminals.

'chassisId/terminal'

When your clock source is connected to a chassis in your session, for example, 'cDAQ1/PFI1'. For more information on terminal see Terminals.

You can also specify multiple destination devices as an array, for example, {'Dev1/PFI1', 'Dev2/PFI1'}.

type

Specify the clock connection type. ScanClock is the only connection type available for clock connections at this time.

Output Arguments

CC

Object representing the clock connection you added.

idx

An index into the array of the session object's Connections property.

Properties

Destination Indicates trigger destination

terminal

Source Indicates trigger source terminal

Terminals available on device or

CompactDAQ chassis

Examples Add External Scan Clock

Create a session and add an analog input channel from Dev1 to the session.

```
s=daq.createSession('ni')
s.addAnalogInputChannel('Dev1','ai0', 'Voltage');
```

Add a clock connection from an external device to terminal PFI1 on Dev1 using the 'ScanClock' connection type:

```
s.addClockConnection('external', 'Dev1/PFI1', 'ScanClock')
```

Export Scan Clock to External Device

To add clock connection going to an external destination, create a session and add an analog input channel from Dev1 to the session.

```
s=daq.createSession('ni')
s.addAnalogInputChannel('Dev1','ai0', 'Voltage');
```

Add a clock from terminal PFIO on Dev1 to an external device using the 'ScanClock' connection type:

```
s.addClockConnection('Dev1/PFI1','external','ScanClock')
```

Save Clock Connection in Variables

To display a clock connection in a variable, create a session and add an analog input channel from Dev1 and Dev2 to the session.

```
s=daq.createSession('ni')
s.addAnalogInputChannel('Dev1','ai0','Voltage');
s.addAnalogInputChannel('Dev2','ai1','Voltage');
```

Add a clock connection from terminal PFI1 on Dev1 to terminal PFI0 on Dev2 with 'ScanClock' connection type, storing it in cc:

```
cc = s.addClockConnection('Dev1/PFI1', 'Dev2/PFI0', 'ScanClock');
```

See Also

daq.createSessiondaq.Session.addTriggerConnection |
daq.Session.removeConnection |

Related Examples

- - "Multiple-Device Synchronization"
 - "Multiple-Chassis Synchronization"

Concepts

•

daq.Session.removeConnection

Purpose

Remove clock or trigger connection

Syntax

s.removeConnection(idx);

Description

s.removeConnection(idx); removes the specified clock or trigger with the index, idx, from the session. The connected device remains in the session, but no longer synchronizes with other connected devices in the session.

Input Arguments

idx

Index of the connection you want to remove.

Examples

Create clock and trigger connections in session s:

```
s = daq.createSession('ni');
s.addAnalogInputChannel('Dev1','ai0','Voltage')
s.addAnalogInputChannel('Dev2','ai0','Voltage')
s.addAnalogInputChannel('Dev3','ai0','Voltage')
s.addTriggerConnection('Dev1/PFI0',{'Dev2/PFI0','Dev3/PFI0'}','StartTrigger');
s.addClockConnection('Dev1/PFI1',{'Dev2/PFI1','Dev3/PFI1'}'ScanClock');
```

View existing synchronization connections in the session:

s.Connections

```
Synchronization summary:

Start Trigger is provided by 'Dev1' on terminal 'PFIO' and will be received.

Scan Clock is provided by 'Dev1' on terminal 'PFI1' and will be received.

index ConnectionType Source Destination.
```

```
1 StartTrigger Dev1/PFI0 Dev2/PFI0
2 StartTrigger Dev1/PFI0 Dev3/PFI0
3 ScanClock Dev1/PFI1 Dev2/PFI1
```

daq.Session.removeConnection

4 ScanClock Dev1/PFI1 Dev3/PFI1

Remove the trigger connection with the index 2 from Dev3/PFI0 to Dev1/PFI0:

s.removeConnection(2);

View updated connections:

s.Connections

See Also

daq.createSessiondaq.Session.addClockConnection |
daq.Session.addTriggerConnection |

Concepts

- •
- •
- •

Purpose

Add digital channel

Syntax

s.addDigitalChannel(deviceID, channelID, measurementType)
ch =

s.addDigitalChannel(deviceID,channelID,measurementType)
[ch, idx] = s.addDigitalChannel(deviceID,channelID,
 measurementType)

Description

s.addDigitalChannel(deviceID, channelID, measurementType) adds digital channel on the device represented by deviceID, with the specified port and single line combination and the channel measurement type to the session, s.

ch =

s.addDigitalChannel(deviceID, channelID, measurementType) creates and displays the digital channel ch.

[ch, idx] =

s.addDigitalChannel(deviceID, channelID, measurementType) creates and displays the digital channel ch that was added and the index, idx, which is an index into the array of the session object's Channels property.

Note To input and output decimal values, use the conversion functions:

- decimalToBinaryVector
- binaryVectorToDecimal
- hexToBinaryVector
- binaryVectorToHex

Tips

- Create a session using daq.createSession before adding digital channels.
- Change the Direction of a bidirectional channel before you read or write digital data.

Input Arguments

deviceID - Specify the vendor-defined ID of the device.

A channel is created for the specified device. To obtain the device ID, call daq.getDevices.

channelID - Specify the port and the line you want to use on this device.

Add the channel with both the port and the line specified, for example 'port0/line0'. You can add a port and line or a range of lines on a port that you will use on this channel. To add more than one line, use 'port0/line0:3' or 'port1/line0, port0/line0, port2/line1'.

measurementType - Specify a string that represents the channel's measurement type. Digital channels can be:

- InputOnly
- OutputOnly
- Bidirectional

Output Arguments

ch - Object representing the digital channels added in the session. Use this object to access properties on the channel.

Examine the channels added to the session using the

```
s = daq.createSession('ni')
[ch, idx]=s.addDigitalChannel('dev1', 'Port0/Line0:1', 'InputOnly
```

idx - An index into the array of the session object's Channels property. The index displays the position of the channel in the session. If you add a channel with channel ID 'port0/line0' as the first channel in your session, the session index will be 1. Use the index to access the Channels properties.

Examine the direction of two digital channels:

```
s = daq.createSession('ni')
[ch, idx]=s.addDigitalChannel('dev1', 'Port0/Line0:1', 'InputOnly')
s.Channels(idx).Direction
```

Properties

Device annel device

Direction Specify digital channel direction

ID ID

Name Specify descriptive name for the

channel

Examples Add Digital Channels

Discover available digital devices on your system, create a session with digital channels.

Find all installed devices and get detailed subsystem information for NI USB-6255:

d=daq.getDevices

d =

Data acquisition devices:

```
2
                      National Instruments USB-6363
   ni Dev2
d(1)
ans =
ni: National Instruments USB-6255 (Device ID: 'Dev1')
  Analog input subsystem supports:
      7 ranges supported
      Rates from 0.1 to 1250000.0 scans/sec
      80 channels ('ai0' - 'ai79')
      'Voltage' measurement type
   Analog output subsystem supports:
      -5.0 to +5.0 Volts, -10 to +10 Volts ranges
      Rates from 0.1 to 2857142.9 scans/sec
      2 channels ('ao0', 'ao1')
      'Voltage' measurement type
   Digital subsystem supports:
      24 channels ('port0/line0' - 'port2/line7')
      'InputOnly','OutputOnly','Bidirectional' measurement types
   Counter input subsystem supports:
      Rates from 0.1 to 80000000.0 scans/sec
      2 channels ('ctr0', 'ctr1')
      'EdgeCount', 'PulseWidth', 'Frequency', 'Position' measurement type
   Counter output subsystem supports:
      Rates from 0.1 to 80000000.0 scans/sec
      2 channels ('ctr0','ctr1')
      'PulseGeneration' measurement type
Create a session with input, output, and bidirectional channels using
Dev1:
```

s = daq.createSession('ni');

```
s.addDigitalChannel('dev1', 'Port0/Line2:3', 'OutputOnly');
                 s.addDigitalChannel('dev1', 'Port2/Line0:1', 'Bidirectional')
                 ans =
                 Data acquisition session using National Instruments hardware:
                    Clocked operations using startForeground and startBackground are disab
                    Only on-demand operations using inputSingleScan and outputSingleScan of
                   Number of channels: 6
                       index Type Device Channel MeasurementType Range Name
                          dio Dev1 portO/lineO InputOnly
                                                                          n/a
                            dio Dev1 portO/line1 InputOnly
                      2
                                                                           n/a
                            dio Dev1 port0/line2 OutputOnly
                      3
                                                                           n/a
                            dio Dev1 port0/line3 OutputOnly
                                                                           n/a
                      5
                            dio Dev1 port2/lineO Bidirectional (Unknown) n/a
                            dio Dev1 port2/line1 Bidirectional (Unknown) n/a
See Also
                 dag.Session.inputSingleScan | dag.Session.outputSingleScan
                 | daq.createSession | decimalToBinaryVector
                 | binaryVectorToDecimal | hexToBinaryVector |
                 binaryVectorToHex
Related
                 • "Acquire Non-Clocked Digital Data"
Examples
                 • "Generate Non-Clocked Digital Data"
Concepts
```

s.addDigitalChannel('dev1', 'Port0/Line0:1', 'InputOnly');

decimalToBinaryVector

Purpose

Convert decimal value to binary vector

Syntax

decimalToBinaryVector(decimalNumber)

decimalToBinaryVector(decimalNumber, numberOfBits)

decimalToBinaryVector(decimalNumber, numberOfBits, bitOrder)

decimalToBinaryVector(decimalNumber,[],bitOrder)

Description

decimalToBinaryVector(decimalNumber) converts positive decimal number to a binary vector, represented using the minimum number of bits.

decimalToBinaryVector(decimalNumber, numberOfBits) converts decimal number to a binary vector with the specified number of bits.

decimalToBinaryVector(decimalNumber, numberOfBits, bitOrder) converts decimal number to a binary vector with the specified number of bits in the specified bit ordering.

decimalToBinaryVector(decimalNumber,[],bitOrder) converts decimal number to a binary vector with default number of bits in the specified bit ordering.

Input Arguments

decimalNumber - Specify the number you want to convert to binary vector.

You can only specify positive integer scalars.

Data Types

single | double | int8 | int16 | int32 | int64 | uint8 | uint16 | uint32 | uint64

numberOfBits - specify the number of bits required to correctly represent the decimal number.

This is an optional argument. Specify the number of bits you want your decimal represented in. If you do not specify the number of bits, the number is represented using the minimum number of bits needed.

Default: If you do not specify the number of bits, the number is represented using the minimum number of bits needed.

bitOrder - Specify the bit order for the binary vector representation.

Specify MSBFirst if you want the first element of the output to contain the most significant bit of the decimal number. By default, the MSB is displayed first.

Specify LSBFirst if you want the first element of the output to contain the least significant bit of the decimal number.

Default: The default for this optional argument is MSBFirst.

Examples

Convert the decimal 6 to a binary vector

```
decimalToBinaryVector(6)
ans =
1 1 0
```

Convert the decimal 6 into a binary vector of specific bits

Convert the number 6 to an 8-bit binary vector with the most significant bit as the first element:

```
decimalToBinaryVector(6, 8, 'MSBFirst')
ans =
     0     0     0     0     1     1     0
```

Convert the decimal 6 into a binary vector with LSB first

Convert the number 6 to a binary vector with the least significant bit as the first element, using the default bit length:

```
decimalToBinaryVector(6, [], 'LSBFirst')
ans =
     0     1     1
```

See Also hexToBinaryVectorbinaryVectorToDecimalbinaryVectorToHex

Related Examples

• "Generate Signals Using Decimal Data"

binaryVectorToDecimal

Purpose

Convert binary vector value to decimal

Syntax

binaryVectorToDecimal(binaryVector)

binaryVectorToDecimal(binaryVector,bitOrder)

Description

binaryVectorToDecimal(binaryVector) converts binary vector to a

decimal.

binaryVectorToDecimal(binaryVector,bitOrder) converts binary

vector with the specified bit orientation.

Input Arguments

binary Vector - Specify the binary vector you want to convert to a decimal.

bitOrder - Specify the bit order for the binary vector representation.

Specify MSBFirst if you want the first element of the output to contain the most significant bit of the decimal number. By default, the MSB is displayed first.

Specify LSBFirst if you want the first element of the output to contain the least significant bit of the decimal number.

Default: The default for this optional argument is MSBFirst.

Examples

Convert the binary vector 1 1 0 to a decimal

binaryVectorToDecimal([1 1 0])

ans =

6

Convert a binary vector 1 1 0 with LSB first

binaryVectorToDecimal([1 1 0], 'LSBFirst')

binaryVectorToDecimal

ans =

6

 $\textbf{See Also} \quad \text{hexToBinaryVectordecimalToBinaryVectorbinaryVectorToHex}$

Related Examples

• "Generate Signals Using Decimal Data"

hexToBinaryVector

Purpose Convert hexadecimal value to binary vector

Syntax hexToBinaryVector(hexNumber)

hexToBinaryVector(hexNumber, numberOfBits)

hexToBinaryVector(hexNumber, numberOfBits, bitOrder)

Description hexToBinaryVector(hexNumber) converts hexadecimal number to

a binary vector.

 $\verb|hexToBinaryVector(hexNumber, numberOfBits)| converts \ hexadecimal$

number to a binary vector with the specified number of bits.

hexToBinaryVector(hexNumber,numberOfBits,bitOrder) converts hexadecimal number to a binary vector with the specified number of

bits in the specified bit ordering.

Input Arguments hexNumber - Specify the hexadecimal you want to convert to binary vector.

Specify a hexadecimal number as a character.

numberOfBits - specify the number of bits required to correctly represent the decimal number.

This is an optional argument. If you do not specify the number of bits, the number is represented using the minimum number of bits needed.

bitOrder - Specify the bit order for the binary vector representation.

Specify MSBFirst if you want the first element of the output to contain the most significant bit of the decimal number. By default, the MSB is displayed first.

Specify LSBFirst if you want the first element of the output to contain the least significant bit of the decimal number.

Default: The default for this optional argument is MSBFirst.

Examples

Convert the hexadecimal 'A1' to a binary vector

hexToBinaryVector('A1')

ans = 0 0 0 0

Convert the hexadecimal 'A1' number into a binary vector of specific bits

Convert 'A1' to an 12-bit binary vector with the most significant bit as the first element:

hexToBinaryVector('A1',12, 'MSBFirst')

ans =

Convert the hexadecimal 'A1' into a binary vector with LSB

0

first

Convert the number 'A1' to a binary vector with the least significant bit as the first element, using the default bit length:

hexToBinaryVector('A1', [], 'LSBFirst')

ans = 0 0 0 1 1

0 0 1

See Also

 ${\tt decimalToBinaryVectorbinaryVectorToDecimal binaryVectorToHex}$

Related **Examples**

• "Acquire Digital Data In Hexadecimal Values"

0

binaryVectorToHex

Purpose Convert binary vector value to hexadecimal

Syntax binaryVectorToHex(binaryVector)

binaryVectorToHex(binaryVector,bitOrder)

Description binaryVectorToHex(binaryVector) converts input binary vector to

a hexadecimal.

 $\verb|binary| VectorToHex(binary| Vector, \verb|bitOrder|) converts input binary|$

vector with the specified bit orientation.

Input Arguments

binaryVector - Specify the binary vector you want to convert to a hexadecimal. The binary vector must be a row vector with Os and

bitOrder - Specify the bit order for the binary vector representation.

Specify MSBFirst if you want the first element of the output to contain the most significant bit of the decimal number. By default, the MSB is displayed first.

Specify LSBFirst if you want the first element of the output to contain the least significant bit of the decimal number.

Default: The default for this optional argument is MSBFirst.

Examples Convert the binary vector 1 0 1 0 0 0 0 1 to a hexadecimal

binaryVectorToHex([0 0 1 1 1 1 0 1])

ans =

3D

Convert a binary vector 1 0 1 0 0 0 0 1 with LSB first

Convert the binary vector to a hexadecimal with the least significant bit first:

```
binaryVectorToHex([0 0 1 1 1 1 0 1], 'LSBFirst')
ans =
    BC
```

See Also hexToBinaryVectorbinaryVectorToDecimaldecimalToBinaryVector

Related Examples

• "Acquire Digital Data In Hexadecimal Values"

binaryVectorToHex

Index

A	В
ac	BiDirectionalBit property 4-2
coupling 4-6	binvec2dec function 8-29
acquiring data	BitsPerSample property 4-4
single point 8-72	block
ActiveEdge property	Analog Input 5-2
Session 2-2	Analog Input (Single Sample) 5-12
ActivePulse property 2-4	Analog Output 5-17
ADCTimingMode property	Analog Output (Single Sample) 5-25
Session 2-6	Digital Input 5-29
addchannel function 8-2	Digital Output 5-34
addline function 8-8	BridgeMode
addmuxchannel function 8-12	session-based
Advantech hardware	property 2-9
properties 3-3	buffer
AMUX-64T	configuration 2-10
adding channels 8-12	BufferingConfig property 2-10
channel indices 8-98	BufferingMode property 2-13
Analog Input (Single Sample) block 5-12	
Analog Input block 5-2	C
analog input object	•
acquisition	channel properties
single point 8-72	AI object 1-6
properties	AO object 1-10
channel 1-6	Channel property 2-15
common 1-3	ChannelName property 2-17
Analog Output (Single Sample) block 5-25	Channels property
Analog Output block 5-17	Session 2-19
analog output object	ChannelSkew property 2-21
output	ChannelSkewMode property 2-22
single point 8-112	cleaning up the MATLAB environment
properties	daqfind function 8-54
channel 1-10	clear function 8-30
common 1-8	ClockSource property 2-25
analogoutput function 8-22	common properties
AutoSyncDSA	AI object 1-3
property 2-8	AO object 1-8
	DIO object 1-12
	Connections property
	Session 2-29

CountDirection property 2-31	disp function 8-60
Coupling property 4-6	DMA
	NI hardware 4-24
D	DurationInSeconds property 2-42
dag.createSession function 8-146	DutyCycle property 2-43
dag.getDevices function 8-149	
daq.getVendors function 8-152	E
dagcallback function 8-32	EncoderType property 2-45
dagfind function 8-34	EventLog property 2-47
daghelp function 8-37	ExcitationCurrent property 2-50
daghwinfo function 8-40	ExcitationSource property 2-51
dagmem function 8-44	ExcitationVoltage property 2-52
dagread function 8-47	external clock
dagregister function 8-48	clock sources 2-25
dagreset function 8-50	ExternalClockDriveLine property 4-8
DataMissedFcn property 2-33	ExternalClockSource property 4-9
dc	ExternalSampleClockDriveLine property 4-10
coupling 4-6	ExternalSampleClockSource property 4-11
dec2binvec function 8-51	ExternalScanClockDriveLine property 4-13
DefaultChannelValue property 2-35	ExternalScanClockSource property 4-14
delete function 8-53	ExternalTriggerDriveLine property 4-15
Destination	ExternalTriggerTimeout
property 2-37	property 2-53
Device property 2-38	extracting data
device-specific properties	event information 8-71
Advantech hardware 3-3	native data 2-110
MCC hardware 3-4	
NI hardware 3-5	F
parallel port 3-7	-
sound cards 3-8	FIFO
digital I/O object	TransferMode 4-24
properties	finding device objects 8-34
common 1-12	flushdata $8\text{-}62$
line 1-12	flushdata function 8-62
Digital Input block 5-29	Frequency property 2-54
Digital Output block 5-34	full duplex
digitalio function 8-56	BitsPerSample property 4-4
Direction property 2-39 to 2-40	functions
disk logging 2-97	addchannel 8-2

addline 8-8	putsample 8-112
addmuxchannel 8-12	putvalue 8-114
analogoutput 8-22	save 8-116
binvec2dec 8-29	set 8-118
clear 8-30	setverify 8-121
daq.createSession 8-146	showdagevents 8-124
daq.getDevices 8-149	size 8-127
daq.getVendors $8\text{-}152$	softscope 8-130
daqcallback 8-32	start 8-139
daqfind 8-34	stop 8-141
daqhelp 8-37	trigger 8-143
daqhwinfo 8-40	wait 8-144
daqmem 8-44	
dagread 8-47	G
daqregister 8-48	
dagreset 8-50	get function 8-64
dec2binvec 8-51	getdata function 8-67
delete 8-53	getsample function 8-72
digitalio 8-56	getvalue function 8-74
disp 8-60	
flushdata 8-62	Н
get 8-64	hardware
getdata 8-67	initializing 8-50
getsample 8-72	holding the last output value 4-18
getvalue 8-74	HwChannel property 2-55
inspect 8-76	HwDigitalTriggerSource property 4-16
ischannel 8-78	HwLine property 2-57
isdioline 8-80	nweine property 2-97
islogging 8-82	_
isrunning $8-84$	I
issending 8-86	ID property 2-59
isvalid 8-88	IdleState property 2-61
length 8-91	Index property 2-63
load 8-93	InitialCount property 2-66
makenames 8-96	InitialDelay property 2-65
muxchanidx 8-98	initializing the hardware 8-50
obj2mfile $8\text{-}100$	InitialTriggerTime property 2-68
peekdata 8-103	InputOverRangeFcn property 2-70
propinfo 8-106	InputRange property 2-72
putdata 8-109	InputType property 2-75

inspect function 8-76	MeasurementType property 2-105
interrupts	mux board
NI hardware 4-24	adding channels 8-12
ischannel function 8-78	channel indices 8-98
IsContinuous property 2-77	muxchanidx function 8-98
isdioline function 8-80	
IsDone property 2-79	N
islogging function 8-82	N
IsLogging property 2-81	Name property 2-107 to 2-108
IsNotifyWhenDataAvailableExceedsAuto	National Instruments hardware
property 2-83	data transfer mechanisms 4-24
IsNotifyWhenScansQueuedBelowAuto	properties 3-5
property 2-84	native data
isrunning function 8-84	getdata 8-67
IsRunning property 2-85	offset 2-110
issending function 8-86	putdata 8-109
IsSimulated property 2-87	scaling 2-113
isvalid function 8-88	NativeOffset property 2-110
IsWaitingForExternalTrigger	NativeScaling property 2-113
property 2-89	NominalBridgeResistance property 2-115
• • •	NotifyWhenDataAvailableExceeds
L	property 2-116
_	NotifyWhenScansQueuedBelow property 2-119
length function 8-91	NumberOfScans property 2-120
Line property 2-90	NumMuxBoards property 4-17
LineName property 2-92	
load function 8-93	0
LogFileName property 2-94	
Logging property 2-95	obj2mfile function 8-100
LoggingMode property 2-97	OutOfDataMode property 4-18
LogToDiskMode property 2-99	OutputRange property 2-122
	outputting data
M	holding the last value 4-18
	single point 8-112
makenames function 8-96	
ManualTriggerHwOn property 2-101	P
maximum samples queued 2-103	parallel port
MaxSamplesQueued property 2-103	device-specific properties 3-7
Measurement Computing hardware	Parent property 2-124
properties 3-4	peekdata function 8-103
	F 10.100 1011 0 100

Port property 2-126 InputRange 2-72 PortAddress property 4-20 InputType 2-75 properties IsWaitingForExternalTrigger 2-89 ActiveEdge 2-2 Line 2-90 ADCTimingMode 2-6 LineName 2-92 AutoSyncDSA 2-8 LogFileName 2-94 BiDirectionalBit 4-2 Logging 2-95 BitsPerSample 4-4 LoggingMode 2-97 BridgeMode LogToDiskMode 2-99 property 2-9 ManualTriggerHwOn 2-101 MaxSamplesQueued 2-103 BufferingConfig 2-10 BufferingMode 2-13 Name 2-107 to 2-108 Channel 2-15 NativeOffset 2-110 ChannelName 2-17 NativeScaling 2-113 Channels 2-19 NumMuxBoards 4-17 ChannelSkew 2-21 OutOfDataMode 4-18 ChannelSkewMode 2-22 OutputRange 2-122 ClockSource 2-25 Parent 2-124 Connections 2-29 Port 2-126 Coupling 4-6 PortAddress 4-20 DataMissedFcn 2-33 Range 2-129 DefaultChannelValue 2-35 Rate 2-130 Destination 2-37 RateLimit 2-132 Direction 2-39 to 2-40RepeatOutput 2-133 EnhancedAliasRejectionEnable 2-46 Running 2-137 EventLog 2-47 RuntimeErrorFcn 2-139 ExternalClockDriveLine 4-8 SampleRate 2-141 ExternalClockSource 4-9 SamplesAcquired 2-143 ExternalSampleClockDriveLine 4-10 SamplesAcquiredFcn 2-144 ExternalSampleClockSource 4-11 SamplesAcquiredFcnCount 2-146ExternalScanClockDriveLine 4-13 SamplesAvailable 2-147 ExternalScanClockSource 4-14 SamplesOutput 2-149 ExternalTriggerDriveLine 4-15 SamplesOutputFcn 2-150 ExternalTriggerTimeout 2-53 SamplesOutputFcnCount 2-152 HwChannel 2-55 SamplesPerTrigger 2-153 HwDigitalTriggerSource 4-16 ScansAcquired 2-155HwLine 2-57ScansOutputByHardware 2-156 Index 2-63ScansQueued 2-157 InitialTriggerTime 2-68 Sending 2-158InputOverRangeFcn 2-70 SensorRange 2-161

ShuntLocation 2-163	AO object 1-10
ShuntResistance 2-164	common
Source 2-165	AI object 1-3
StandardSampleRates 4-22	AO object 1-8
StartFcn 2-166	DIO object 1-12
StopFcn 2-168	device-specific
Tag 2-170	Advantech hardware 3-3
TerminalConfig 2-173	MCC hardware 3-4
ThermocoupleType 2-176	NI hardware 3-5
Timeout $2 ext{-} 178$	parallel port 3-7
TimerFcn 2-180	sound cards 3-8
TimerPeriod $2 \text{-} 182$	line
TransferMode 4-24	DIO object 1-12
TriggerChannel 2-183	propinfo function 8-106
TriggerCondition 2-185 2-189	putdata function 8-109
TriggerConditionValue 2-190	putsample function 8-112
TriggerDelay 2-192	putvalue function 8-114
TriggerDelayUnits 2-194	
TriggerFcn 2-195	Q
TriggerRepeat 2-197	
TriggersExecuted 2-199	queuing data for output
TriggersPerRun 2-201	maximum number of samples 2-103
TriggersRemaining 2-202	
TriggerType $2\text{-}203\ 2\text{-}206$	R
Type 2-207 2-209	RO
Units $2\text{-}211$	session-based
session-based 2-210	property 2-128
UnitsRange 2-212	Range property
UserData 2-214	Session 2-129
Vendor $2\text{-}215$	Rate property
property	Session 2-130
R0	RateLimit property
session-based 2-128	Session 2-132
RTDConfiguration	RepeatOutput property 2-133
session-based 2-135	resetting the hardware 8-50
RTDType	RTDConfiguration
session-based 2-136	session-based
property types	property 2-135
channel	RTDType
AI object 1-6	J F -

session-based	ExcitationSource 2-51
property 2-136	ExcitationVoltage $2\text{-}52$
Running property 2-137	Frequency 2-54
RuntimeErrorFcn property 2-139	ID 2-59
	IdleState 2-61
S	InitialCount 2-66
	InitialDelay 2-65
SampleRate property 2-141	IsContinuous 2-77
SamplesAcquired property 2-143	IsDone 2-79
SamplesAcquiredFcn property 2-144	IsLogging 2-81
SamplesAcquiredFcnCount property 2-146	IsNotifyWhenDataAvailableExceedsAuto 2-83
SamplesAvailable property 2-147	IsNotifyWhenScansQueuedBelowAuto 2-84
SamplesOutput property 2-149	IsRunning 2-85
SamplesOutputFcn property 2-150	IsSimulated 2-87
SamplesOutputFcnCount property 2-152	MeasurementType 2-105
SamplesPerTrigger property 2-153	NominalBridgeResistance 2-115
save function 8-116	NotifyWhenDataAvailableExceeds 2-116
ScansAcquired property	NotifyWhenScansQueuedBelow 2-119
Session 2-155	NumberOfScans 2-120
ScansOutputByHardware property	Sensitivity $2 ext{-} 159$
Session 2-156	Terminal $2 ext{-} 172$
ScansQueued property	Terminals 2-174
Session 2-157	ZResetCondition 2-216
Sending property 2-158	ZResetEnable 2-217
Sensitivity property 2-159	ZResetValue 2-218
SensorRange property 2-161	set function 8-118
sensors	setverify function 8-121
range 2-161	showdagevents function 8-124
Session	ShuntLocation
channels 2-19 2-29 2-129 to 2-130 2-132	property 2-163
2-155 to 2-157 2-215	ShuntResistance
counter channels 2-2 2-6 2-176	property 2-164
Session properties	Simulink block
ActivePulse $2\text{-}4$	Analog Input 5-2
CountDirection 2-31	Analog Input (Single Sample) 5-12
Device 2-38	Analog Output 5-17 5-25
DurationInSeconds 2-42	Digital Input 5-29
DutyCycle 2-43	Digital Output 5-34
EncoderType $2\text{-}45$	single-point
ExcitationCurrent 2-50	acquisition 8-72

output 8-112 TriggerFcn property 2-195 size function 8-127 TriggerRepeat property 2-197 softscope function 8-130 triggers software clock session-based 8-201 MCC hardware 2-25 synchronizing for AI and AO 2-101 sound cards TriggersExecuted property 2-199 device-specific properties 3-8 TriggersPerRun standard sample rates 4-22 property 2-201 Source TriggersRemaining property 2-165 property 2-202 StandardSampleRates property 4-22 TriggerType property 2-203 2-206 start function 8-139 Type property 2-207 2-209 StartFcn property 2-166 stop function 8-141 U StopFcn property 2-168 Units synchronizing triggers 2-101 session-based property 2-210 T Units property 2-211 Tag property 2-170 UnitsRange property 2-212 Terminal property 2-172 UserData property 2-214 TerminalConfig property 2-173 Terminals property 2-174 ThermocoupleType property Vendor property Session 2-176 Session 2-215 Timeout property 2-178 TimerFcn property 2-180 TimerPeriod property 2-182 TransferMode property 4-24 wait function 8-144 trigger connection 8-201

Z

ZResetCondition property 2-216

ZResetEnable property 2-217

ZResetValue property 2-218

Index-8

trigger function 8-143

TriggerChannel property 2-183

TriggerDelay property 2-192
TriggerDelayUnits property 2-194

TriggerCondition property 2-185 2-189

TriggerConditionValue property 2-190